



ALADDIN HW3.0

User manual

Rev. 19 EN - 2019

Caution: Federal law restricts this device to sale by or on the order of an **optometrist**, **optician**, or an **ophthalmologist**.

The manufacturer has a policy of continuous improvement of its products, so it is possible that some instructions, specifications and pictures in this manual may differ slightly from the product you purchased. The manufacturer also reserves the right to make any changes to this manual without notice.

The original text of this manual is in English.

Accessibility and scope of the manual

Thank you for choosing this product.

Please read the information in this manual carefully. You must be familiar with its contents in order to work with the device.

Keep these instructions in a safe place close to the device. The manual must be at hand at all times. For a correct use of the instrument, read the instructions carefully.

The purpose of this manual is to inform the user as to all the device's functions, settings, safety, installation, maintenance, cleaning and storage instructions.

ALADDIN HW3.0

Product cod. 1240212

SW v.: 1.9.x



VISIA imaging S.r.l. Via Martiri della Libertà 95/e 52027 San Giovanni Valdarno (AR) Italy

Distributor

Topcon Europe Medical B.V. Essebaan 11 2908 LJ Capelle a/d IJssel Niederlande www.topcon.eu medical@topcon.eu

Contents

1	INT	INTENDED USE		
	1.1	Intended users	. 7	
	1.2	Places of use	. 7	
	1.3	Contraindications	. 7	
	1.4	Description of functionalities	. 7	
	1.5	Essential Performance	. 8	
	1.6	Interaction with the patient	. 8	
2	PRE	CAUTIONS	. 9	
	2.1	Electromagnetic compatibility	10	
	2.1	1 EM Emission	10	
	2.1	2 EM Immunity	10	
3	SYN	/BOLS	12	
	3.1	Sample labelling on the device	13	
4	SAF	ETY INSTRUCTIONS	14	
	4.1	General	14	
	4.2	Electrical safety	15	
	4.3	LED emission safety	15	
	4.4	Installation with external devices or IT Network	15	
	4.5	Transport and packaging	16	
	4.6	Cleaning	16	
	4.7	Checking the measurements	17	
	4.8	Privacy & Cybersecurity	18	
	4.8	1 Privacy	18	
	4.8	2 Privacy & Security options	20	
5	PRO	DDUCT WARRANTY AND RELIABILITY	21	
6	LEG	AL PROVISIONS	22	
7 MAIN COMPONENTS		23		
8	INS	TALLATION /UNINSTALLATION OF THE SYSTEM	24	
	8.1	Installing the system	24	
	8.2	Uninstalling the system	26	
9	STA	NDARD ACCESSORIES AND EQUIPMENTS	29	
1() SET	TING UP THE INSTRUMENT	30	
	10.1	Connection modes	30	

10.2	Pow	er on procedure	30
11 OP	ERATI	NG INSTRUCTIONS	31
11.1	Use	r Login/Logout	31
11.	1.1	User Logout	32
11.	1.2	Auto LogOFF/ Screen Saver	32
11.	1.3	Technician Support staff Login	32
11.2	Ove	rview	33
11.3	Che	cking the calibration	33
11.4	Pati	ent entry/selection	37
11.4	4.1	Creating a new patient	37
11.4	4.2	Selecting or modifying a patient	41
11.4	4.3	Selecting a patient from Server	43
11.5	Acq	uisition: general instructions	46
11.	5.1	Positioning the patient	46
11.	5.2	Description of the Acquisition screen	47
11.6	Bior	netry: acquisition procedure	50
11.7	Full	biometry acquisition (K-AL-CCT-ACD-LT)	52
11.8	Acq	uisition of axial length measurements (AL)	53
11.9	Acq	uisition of the anterior segment sections (CCT-ACD-LT)	54
11.10	K	eratometry acquisition (KER)	54
11.11	Р	upillometry	56
11.12	R	eport printing	58
11.	12.1	Available Printers	59
11.	12.2	Custom Reports	60
11.13	D	ata Exportation	61
11.14	IC	DL calculation	63
11.	14.1	Data	64
11.	14.2	Spherical IOL calculation	66
11.	14.3	Toric IOL calculation	68
11.	14.4	Post Refractive IOL calculation	71
11.	14.5	Barrett Calculator	73
11.	14.6	Olsen Calculator	83
11.15	D	ata saving	87
11.16	R	X/AL Trend Function	. 88
11.	16.1	RX/AL Trend Action Window	88

	11.16.2	New Exam Action	
	11.16.3	Data Review	
	11.16.4	RX/AL Trend Option	
12	MEASU	IREMENTS	
12	2.1 То	pographic map (KER)	
	12.1.1	Keratometry	100
	12.1.2	Keratorefractive indices	100
	12.1.3	Keratoconus	101
	12.1.4	Pupil	102
	12.1.5	Profile	103
12	2.2 Ze	rnike	104
12	2.3 Ax	ial Lenght (AL)	107
12	2.4 An	terior Segment Sections (ANT, CCT-ACD-LT)	108
12	2.5 Pu	pillometry (PUP)	109
	12.5.1	Graphs	110
12	2.6 W	nite To White (WTW)	114
13	SETTIN	GS	116
13	8.1 Ge	neral	116
13	3.2 M	easurements	117
	13.2.1	Мар	117
	13.2.2	Acquisition	
	13.2.3	Pupillometry	118
	13.2.4	Map Option	119
	13.2.5	Topography Map Color scale description	119
13	8.3 Su	rgeons	120
13	8.4 10		121
	13.4.1	General	121
	13.4.2	Preset	122
	13.4.3	IOL list	124
13	3.5 Co	nnectivity	132
	13.5.1	Network folder configuration	132
	13.5.2	XML Export	
	13.5.3	IMAGEnet i-base software	133
	13.5.4	IMAGEnet 6 Server software	
	13.5.5	Export to External Software settings	135

5

	13.	5.6	DICOM	136
13	3.6	Adm	in	138
	13.0	5.1	Report	139
	13.0	5.2	Remote Assistance	139
	13.0	6.3	Updating the integrated software	140
	13.0	5.4	Backup & Restore	145
	13.0	6.5	Shut down	149
13	3.7	Priva	acy & Security Settings	150
	13.	7.1	Password protected Login	150
	13.	7.2	Screen Saver / Auto Log OFF	151
	13.	7.3	Hide Patient Details in reports and exported exam packages	151
	13.	7.4	Anonymize filenames of exported reports and exam packages	152
	13.	7.5	Password protect exported PDF reports	152
	13.	7.6	Hide Patient Names and Disable Actions when in assistance mode	152
14	CHA	ANGIN	NG THE FUSES	154
15	TEC	HNIC	AL SPECIFICATIONS	155
16	PER	FORN	MANCE TESTS	158
17	DEC	CLARA	TION OF CONFORMITY	161
18	APF	PENDI	X A: INSTALLING AN EXTERNAL PRINTER	162
18	3.1	Gett	ing drivers and transferring them to ALADDIN	162
18	3.2	Disa	bling the Write Filter	162
18	3.3	Insta	alling a local printer (USB)	163
18	3.4	Insta	alling a network printer (LAN)	166
18	3.5	Re-E	nabling the Write Filter	167
19	APF	PENDI	X B: BIBLIOGRAPHY	168

1 INTENDED USE

ALADDIN HW3.0 is intended for biometric determination of the following ocular measurements: axial length, corneal radius, corneal cylinder axis, anterior chamber depth, central corneal thickness, crystalline lens thickness, white-to-white (WTW) and pupil diameter of the human eye. ALADDIN HW3.0 also measures corneal topography.

For patients who are candidates for intraocular lens (IOL) implantation, ALADDIN HW3.0 also aids in the calculation of the appropriate IOL power and type to be implanted.

ALADDIN HW3.0 is intended for use by physicians and eye-care professionals and may only be used under the supervision of a physician.

1.1 Intended users

ALADDIN HW3.0 is intended for use by physicians and eye-care professionals (opticians, ophthalmologists) and may only be used under the supervision of a physician.

For surgery and intraocular lens implantation, the device can only be used under medical supervision. For the other applications, the device must be used by qualified personnel.

1.2 Places of use

The intended places of use are: health care centers, doctors' surgeries, operating theatres.

1.3 Contraindications

Patient could have a dazzle effect, after the exam, dues to the device lights, but it disappears in few minutes.

1.4 Description of functionalities

The ALADDIN HW3.0 is a combined device for the biometric measurements of ocular structures. The measurements assist in the determination of the appropriate power and type of intraocular lens.

Below a summary of the device functions:

OPTICAL BIOMETRY - Measurement of the following ocular structures by means of the low-coherence optical interferometry method:

- Axial Length (AL): distance between the cornea and the inner limiting membrane.
- Anterior Chamber Depth (ACD): distance between the anterior surface of the crystalline (anterior capsule) and the outermost stratum of the cornea (epithelium), measured along the central axis where the latter is biggest.
- **Central Corneal Thickness (CCT)**: distance between the outer stratum of the cornea (epithelium) and the inner stratum (endothelium) on the central axis.
- Lens Thickness (LT): distance from the anterior surface of the crystalline lens (anterior capsule) and the posterior surface, measured along the central axis.

TOPOGRAPHY - Acquisition of the topographic map of the eye, through the reflection of 24 rings of the Placido disk at a distance of 80 millimeters from the eye. The following measurement are performed:

- **Keratometry (KER)**: measurement of the curvature of the anterior surface of the cornea. Includes all the basic functions of corneal topography (ketorefractive parameters and Zernike analysis).
- White-to-White (WTW) or corneal diameter: horizontal distance between the borders of the corneal limbus.

PUPILLOMETRY - It is performed with image analysis of the sequence of acquisition. The acquisition is performed with LEDs of different wavelengths. The device uses infrared LEDs to dilate the pupil and white LEDs to reproduce photopic light conditions and to contract the pupil. There are three modes:

- DYNAMIC Pupillometry acquisition in dynamic controlled light condition: mesopic → photopic → mesopic.
- PHOTOPIC Pupillometry acquisition in static photopic controlled light conditions.
- MESOPIC Pupillometry acquisition in static mesopic controlled light conditions.

IOL CALCULATION - Intraocular lens power calculation and suggestion made using the following formulae: Holladay 1, Haigis, Hoffer Q, SRK / T, SRK II, Camellin-Calossi, Shammas No history, Barrett Universal II. There are three modes:

- **IOL calculation**:Calculating the total spherical power for an intraocular lens to be implanted.
- **Toric IOL calculation**: Calculating the spherical equivalent power, cylinder power and placement axis for a toric intraocular lens to be implanted.
- **Post refractive IOL calculation**: Calculating intraocular lenses for patients who have undergone refractive surgery to correct myopia or hypermetropia.

The ALADDIN HW3.0 has an onboard PC with a dedicated software providing all the functionalities described above.

Please refer to the literature references listed in **19. Appendix: bibliography.**

1.5 Essential Performance

- Keratometry (KER) Measurement
- Axial Length (AL) Measurement
- Anterior Chamber Depth (ACD), Central Corneal Thickness (CCT), Lens Thickness (LT) Measurement

1.6 Interaction with the patient

Patients do not control the device. They are positioned with their chin on the chinrest and their forehead on the forehead-rest and they are asked to stay perfectly still and look at the fixation point with one eye. The device is entirely controlled by specialized personnel.

The parts in contact with the patients are the followings:

- ABS chinrest
- Teflon forehead rest

2 PRECAUTIONS

This electronic instrument is a precision tool and it is intended to be used in professional healthcare environment including hospitals, physician's offices, surgical centers and limited care facilities, where equipment and systems are administered by healthcare professionals. Make sure to use it and keep it in a suitable place, at a normal temperature, humidity and atmospheric pressure out of direct sunlight.

- To ensure proper functioning, install the instrument in a vibration-free location.
- Connect all cables correctly before use.
- Use the recommended network voltage.
- When the instrument is not in use, turn off the power supply and protect it from the sun and from dust.
- To obtain accurate and reliable measurements, keep the measuring cone clean and free of dust.

This product conforms to the EMC standard (IEC 60601-1-2 4th Edition).

- ELECTRICAL MEDICAL DEVICES require special EMC precautions and must be installed and activated in accordance with the EMC instructions provided in the accompanying documentation.
- Use of accessories and cables other than those supplied with the instrument, except cables sold by the equipment manufacturer as spare parts, may lead to an increase in emissions and reduce the device's or system's immunity.
- The eventual cables connected to USB and LAN ports must be less than 3 meters length.
- The device should not be used adjacent to or stacked with other equipment; if adjacent or stacked use is inevitable, the equipment should be observed to verify normal operation in the configuration in which it will be used.
- Portable RF communications equipment should be used no closer than 30 cm to any part of the equipment, including specified cables. Otherwise, degradation of the performance of this equipment could result.
- Failure to follow warnings related to electromagnetic compatibility, can compromise the essential performance or basic safety of the device, affect the proper operation of the camera and software, false alarms, interruption of operations or wrong measurement.
- In presence of high level EMC disturbance (i.e. Electrostatic Discharge, Electrical fast transient/Burst) the device could stop working. In case, don't use on patient until the source of disturbance is removed.
- Care is exercised in design and manufacturing to minimize damage to devices under normal use. However, electronic devices are susceptible to many environmental stresses. The device may be affected by electrostatic discharge (ESD). In an environment likely to cause ESD, such as air conditioning, humidification, non-conductive floor coverings, synthetic clothing, discharge any charge collected on your body before touching the device.

The FDA labelling for some IOLs contain sizing based upon white-to-white measurements derived from studies in which this measurement is done with callipers. It is unknown whether the white-to-white measurement from this device yields results systematically biased compared to those from calliper measurements. Thus, sizing based upon white-to-white measurements from this device may not be consistent with those based upon measurements with callipers.

2.1 Electromagnetic compatibility

2.1.1 EM Emission

The ALADDIN HW3.0 device fulfills requirements of IEC 60601-1-2 4th Edition:

Emission limit	Standard	Compliance
Conducted and radiated RF	CISPR 11	Class B
Harmonic distortion	IEC 61000-3:2	Class A
Voltage fluctuations/ flicker emissions	IEC 61000-3:3	Compliant

2.1.2 EM Immunity

The ALADDIN HW3.0 device fulfills requirements of IEC 60601-1-2 4th Edition:

Immunity test	Standard		Test level										
Electrostatic discharge (ESD)	IEC 61000-4-2	\pm 8kV contact \pm 15kV air											
Electrical fast transient/Burst	Electrical fast transient/Burst IEC 61000-4-4		±2kV 100kHz repetition frequency										
Surges	Surges IEC 61000-4-5		n mode tial mode										
Rated Power frequency magnetic field	IEC 61000-4-8	30 A/m											
	IEC 61000-4-6	Level	Frequency	Modulation									
		3V	150kHz÷80MHz	1kHz 80% AM									
Conducted disturbances induced		6V	6,765MHz÷6,795MHz	1kHz 80% AM									
by RF fields		6V	13,553MHz÷13,567MHz	1kHz 80% AM									
											6V	26,957MHz÷27,283MHz	1kHz 80% AM
		6V	40,66 MHz ÷ 40,70 MHz	1kHz 80% AM									
Radiated RF EM fields	IEC 61000-4-3	Field (V/m)	Frequency	Modulation									
,	120 01000 4 3	3	80MHz÷2700MHz	1kHz 80% AM									

27	380MHz÷390MHz	18Hz 50% PM
28	430MHz÷470MHz	18Hz 50% PM
9	704MHz÷787MHz	217Hz 50% PM
28	800MHz÷960MHz	18Hz 50% PM
28	1700MHz÷1990MHz	217Hz 50% PM
28	2400MHz÷2570MHz	217Hz 50% PM
9	5100MHz÷5800MHz	217Hz 50% PM

Recommended separation distances between portable and mobile RF communication equipment and the device

The ALADDIN HW3.0 device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of device can help prevent electromagnetic interference by keeping a minimum distance between portable and mobile RF communication equipment (transmitters) and the device as recommended below, according to the maximum output power of the communication equipment.

Rated maximum output	Separation distance according to transmitter frequency (m)			
of transmitter (W)	150kHz to 80MHz d = 1.2 ·√P	80MHz to 800MHz d = 1.2 .√P	800MHz to 2GHz d = 2.3 ·√P	
0.01	0.12	0.12	0.23	
0.1	0.38	0.38	0.73	
1	1.2	1.2	2.3	
10	3.8	3.8	7.3	
100	12	12	23	

For transmitters rated at a maximum output power not listed above the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note:

(1) At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies

(2) These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

3 SYMBOLS

Symbols	IEC publications	Description
Ŕ	IEC 60417-5840	CLASS I DEVICE ACCORDING TO EN 60601-1 APPLIED PART TYPE B
CE		PRODUCT COMPLIANT WITH DIRECTIVE 93/42/EC
Туре А	EN ISO 19980	CORNEAL TOPOGRAPHY ACCORDING TO ISO 19980:2005
~	IEC 60417-5032	ALTERNATE CURRENT
REF	EN ISO 15223-1	REFERENCE OR MODEL NUMBER
8	ISO 7010-M002	FOLLOW THE INSTRUCTIONS FOR USE
	ISO 7010-W001	GENERAL WARNING
	ISO 7010-W001	CAUTION (GENERAL WARNING): TO AVOID INJURY CAUSED BY ELECTRIC SHOCK, DO NOT OPEN THE COVER. ASK YOUR DEALER FOR SERVICE.
	ISO 7010-W001	GENERAL WARNING: BE CAREFUL NOT TO HIT THE PATIENT'S EYES OR NOSE WITH THE INSTRUMENT DURING OPERATION. THE PATIENT MAY BE INJURED.
A AA	EN ISO 15223-1	MANUFACTURER
X	EN ISO 15223-1	TEMPERATURE LIMITATION Indicate the temperature limits to which the medical device can be safely exposed.
<i>%</i>	EN ISO 15223-1	HUMIDITY LIMITATION Indicate the range of humidity to which the medical device can be safely exposed.
\$•\$	EN ISO 15223-1	ATHMOSPHERIC PRESSURE LIMITATION Indicate the range of atmospheric pressure to which the medical device can be safely exposed.

Ť	EN ISO 15223-1	KEEP DRY Indicates a medical device that needs to be protected from moisture.
		HANDLE WITH CARE
ISO 780		THIS WAY UP Indicates correct upright position of the transport package.
С С	IEC 60417-5009	STAND-BY SWITCH
	 This symbol is solely applicable for EC member countries. To avoid potential negative consequences for the environment and possibly human health, this instrument should be disposed of (i) for EU member countries – in accordance with WEEE directive 2012/19/UE (Directive on Waste Electrical and Electronic Equipment) or (ii) for all other countries – in accordance with local disposal and recycling laws. 	

3.1 Sample labelling on the device



4 SAFETY INSTRUCTIONS

4.1 General

- ALADDIN HW3.0 should be used only for its intended purposes as detailed in this manual.
- It must be installed by qualified personnel.
- The device must be used in the environmental conditions as specified in this document.
- The least favorable environment is defined as the maximum values of temperature for the unit to be operating in, while the unit is consuming the maximum current. The environmental value is stated as +40°C. The maximum current absorption occurs during full biometry acquisition.
- The maximum temperature of applied parts (chinrest and headrest) can exceed 41°C when the device is
 used at environmental temperature close to 40°C. The device temperature doesn't exceed 48°C anyway.
 Considering the examination duration, the patient condition and the parts that are in contact with the
 patient, there aren't any known contraindications about to the contact with the device.
- If the device has just been delivered or has been subjected to thermal shock, wait at least one hour before making measurements on patients.
- Keep this manual at hand close to the device at all times.
- The physician or device user must inform the patient of the pertinent safety instructions and ensure that they are adhered to.
- Connect the device to the supply mains using one of the cables supplied with the device
- A Position the unit so that it is not difficult to disconnect the plug for connection to the supply main.
- Perform all the control functions (detailed in the relative section in this document) before carrying out measurements on patients. In addition, if software interface shows an "Initializing error" warning, don't go on with measurements. Also "Low repeatability of measure" warning originates a wrong IOL calculation.
- Only personnel with the appropriate training and experience may use the device and interpret the results.
- Use of the device requires training and professional skills. The scholastic and cultural preparation of the user and the User Manual reading are enough as training.
- Turn off the device if it is not going to be used for a long period of time.
- A If external forces act on the device (e.g. if it is knocked or dropped), it must be thoroughly checked before proceeding to examine patients. To do this, refer to the relative section in this manual. If necessary, send the device in for repair.
- Use only original ALADDIN HW3.0 accessories and spare parts specific for this device.
- Remove all the covering (dust sheet) from the device before turning it on.
- Do not use the device close to highly inflammable materials or in areas with an explosion hazard.
- Unauthorized installation of software in the device is forbidden.
- After the examination, the patient may be slightly dazed. It is recommended to advise the patient to wait a few minutes before driving or performing actions that require perfect vision.

• A When operating the chinrest up/down switch, be careful not to pinch the patient's hand. The patient may be injured.

4.2 Electrical safety

- To avoid risk of electric shock, this device must be connected to supply mains with protective earth.
- ALADDIN HW3.0 has an on-board power supply unit installed. For connection to the mains, use only the manufacturer-approved cables provided with the device.
- Before performing maintenance on the device, turn it off and disconnect the power cable.
- Maintenance activities must be done in absence of the patient.
- Fuse change can be done by the user by following the instructions and safety precautions described in this Manual.
- Do not touch the LAB/USB ports contacts and the patient at the same time.

4.3 LED emission safety

CAUTION - The light emitted from this device is potentially hazardous. The longer the duration of exposure, the higher the risk of ocular damage. Exposure to the device light when using it at maximum intensity will infringe the safety indication after a 60-minutes use.

ALADDIN HW3.0 has a series of LEDs of various types and powers installed. All the characteristics are detailed in the Technical Specifications section in this manual.

The LED groups comply with the emission limits for the Group 2 instruments according to the standard ISO 15004-2.

4.4 Installation with external devices or IT Network

ALADDIN HW3.0 complies with the CE marking requirements.

- A Before connecting an external device, such as a computer, printer, monitor, keyboard, mouse or other devices, make sure that they comply with the EN 60950-1 standard and have the CE marking.
- An external keyboard or another input device compatible with "*keyboard wedge interface*" (PS/2) such as barcode or card readers can be connected to the device to input text.
- When ALADDIN HW3.0 is installed in rooms for medical use, the PC and the connected printer must be powered by means of an IEC 60601-1 compliant insulating transformer.
- If ALADDIN HW3.0 is installed in rooms for medical use without a computer, it is not necessary to use an insulating transformer.
- Do not use mobile phones or other devices not compliant with the requirements of class B EMC close to ALADDIN HW3.0.
- Every external device that has to be connected to ALADDIN HW3.0 must have a connection cable (USB or LAN) with a maximum length of 3 m.
- After connection of external devices to the USB or LAN, the end installator must check that the system maintains basic safety and essential performance of the product in compliance with IEC 60601-1.

- The purpose of ALADDIN HW3.0 connection to an IT network is report printing and remote technical assistance.
- The ALADDIN HW3.0 USB port must be connected to printer with USB or LAN interface. Ask Topcon technical assistance for printer driver installation.
- The ALADDIN HW3.0 can be connected to a Local Area Network (LAN) through the LAN connector. The network must have Ethernet protocol (IEEE 802.3). Ask Topcon technical assistance and the system administrator for ALADDIN HW3.0 and network settings.
- The purpose of ALADDIN HW3.0 connection is saving PDF report on an external network folder or technical service intervention on the machine.
- Connection of ALADDIN HW3.0 to a computer network that includes other equipment could result in previously unidentified RISKS; identify, analyse, and control such RISKS (refer to IEC 60601-1:2005).
- Subsequent changes to a computer network could introduce new RISKS and require new analysis.
- Changes to the computer network include:
 - Changes in computer or data network configuration
 - Connection of additional items to computer network
 - Disconnecting items from computer network
 - Update of equipment connected to computer network
 - Upgrade of equipment connected to computer network
- The term computer network used here corresponds to the term network/data coupling in IEC 60601-1:2005.
- Do not change the Ethernet settings "Local Area Connection 2" adapter.
 If specific configuration for LAN network is needed in the settings can be modified accordingly on "Local Area Connection" Ethernet adapter.

4.5 Transport and packaging

- The device must be transported and stored in its original packaging.
- For the storage and transport conditions, refer to the relative section in this document.
- Carefully keep the original packaging in order to use it if you need to transport the device.
- To move the device for short distances (without packaging) and to insert it in and remove it from the original packaging, grip the device with both hands, one on the front headrest arch and the other in the recess on the rear of the device (where the locking system is).
- Completely unscrew the two locking screws and the semi-lock (Figure 2) before use.
- Lower the instrument to its minimum height using the joystick, then lock ALADDIN HW3.0 using the semi-lock and the two locking screws for transportation (Figure 2).

4.6 Cleaning

- Regularly clean dust off the device using a soft cloth. For more persistent superficial dirt, use a soft cloth dampened with water or alcohol at maximum 70%.
- A Be careful not to get the device wet and clean it only as indicated to prevent damaging it. Never use solvents or other abrasive agents.
- The device comes with a dust cover to be used to protect it. Cover ALADDIN HW3.0 if it is not going to be used for a long period of time.

• Before turning on the device, remove the cover. Never put the cover on when the device is on.

4.7 Checking the measurements

- The calibration must be checked when the device has been transported from one place to another and when it has suffered an impact or thermal shocks.
- Check the measurements every day when turning on the device using the instrument provided.
- The user of the device must check that the measurements provided by the device are plausible.
- It is advisable to visually check all the light sources before examining patients to make sure that they come on properly.
- If the device frequently emits error signals, turn it off and contact technical support to have the device checked.
- In patients with blue eyes, acquisition of pupillometry in mesopic lighting conditions can be difficult to accomplish. In this case, we suggest acquiring the mesopic data through dynamic pupillometry.
- Contact lenses must not be worn by the patient during data acquisition.

4.8 **Privacy & Cybersecurity**

- When performing the installation of a new unit the user MUST set his own credentials to prevent unauthorized physical access to the device. To set up the access login to the instrument operation and on-board data refer to section 13.7.
- Make sure the USB devices you intend to connect to the instrument are secured against malware/viruses.
- A Patient data on USB devices can become corrupted when inserting into computers for backup or transfer.
- The use of antivirus software on computers is recommended and it is responsibility of the user.
- To protect data exported to USB from unauthorized access, use dedicated USB data for storage.
- Installation of any unapproved software, including drivers, could degrade the performances of the instrument and may void the instrument warranty.

4.8.1 Privacy

It is recommended that the operator understand the characteristics of this device in relation to Data Protection.

4.8.1.1 Definitions and application within this device

Personal Data: means any information relating to an identified or identifiable natural person ('data subject'); an identifiable natural person is one who can be identified, directly or indirectly, in particular by reference to an identifier such as a name, an identification number, location data, an online identifier or to one or more factors specific to the physical, physiological, genetic, mental, economic, cultural or social identity of that natural person.

Personal Data used with the device is:

- Subject identification data: subject name, subject date of birth, subject gender, subject ID (according to the practise ID application policy). Refer to section 11.4 for further details on the management of subject identification data.
- Subject physical/physiological data: subject eyes biometric measurements, as acquired and stored in each examination belonging to the subject. Refer to section 1.4 for further details on the types of biometric measurements performed by the device.

Controller: means the natural or legal person, public authority, agency or other body which, alone or jointly with others, determines the purposes and means of the processing of personal data.

The controller is the representative of the practise, clinic or hospital where the device is installed.

Processor: means a natural or legal person, public authority, agency or other body which processes personal data on behalf of the controller.

The processor is/are the person/s designed by the practise, clinic or hospital to operate the device and use the results.

4.8.1.2 Device characteristics to support Data Protection

Data Protection

The device provides by default encryption of data stored, to avoid unauthorized access to Personal Data.

The device provides by default technical measures to provide integrity, availability and resilience of subject data by means of encryption, data integrity checksum and redundancy.

The data can be accessed in relation to subject identification only within the application itself which is protected by login that can be configured by the controller to grant confidentiality.

The subject data it is used within the device on-board application only for storage and identification purpose. The controller can apply an explicit pseudonymisation rule to handle subject identification by using only ID subject identification and not using subject name, DoB, gender in the device.

Right of Access by the data subject

The controller and processor have access to data related to the subject which can be stored on the device and can make the data available to the subject in the normal operating mode. The data can be also accessed in the form of printable or electronic file (PDF or XML file).

Right to Rectification or Erasure

The controller and processor can, upon request from the data subject, rectify the subject identification data.

The controller and processor can, upon request from the data subject, erase permanently the subject identification data and subject physical/physiological data. After erasure no trace of the subject data is present on the device.

Refer to section 11.4.2.2 for details on how to rectify or erase patient data.

Right to restriction of processing

The device application doesn't automatically process or transmit the subject data.

Right to data portability

The controller and processor can provide to the data subject the data in the form of exported or printed reports or XML data file.

Right to object

No automated processing of data subject is performed in this device. The processing is performed on decision of the processor.

Automated individual decision-making, including profiling

No automated decision-making or profiling is performed in this device.

4.8.1.3 Privacy and Data Protection during assistance/support services

The device provides functionalities to assure Data Protection also during assistance/support services, to avoid disclosure of Personal Data to unauthorized personnel and assure data integrity and availability.

Refer to section 13.7 for specific options.

4.8.2 Privacy & Security options

When the on board application is booted, the following prompt is shown.



By pressing the OK button you get access to the Privacy & Security settings of the application, where you can define a login password and customize the level and type of data protection.

By pressing Cancel you choose to not proceed with the customization of the data protection options, the message will be prompted at the next boot.

The Privacy and Security options can be access at any time from the settings. For the details on each option refer to section 13.7.

PRIVACY AND SECURITY			
In this section you can configure your data protection options within this device. The main option is the definition of the user password to access the device. For further details on the data protection options please refer to sections 4.8 and 13.7 of the user manual.			
Use Password Login to operate the device Configure Password A PASSWORD NOT CONFIGURED			
Auto LOG OFF / Screen Saver Timeout [min]:			
Hide Patient Details in reports and exported exam packages			
Anonymize filenames of exported reports and exam packages			
Password protect exported PDF reports (use login password)			
Hide Patient Names and Disable Actions when in assistance mode			
Close			

- 5 PRODUCT WARRANTY AND RELIABILITY
- A The product warranty is valid only if all the instructions detailed in this document are followed.
- The product warranty is forfeited in the event of loss or damage due to improper or incorrect use of the device.
- The product warranty is valid only if it is equipped with its original accessories.
- A If the device is opened by unauthorized personnel, the manufacturer is relieved of all responsibility and the warranty shall become null and void.
- **N.B.**: Modifications or repairs to the product, especially where they require opening the device, may only be carried out by technical personnel authorized by the manufacturer.

6 LEGAL PROVISIONS

93/42/EEC – 2007/47/EC:	\rightarrow	Class IIA medical device
IEC 60601-1:	\rightarrow	Class I type B continuous operation
IEC 60601-1-2:	\rightarrow	EMC
ISO 15004-2:	\rightarrow	Group 2
UNI EN ISO 19980	\rightarrow	Туре А

ETL CLASSIFIED	ETL Standards:		
	Medical Electrical Equipment – Part 1: General Requirements For		
	Basic Safety and Essential Performance [AAMI ES60601-		
	1:2005+A1]		
C US	Medical Electrical Equipment – Part 1: General Requirements For		
	Basic Safety and Essential Performance [CSA C22.2#60601-1:2014		
Intertek	Ed.3]		

7 MAIN COMPONENTS

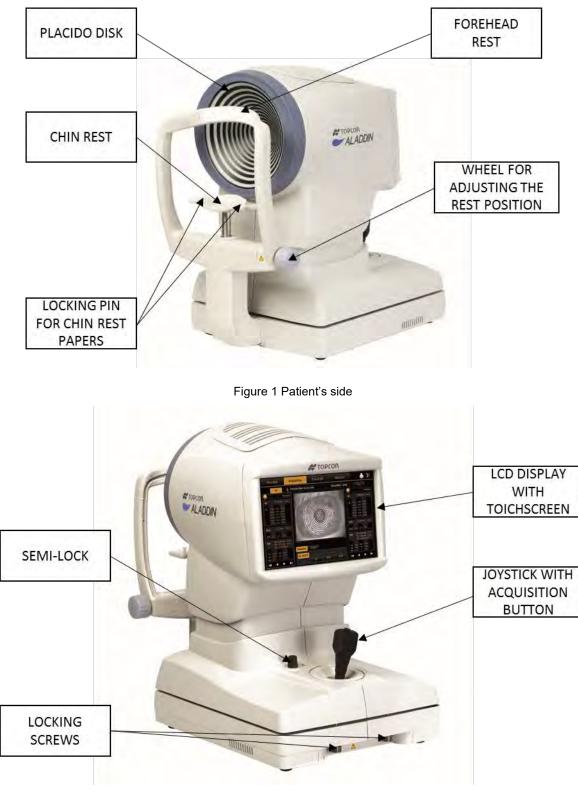


Figure 2 User's side

NB: The parts in contact with the patient (applied parts) are the forehead rest in Teflon and the chin rest in acrylonitrile butadiene styrene resin (ABS)

8 INSTALLATION /UNINSTALLATION OF THE SYSTEM

ALADDIN HW3.0 is packed for shipping in a double cardboard box on a dedicated pallet with specially shaped cardboard parts inside to guarantee instrument safety during shipment.

Keep the original packaging for future use. The system must always be moved/shipped in its original packaging, which is specifically designed for damage protection.

8.1 Installing the system

Before installing the system, read the "Safety Instructions" in this manual.



Figure 3

Figure 3 shows the complete packaging of the instrument.

Cut the extensible film and the packing straps. Open the external box, and remove the wood panel as shown in Figure 4.



Figure 4

Remove the manual and the accessories from the dedicated spaces between the two pieces of cardboard (see Figure 5).



The accessories are:

- "Topcon" box:
 - o calibration checking device
 - o chin rest paper
 - o chin rest pins
 - o touchscreen pen
 - o silicon cloth
- Power cable
- "Topcon" ALADDIN HW3.0 dust cover
- ALADDIN HW3.0 user manual

Open the internal box and remove the specially shaped cardboard that holds the instrument. The instrument can now be taken out of the package. The steps are illustrated in Figure 6.



Rev. 19 EN of 10/07/2019

ALADDIN - User manual



Figure 6

Be careful when taking ALADDIN HW 3.0 out of the box gripping it by the chin rest arch and the base beside the joystick.

Remove the Nylon cover.

Place the instrument on a flat surface.

Completely unscrew the two locking screws and the semi-lock (Figure 2).

Connect the power cable provided. The instrument is now ready for use.

8.2 Uninstalling the system

Take the original packaging.



Set the instrument to the minimum height using the joystick. Lock the device using the instrument semi-lock and the two locking screws for transportation (Figure 2).





Figure 7

Place the Nylon cover over the instrument and insert it in the box, as shown in Figure 7. Follow the sequence of steps shown in Figure 8.







Put the accessories in the dedicated spaces. Position the wood panel with the shock absorbers in the lower part. Close the external box with strong packing tape or use extensible film and packing straps.

9 STANDARD ACCESSORIES AND EQUIPMENTS

Calibration checking device The calibration checking device shows the serial number of the instrument with which it is associated. To properly check calibration, the calibrator provided with the instrument must always be used.	
Power cable	99
Manual	ALADDIN User manual User manual
Protective cover	TOPCON
 Accessories box containing: Touchscreen pen Silicon cloth Chin rest paper Chin rest pins 	

10 SETTING UP THE INSTRUMENT

10.1 Connection modes

Before connect device to external devices or IT Network, please read paragraph 4.4.

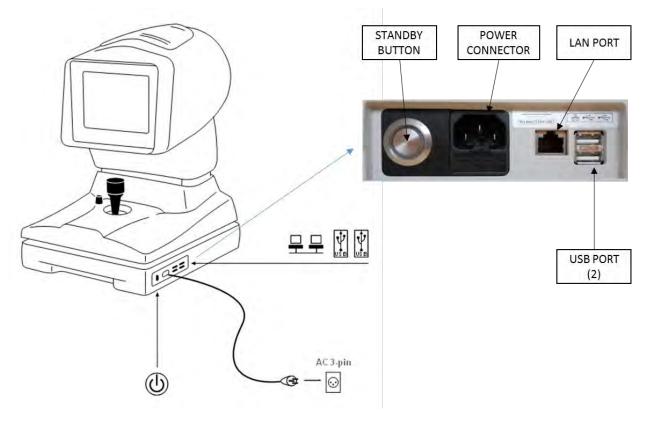


Figure 9

10.2 Power on procedure

Make sure the power cord provided is connected to the mains.

Press the stand-by button (see Figure 9) and wait for the system loading until the screen showed in Figure 16 appears.

11 OPERATING INSTRUCTIONS

ALADDIN HW3.0 is designed to work in stand-alone mode. For this reason, all the software functions are automatically loaded when the device is turned on, enabling the user to control the device and guiding him or her through the various phases:

- Entry of patient data
- Acquisition of the various possible modes
- Display of measurements
- Selection of intraocular lenses

More information for each function and setting are provided in the following paragraphs.

To interact with the software, the LCD display with touchscreen is used. To activate the button or the desired function, simply touch the screen close to the command. The screen is highly sensitive. Minimum pressure is required, indeed advised.

11.1 User Login/Logout

The user MUST configure the login password to access the device operations and stored data, in order to assure confidentiality and integrity of personal data and avoid unauthorized access.

It is recommended to make a robust password choice. It is recommended to keep record of the password in a safe place.

It is not necessary to disclose the user password to technical support staff. Refer to section 11.1.3 for details.

Refer to section 13.7 on how to configure and manager the login access. If the user login password is configured and enabled, the login screen will be prompted at the start-up of the Aladdin on-board application. The login screen is prompted also to allow access to the privacy and security settings (section 13.7).



Type the user password as configured to access the operations and data of the device.

If the password is forgotten or access problems are encountered please contact assistance. It is possible to get remote support by pressing on the button "Remote Assistance" if the device is connected to internet. Refer to section 13.6.2 for further details on TeamViewer remote support.

11.1.1 User Logout


It is possible, at any time, to lock the device operations by accessing the main views of the applications and

pressing on the lock button. The login screen will be prompted again.



11.1.2 Auto LogOFF/ Screen Saver

The device locks automatically if no user actions occurs for a configurable amount of time. The Aladdin screen saver is then shown.

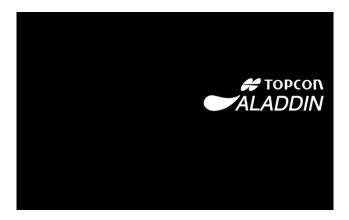

To unlock the touch the device screen. If the login password is configured and enabled the user login password is requested.

To configure the screen saver timeout or enable/disable it refer to section 13.7.2.

11.1.3 Technician Support staff Login

Technical staff can log into the device and on-board application without knowing the user password. The authorized technical staff knows the technician password which allows to do service on the device having data protection features applied according to privacy options of section 13.7.

To let the technical staff log in properly press on the lock button before.

11.2 Overview

Every working environment has the same screen layout. In Figure 10 a sample of the acquisition screen is shown.

The active w	orking area is enlightened in orange.	Main menu	Print	Settings
Ma	in Acquisition	IOL Calculation	Measurements	
	OD S IPPCON DEMC		05/11/2014 - 12:0	os
AL V AL ACD ACD LT LT KER	Phakic		AL V AL ACD ACD ACD IT IT IT	Phakic
K2 Axis CCT ↓	K-AL-AN		CCT-ACD-LT	
Righ	t eye data		Left	eye data
	The measurement select enlightened in oran		Live images from camera	



11.3 Checking the calibration


- The calibration must be checked:
 - every day before starting patient examinations;
 - when the device has been transported from one place to another;
 - when it has suffered an impact or thermal shocks.

To check the calibration, turn on the device, and when asked to check the calibration, press START. The test patient is automatically created.

New	Searc	h	Sea	rch Server	Acq	uisition	03/12/201 13:46	8 🤷 🧎
ATIENT DETAILS								
Last Name				ID				
								Ök
Name		WARNING	5: Check ca	alibration before s	tarting the m	easurements	1	Clear
Data of Rith	(dd/mm/yyyy)	on patient	ts					
Date of Birth	ICIC/mm/www							
	(44,1,111,1,1,1,1,1,1,1,1,1,1,1,1,1,1,1,			Start C	ancel			
				Start C	ancel			
		1	2		ancel 5 6	78	9	0 ESC
TO	Φροοη			3 4		78 YU	9 1 0	
	ορςοη	Q	2 W	3 4 E R	5 6 T	ΥU	ΙΟ	P
₩ TO ALA	ορςοη			3 4	5 6 T			

Figure 11

The next screen explains how to check calibration:

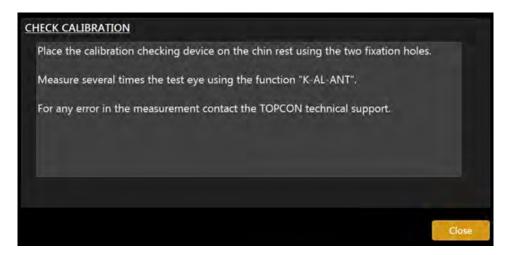



Set the calibration tool supplied with the device in the special holes in the chin rest and press until the tool is blocked on the device. Check that the calibration tool is perfectly aligned with the device. If the calibration tool is positioned correctly, all the rings of the Placido disk should be seen reflected in the center on the surface of the hemisphere (Figure 13).

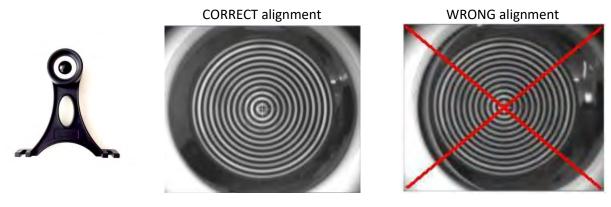


Figure 13

Now check and several times acquire the calibration checking device using the complete acquisition (K-AL-ANT). For details on this acquisition mode, see chapter 11.5.

If the calibration is ok, the "Valid" word will be display for all the measurements (Figure 14). If the measurements are incorrect, the words "Repeat" or "Not Valid" will be displayed besides the wrong measurement (Figure 15).

See table below for instructions.

MESSAGE	MEANING	WHAT TO DO
VALID	Good calibration	Click on "Main" to start a new examination. When asked, press Yes to save the "Calibration Check". NB: the "Valid" word must be displayed for ALL the measurements.
REPEAT	Acquisition failed: cannot check calibration	 Acquire again until three times following these suggestions: improve the environmental conditions (less light and no reflections on the sphere); clean the sphere of the calibration tool;
NOT VALID	Incorrect calibration	 make sure that the calibration tool is position correctly. If the calibration check is still not valid, do not take any patient measurement and contact Technical Support.

ALADDIN - User manual

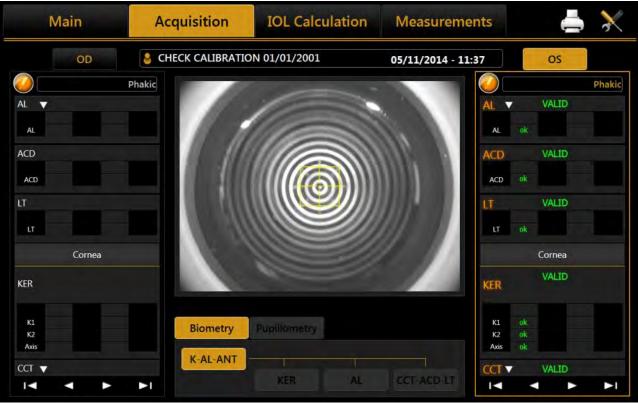


Figure 14

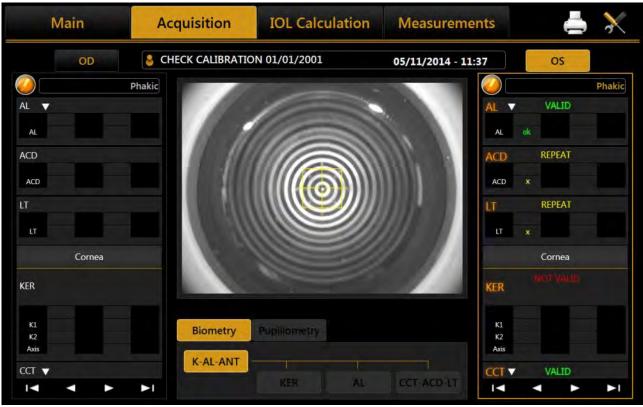


Figure 15

11.4 Patient entry/selection

When the instrument is turned on, the software displays the following screen. Before start examination, it is required to create a new patient or select a previously saved one from database.

New	Search	Sear	rch Ser	ver	Á	cquisi	tion	19/0	6/2016 2:	1:25 🔀
PATIENT DETAILS										
Last Name			ID				-			Ok
Name			Gender							lear
Date of Birth (do	l/mm/vvvv)			Vale emale						
	1	2	3 4	4 !	5 (5 7	7 8	3 9	0	ESC
🖨 ТОР		Q W	E	R	T	Y	U	I	0	P
ALAD	DIN p	A S	D		G		J	К	L	J.
	Z	x	c	v	B	N	м	i,		~

Figure 16

Figure 16 shows the section for creating a new patient, entering Last Name, Name and Birth Date as required fields (Gender and ID are optional). You can set from the settings environment to have only the ID as required field.

11.4.1 Creating a new patient

To create a new patient, select the **"New"** tab and enter the data using the on-screen keyboard. Once you have entered the new patient data, click on the **"Ok"** button or select the **"Acquisition"** tab to confirm the information and continue with the examination. If you want to empty all the fields click on the **"Clear"** button. Before going into the acquisition environment, additional information on the patient is required, in particular the presence and type of crystalline and the nature of the vitreous body (Figure 18).

An external keyboard or another input device compatible with "*keyboard wedge interface*" (USB Type A) such as barcode or card reader can be connected to the device to input text. The user must assure that the desired textbox is under focus before the input action.

Before connecting an external device, such as a computer, printer, monitor, keyboard, mouse or other devices, make sure that they comply with the EN 60950-1 standard and have the CE marking.

11.4.1.1 Entering special characters

A special character can be entered simply by touching and holding for a few seconds the corresponding letter as shown in Figure 17:



Figure 17

11.4.1.2 Selecting crystalline and vitreous body type

Once the patient identity record has been created, it is possible to select the type of cristallyne and vitreous body for each patient's eye, by pressing the "Acquisition" button (please see the following Figure 18)



Figure 18

For each eye, select the type of crystalline currently present:

- Phakic: the patient has a natural crystalline lens.
- Aphakic: the patient does not have any crystalline lens from birth or as a result of surgery.
- **Pseudophakic:** the patient has an intraocular lens substituting the crystalline. In this case, it is very important to also detail the type of material used by the surgeon:

- o **Unknown**
- o Silicon
- o **PMMA**
- Acrylate
- o Memory

The measured axial eye length depends on the measuring mode selected. Depending on the measuring mode selected, ALADDIN HW3.0 corrects the measurement with a constant defined as follows.

ALADDIN HW3.0 device takes into consideration two conditions of the eye that can alter the measurement of axial length:

- Vitreous body filled of silicone oil
- Implant of intra ocular lens

The difference of the measurement is caused by a different group refraction index considered in the formula. According to bibliographic data, the calculations have been performed to assess the amount of correction that must be applied to correct the measurement in these special cases.

The correction data have been compared with predicate device assumptions and a table of corrections has been elaborated as follows:

The correction values (in mm) of the natural vitreous body

Phakic	0
Aphakic	0.21
Pseudophakic Unknown material	0.11
Pseudophakic Silicone IOL	0.12
Pseudophakic PMMA IOL	0.11
Pseudophakic Acrylic IOL	0.1
Pseudophakic Memory IOL	0.11

For the vitreous body you can choose between:

- **Natural:** the vitreous body has never been operated or treated such as to alter its composition.
- Silicon Oil: the vitreous body has been filled, even only partly, with silicon oil.

The correction values (in mm) of the vitreous body filled by Silicon Oil

Phakic	-0.74
Aphakic	-0.86
Pseudophakic Unknown material	-0.75
Pseudophakic Silicone IOL	-0.74
Pseudophakic PMMA IOL	-0.75
Pseudophakic Acrylic IOL	-0.76
Pseudophakic Memory IOL	-0.75

All this information is required because, on the basis of the artificial materials and their optical properties present inside the eye, the instrument always corrects the measurements obtained to the most precise value possible.

Once this information has been entered, you can access the acquisition environment.

For more details on the acquisition environment see the dedicated section.

The vitreous body nature is expressed, if different from natural, in the acquisition view as well as in the output reports, as shown in the following figures. The lens nature is always reported.



11.4.2 Selecting or modifying a patient

On the input screen, click on the **"List"** tab to access all the patients included in the local database (see Figure 19).

New	Search	Search Serv	/er	Acquisition	14/06/20	14/06/2016 11:55 እ	
Last Name	PATIEN Last Na		<u>EXAM L</u>	<u>IST</u>			
🖲 Surname		SAMPLE	26/0	5/2016 14:21:00		New Exam	
	Name		28/07/2014 19:20:00				
SAMPLE TOPCON1	Deter	TOPCON1					
TOPCON DEMO	Date o	20/08/1984				Open	
TOP-CON DEMO	ID					Delete	
		p1t9f31MIM				Delete	
	Gender ©			•		Post Op	
	1	2 3 4	5	6 7 8	9	0 ESC	
	C	Q W E	R	TYU	D C	D P	
	А	S D	F G	H J	KL	L.	
Edit Pag	e1/1 Z	хсч	/ B	N M	Ú)	~	

Figure 19

On this screen you can select a previously created patient and the examinations associated with him/her.

The list can be viewed by **patient ID** or by **Last Name** (and name) selecting the corresponding radio button.

If you type into the **"Last Name"** field, a search is done in the local database for patients with the corresponding surname or whose surname contains the selected key, same for patient ID.

By pressing the button on the right, the patient list is ordered alphabetically (A to Z) or by last exam date (most recent first).

11.4.2.1 Open an examination or acquire data for the selected patient

In the left column, clicking on a patient in the **"Exam List"** frame displays the list of associated examinations. In this list, you can access examinations or delete them, using the "Open" or "Delete" buttons. After having selected a patient, another examination can be carried out by pressing the **"Acquisition"** tab or pressing on "New Exam" button.

11.4.2.2 Delete or edit the selected patient

From the list of patients, select the exam you want to delete and press the "Delete" button. The program will ask you to confirm the choice.

Press "Edit" to change the name, surname or date of birth. This takes you back to the initial "**New**" tab. From here, you can edit the information you need to change and press **"Ok"** or **"Cancel"** to confirm or cancel the changes.

ALADDIN - User manual

11.4.2.3 Insert the Post-Op (after surgery) refraction data

Through this function the user can update the data related to a single exam of the chosen patient. This means that if the patient has already undergone surgery, the new refractive status can be recorded as a main factor to personalize constants of the implanted IOL.

R	👌 D. TOPCO	N 01/01/1950	0	10	0/02/2015 - 17	7:55	L
POST OP DAT	<u>A</u>			PRE OP DA	<u>TA</u>		
Surgeon	SurgeonPP Spherical		▼ orical	AL (mm)	23.73	ACD (mm)	3.14
Manufacturer	Oculentis		•	K1 (mm)	8.00	K2 (mm)	8.28
Model	L-303		•	Wtw (mm)	11.69		
IOL (D)	21	Sphere	-0.75				
Target	-0.5	CYL	0				
		Axis	0				
					Reset	Cancel	Close

Figure 20

Opening the Post-Op section, the screen shown Figure 20 will be displayed.

In this section you can insert the Post-Op Data (IOL information plus actual refraction) in the meantime looking at the Pre Operative Data.

11.4.3 Selecting a patient from Server

Once enabled, Aladdin **IMAGEnet i-base's** integration from Aladdin's settings panel (refer to IMAGEnet i-base configuration), it's possible to select a new patient from the patient list retrieved from IMAGEnet i-base (Figure 21).

In the same way, ALADDIN HW3.0 can be activated to search patients from DICOM services (refer to DICOM configuration section):

- DICOM Patient Root Query: search patient's details on enable patient's archive server
- DICOM Modality Worklist: get the list of patients and tasks in the waiting room

The user can search for a patient either by **surname**, by **id** or by **date of birth** (i-base only). A list of patients will be created corresponding to the search criteria (Figure 23). Once selected a patient, the user can create a new examination in the standard mode by clicking on the Acquisition or OK button button.

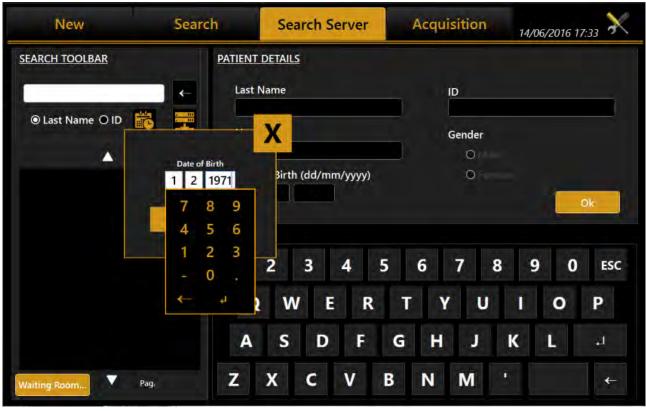


Figure 21

The user can search from IMAGEnet i-base and/or DICOM sources at the same time by enabling/disabling the corresponding options using the server selection button.

New	Search		Sea	rch S	erver		Acqui	sition	25,	/07/2017 0	9:20 🕅
SEARCH TOOLBAR	P.	ATIENT D		i			1	ID			
 Last Name ID SERVER i-base Dicom 		Name Date o	f Birth	(dd/mn	י/yyyy)			Gender O			Ok
		1	2	3	4	5	6	7	8	9 0	ESC
		Q	W	E	R	Т	Y	U		0	Ρ
		Α	S	D	F	G	H	J	К	L.	ц.
Waiting Room	Page 1 / 1	z 2	x	С	V	B	N	М	'*?		÷

Figure 22

New	Search		Sea	arch S	erver		Acqui	isition	25	/07/2017 0	9:19 እ
SEARCH TOOLBAR	<u>P</u> /	ATIENT	DETAIL	<u>s</u>							
*	~	Last N						ID			
● Last Name ● ID		Patie	nt					A-1000	2		
o Last Marrier O Ib		Name						Gender			
A		Test						ОM	ale		
PETCT		Date	of Birth	(dd/mn	n/yyyy)			O Fe	emale		
Patient Test											Ok
DEMO CA800											
DOE JOHN											
Patient Test		1	2	3	4	5	6	7	8	9 0	ESC
DOE MIDDLE JOHNI	NNN	-	2			<u> </u>	•		•	5	1000
		Q	N	/ E	R		Y	U	I	0	Р
		Α	S	D	F	G	(R)	J	К	L.	Ļ
Waiting Room	Page 1 / 1	z	x	с	v	B	N	м	'*?		\leftarrow

Figure 23

ALADDIN - User manual

11.4.3.1 Start an exam from the Waiting Room

If DICOM Modality Worklist service is configured, ALADDIN HW3.0 is able to search for pending patient's examinations in the waiting room. Pressing on the "**Waiting Room...**" button (Figure 22) shows a list of the pending worklists for the current day. The list can filtered by one or more of the other criteria:

- Patient Name
- Patient ID
- Examination date range
- Scheduled Station Name (default is "Aladdin")*
- Modality (default is "OT")*
- * = contact DICOM services administrator for details on these settings

			DICOM Workli	st Query				
Patient II	0			Start Date	-	Scheduled Station	n Name	
Patient N	lame			Stop Date	20	Aladdin Modality		
					20	от		
Patient ID SC-I1 SC-I1 pidP645	Patient Name SAMPLE^PATIENT3 SAMPLE^PATIENT2 SAMPLE^PATIENT	0	Patient Birthdate 7/16/1980 12:00:00 AM	Start Date 6/20/2016 10:50:42 AM 6/20/2016 10:50:42 AM 6/20/2016 10:50:42 AM	(Title Modality	Procedure ID 0000018705	Proce CSPIN
C .				_				3
Reset				Upd	ate Worklist	Start Work	Close	e

Figure 24

Each time the filtering criteria are changed, press "Update Worklist" to update the list of matching items.


Once the desired work is selected, press "**Start Work**" to start a new exam relative to the selected work.

11.5 Acquisition: general instructions


11.5.1 Positioning the patient

To obtain correct measurements it is necessary an adequate positioning of the patient relative to the device. A steady head position and the correct device-to-patient distance are obtained by resting the patient's head well against the chin rest and forehead band. The patient must look steadily at the fixation point in the center of the Placido disk.

A correct alignment with the patient's pupils can be visually checked by the operator referring to the two lines on the forehead supports (see the red arrows in Figure 25).

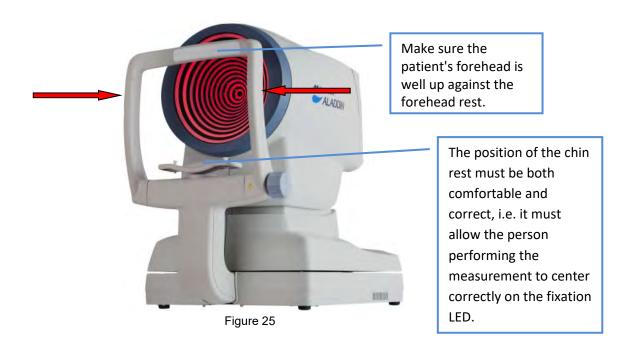

The working distance from the device to the eye is 80.0 mm. The software guides the operator in order to reach the ideal focus condition (see paragraph 11.6)

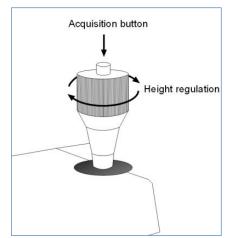

The joystick illustrated in Figure 26 is the only part the user has to physically control during acquisition. The button on the top marked "Acquisition button" starts the acquisition of the various measurements.

The thumb wheel marked "Height Regulation" allows you to adjust the instrument's height according to the patient's position. If the adjuster on the joystick is not enough to achieve the correct

position, move the knob on the chin rest for adjusting the height

Figure 26

11.5.2 Description of the Acquisition screen

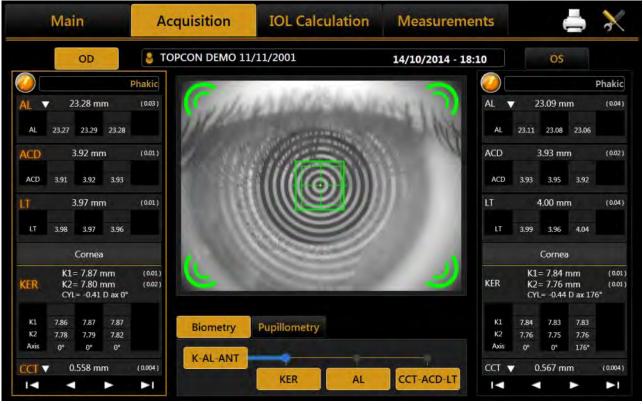


Figure 27

Figure 27 shows the acquisition screen from where all the operations to acquire the required measurements are performed.

The acquisition window has the following commands:

- **OD** and **OS**: indicate the eye being acquired (the one highlighted in yellow); they are normally selected automatically, depending on the position into which you move the instrument.
- *Biometry:* gives access to the biometric measurement section
- *Pupillometry:* gives access to the pupillometry section
- The buttons at the bottom of the data frame for each eye serve to scroll the measurements, as some of these are hidden if more than four acquisitions are made per eye.
- Modify the nature of the lens and vitreous body. The Axial Length will automatically be corrected depending on the new refractive index

When the acquisition is done, and the eye is selected: the correctly acquired measurements are displayed in white and incorrect measurements are displayed in red. See Figure 28.

11.5.2.1 Description of results

For each Biometry and Keratometry result a dedicated section is present. In each section the total result is shown together with the standard deviation between the single results (if more than one) and eventual warning or error signs (described in the following section).





Sometimes the measurement is taken in one of these conditions: bad focus, closed eyelid, tear film irregularity, high standard deviation in multiple measurements, movement, measurement not in range; in this case, a warning sign appears above the measurement.

ATTENTION: When the symbol A is shown above a measurement, it means that the software recognized an error during the acquisition, which could be: bad focus, closed eyelid, tear film irregularity, high standard deviation in multiple measurements, movement, and measurement not in range.



Figure 28

The possible problems in acquisition are found in the software with the following methods:

ALADDIN - User manual

Error cause	Identification
Closed eyelid	Missing ring reflection on the eye Placido Image, on the upper
	hemisphere of the cornea
Movement	Interlace pattern shown in acquired image
BrokenTear film	Missing ring reflection on the eye on Placido Image
Bad focus	Defocus of ring reflection on acquired image
High standard deviation on repetition	Big difference between repeated acquisitions
Measurement not in range	Output out of instrument range of measurement

If a warning sign is shown above a measurement, it is recommended to make further acquisitions until reliable data is obtained.

It is very important that the main types of measurement (KER-AL-ACD) taken are shown without a warning sign, otherwise it will not be possible to proceed to IOL Calculation with the current data (valid ACD is needed only if using *Haigis formula*).

Valid LT and CCT aren't needed to access IOL calculation as they are not used in any of the available formulas.

As shown in Figure 28, accessing IOL Calculation, an error window will be shown warning the user to reacquire the measurements or to manually enter a new set of data.

If the user proceeds with manual input (Figure 29), the software will pre-populate all the fields with acquired data, even those with error or warning signs. It is also possible to enter arbitrary data, if possible taken with other instruments.

						<u>0</u>
BIOMETRY			7	8	9	<u>KERATOMETRY</u>
SOURCE N			4	5	6	SOURCE MANUAL
AL		mm	4	2	•	K1 mm m
ACD LT		mm mm	1	2	3	K2 mm 90
сст		mm		0		Index 1.3375 V
WTW		mm	<i>←</i>		J	
Acquisition	Optical					
						Cancel Save

Figure 29

11.6 Biometry: acquisition procedure

There are two acquisition mode:

MODE	BUTTON	DESCRIPTION
FULL BIOMETRY ACQUISITION	K-AL-ANT	Perform the AL (Axial Length), CCT-ACD-LT sequence (Central Corneal Thickness, Anterior Chamber Depth and Lens Thickness) and KER (Keratometry) acquisitions in a full sequence measurement.
	KER	Perform KER (Keratometry) acquisition
SINGLE ACQUISITION	AL	Perform AL (Axial Length) acquisition
	CCT-ACD-LT	Perform CCT-ACD-LT sequence (Central Corneal Thickness, Anterior Chamber Depth and Lens Thickness)



WARNING: a valid LT result cannot be obtained without acquiring the patient's AL.

WARNING: a standalone CCT-ACD-LT measurement cannot be performed without a valid AL result

Backlighting of the Placido disk is automatically activated when you enter the acquisition environment. If the instrument is not used for a few minutes, the cone turns off; to turn it on again, just press the joystick button.

To acquire the image or measurements in general, whatever mode you are in, simply proceed as follows:

- 1. Align the live image in the center and focus, then press the joystick button to start the acquisition (you need to perform this step only the first time you acquire a patient eye).
- 2. Move the instrument forwards and backwards (following the indications of the red and blue arrows on the screen) to find the ideal focus. While you find the ideal focus achieve the central alignment by centering the two squared aims with vertical or horizontal movements.
- 3. When the green indicators are displayed and the two squares are centered (both green), press the joystick button again and the system will automatically capture the required image and/or measurements.
- 4. Don't move the joystick in the few seconds during the acquisition.

Focusing and centering guidance system is composed of two aspects:

- Centering
- Focusing distance

Centering ideal conditions are achieved by **centering the two squared aims** by means of horizontal and vertical movements, reaching the condition of following picture:

	 -

The two squares assume different colors depending on two aspects: focusing position and centering in tolerance.

Focusing ideal conditions are achieved by **following the 4 indicators** at the corners of the viewport, which explain the needed movement in the "forward/backward" direction.

	The red arrows indic	ate to move the instrument forwards towards the patient's eye.
	TH	ne centering condition is displayed in red color
	The blue arrows indi	cate to move the instrument backwards away from the patient.
	T	ne centering condition is displayed in blue color
	The green signs indic	ate that the ideal focus has been reached.
	Press the joystick but	tton to start the automatic acquisition procedure.
	In this situation the o	entering aims assume different colors, as follows:
6		ne centering condition is out of tolerance and a better center alignment required
		ne centering condition is displayed in green color, press the joystick utton to start acquisition

At the end of the acquisition the measurements are displayed in the corresponding side windows (see Figure 27).

During the acquisition procedure, pulsating dots appear above the "KER", "AL" and "CCT-ACD-LT" buttons to guide the user through the acquisition steps. They are explained in the following table:

AL	The system is waiting for the user click on acquisition button. Follow the guide for right centering and focusing, then click the joystick.
AL	The system is acquiring. Wait until it has finished.
AL	The system has finished the acquiring procedure.
AL	The system is acquiring or waiting for user input in a previous step of the acquisition sequence.

11.7 Full biometry acquisition (K-AL-CCT-ACD-LT)

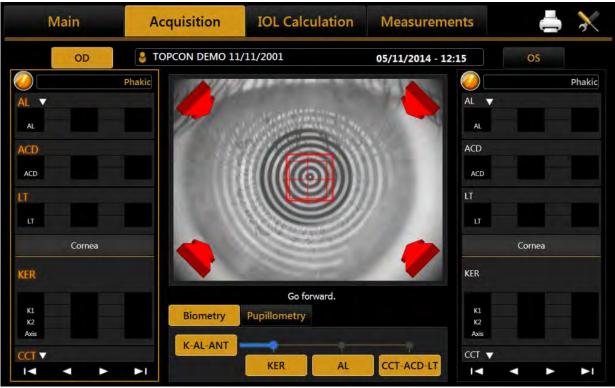


Figure 30

This is a special mode which successively performs all the measurements described in detail in the following paragraphs, specifically:

- Keratometry
- Axial Length
- Anterior Chamber Depth, Central Corneal Thickness and Lens Thickness
- Identification of Mesopic and Photopic pupil

Below the acquisition steps:

- 1. Press the joystick button a first time to start the procedure (Figure 30).
- 2. Find the ideal focus and achieve the central alignment as described in paragraph 11.6.
- 3. Press the joystick button to perform the Keratometry acquisition.
- 4. Next the system will automatically perform also the Axial Length acquisition (the progression of acquisition steps is shown in Figure 31, as described in paragraph 11.6).
- 5. Find the ideal focus and achieve the central alignment as described in paragraph 11.6.
- 6. Press the joystick button to perform the Anterior Chamber Depth, Central Corneal Thickness and Lens Thickness acquisition.
- 7. Wait few seconds for measurements calculation.
- 8. The measurements are displayed in the corresponding side windows (see Figure 27).
- 9. A new acquisition can be done.



Figure 31

11.8 Acquisition of axial length measurements (AL)

Interpretation of axial length measurements. As a rule, an interference signal is produced if the measuring light is reflected by the retinal pigmented epithelium of the eye. This signal is utilized for axial length measurements.

Note: Ultrasonic biometrical instruments measure the axial length as the distance between the cornea and the inner limiting membrane, because the sound waves are reflected at this membrane. To ensure that the measured values obtained with the ALADDIN HW3.0 are compatible with those obtained through acoustic axial length measurement, the system automatically adjusts for the distance difference between the inner limiting membrane and the pigmented epithelium. The displayed axial length values are thus directly comparable to those obtained by immersion ultrasound, and no re-calculation or correction factors are necessary. Deviations may nevertheless occur between the displayed axial lengths and ultrasonic readings (particularly in the applanation procedure).


By selecting this mode, the acquisition environment shown in Figure 32 appears.



Figure 32

The side columns show the measurements performed for the two eyes (OD = right, OS = left).

For each acquisition, six measurements of the axial length are performed. The information displayed is the same as in the "full sequence" acquisition. The six measurements are reported as a group result in the collapsed view, to expand the view and see each single measurement press the arrow next to AL title to expand/collapse the view.

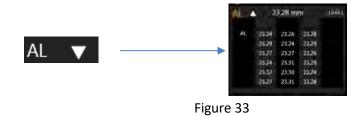

11.9 Acquisition of the anterior segment sections (CCT-ACD-LT)

The AL measurement must be performed before CCT-ACD-LT measurement in order to have a LT result without warning.



A standalone CCT-ACD-LT measurement cannot be perform without a valid AL result.

With this type of acquisition measurements about the sections of the anterior segment of the eye are taken using the interferometer system.

By selecting this mode, the acquisition environment appears as shown in Figure 34.

Press the joystick button and move the device according to the instructions of the automatic guidance system (red/blue arrows and centering aims) until ideal conditions have been reached (green icons and centered aims).

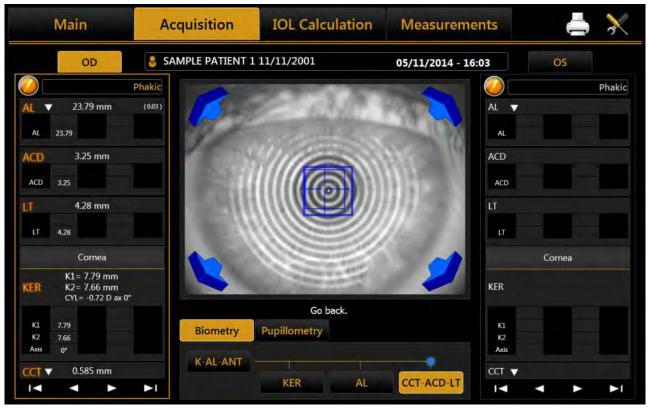


Figure 34

11.10 Keratometry acquisition (KER)

Keratometry is used to measure the corneal curvature. It is based on the reflection of the Placido disk on the eye at a controlled working distance for high measurement precision.

ALADDIN HW3.0 allows the user to acquire the corneal topography of the eye. The "Corneal Map" is obtained from the reflection of 24 rings of the Placido disk at a distance of 80 millimeters from the patient's eye. The position of the device, in relation to the patient's eye thus found, serves as a starting point for fine adjustments to be made in the respective measurement mode.

By selecting this mode, the acquisition environment shown in Figure 35 appears.

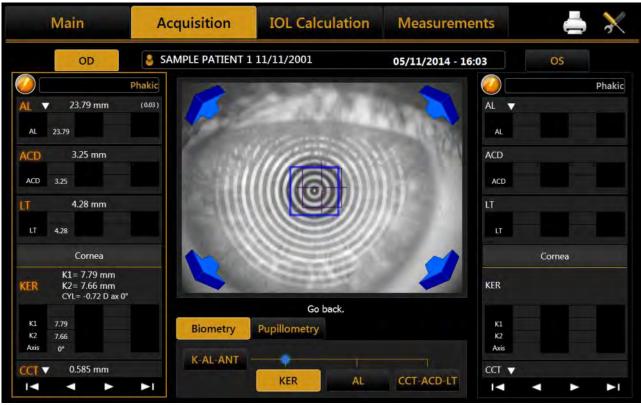


Figure 35

In this mode, the topographic map of the cornea is acquired.

Knowing the distance of the corneal apex, with a precision of microns, at the time of acquisition of the topographical image, the software applies to each of the 256 zero crossing, identified for each of the 24 RINGS, a correction factor given by the ratio between correct mean value and mean radius of the ring.

Concerning the calculation, the software performs the standard calculation of 6,144 zero crossing points, identified at the 24 RINGS along the 256 semi-meridian.

In order to increase the measurement precision, interferometry is used to evaluate the corneal distance.

The keratometry data is evaluated in the left column, which is referring to the right eye. This section has the same interactions as the "K-AL-CCT-ACD-LT" acquisition.

However, since it is not possible for the human eye to stay still, the images acquired manually in a range close to the optimum focus (which is the optimal operating range of the device) can be out-of-focus.

11.11 Pupillometry

Main	Acquisition	IOL Calculation	Measurements	📥 💥
OD	SAMPLE PATIENT 1	11/11/2001	05/11/2014 - 16:30	OS
AL ▼ 23.79 mm	Phakic (0.03)	Starte III	AL	Phakic
AL 23.79	-	A second	AL	
ACD 3.25	3	.	ACD	
LT 4.28		Alexand State	53) u	
Cornea K1= 7.79 mm K2= 7.66 mm CYL= -0.72 D ax 0°		a start	KER	Cornea
K1 7.79 K2 7.66 Axis 0°	Biometry	Pupillometry	K1 K2 Axis	
(CCT ▼ 0.585 mm) I ◄ ◄ ►	Dynamic ►I	Photopic	Mesopic CCT	



By selecting this mode, the acquisition environment shown in Figure 36 appears on the screen. In order to acquire the pupillometry, first of all you need to center the blue rectangle, which is overlaid in the image on the reflection of the four LEDs, as shown in Figure 37.

Press the joystick button to start the acquisition and press the button again to stop the acquisition. As already mentioned in the introductory paragraphs, three types of acquisition can be performed:

- Dynamic pupillometry
- Photopic controlled light conditions (Photopic)
- Mesopic controlled light conditions (Mesopic)


In the case of the dynamic pupillometry, recording of the state of the pupil is started, first in mesopic conditions, then photopic and then mesopic again. The data on the diameters measured are recorded and shown in the **"Measurement"** section.

For the dynamic acquisition, a sequence of images is recorded and allows you to "review" the evolution of the pupil through the various different light conditions to which it is subjected. In the pupillometry acquisition in static controlled light conditions: photopic and mesopic, certain frames are saved, which you can display by scrolling the associated gallery in the Pupil \rightarrow Measurements section.

WARNING: With blue eyes, acquisition of pupillometry in mesopic lighting conditions can be difficult to accomplish. In this case, we suggest acquiring the mesopic data through dynamic pupillometry.



Figure 37

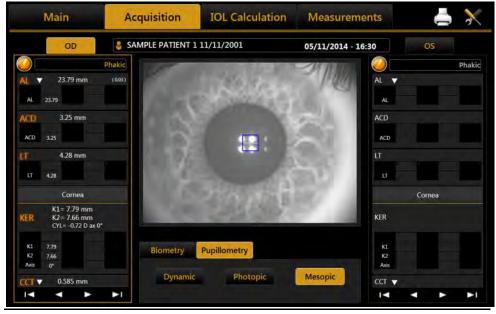

Acquisition in mesopic controlled light conditions

Figure 38

button.

11.12 Report printing

After every measurement, you can print the corresponding report or print all the measurements made in the

current exam. In the top-right corner of the screen, press on the

As shown in Figure 39 you can now select from the left column which report to print and also with which surgeon preset (from the **"Surgeon"** box).

You can print directly to an external printer or a USB drive and also print to both simultaneously ("Both" option).

REPORT SELECTION		SURGEON				
 ✓ Aladdin □ Measurements □ Pupil 	OLSEN IOL CALCULATOR Olsen IOL Olsen Toric IOL	Surgeon Generic				
□ IOL Custom reports □ CustomReportTest	BARRETT IOL CALCULATOR Barrett Universal II Formula Barrett Toric Calculator Barrett True K Barrett True K Toric Calculator Barrett Rx Formula	OUTPUT DEVICES HP Universal Printing PCL 6 USB Drive Export Both Network Folder \ 				
ALADDIN REPORT OPTIONS						
3 zones Anels None O Sim-K	 Pupil O Meridians O Emimeridians 					
Screenshot	Cancel Print					

Figure 39

In addition to these two options you can select a network folder as destination of report to print. For network folder settings, see 13.5.1.

The "Aladdin Report Options" allow to define which overlay will be printed in the topography maps images of the "Aladdin". Refer to sections 12.1 and to 13.2 for details about this options.

By clicking on the Screenshot button (Figure 39) you can open the print screen preview and, as shown in Figure 40, send to print the document as represented in the preview.

Main	Acquisition	IOL Calculation	Measurements	📥 🔀
OD	PRINT SCREEN PREVIE	w	-	OS
AL 🔻 23.48 mm	Main Main Acc	uisition IOL Calculation Me	asuremanta 📥 🔭	
AL 23.48	Parte	con Demo 10/12/1980 10/05	/2015 - 13.49 000 Phake	
ACD 3.98 mm	(C), w 23.48, mm (2.44)			
ACD 3.98	4610 828 mm		42	
[T] 3.85 mm	JI 385 mm			
LT 3.85	K0-49-		Consu ALP	
Cornea	KEN KE - LED B (KP)			Cornea
K1= 7.51 mm KER: K2= 7.29 mm CYL= -1.35 D ax	60 9.00 60 9.00 60 10 60 10		es es es	
K1 7.51				
K2 7.29			Cancel Print	
Axis 0*				
(CCT ▼ 0.523 mm				

Figure 40

11.12.1 Available Printers

The printing form is showing a list of available printers. The available printers are the one installed on the Operating System. Refer to <u>"Appendix: Installing an external printer"</u> or ask to your technical support in order to have the desired printer installed.

The application pre-selects always the last used printer.

REPORT SELECTION		SURGEON
 Aladdin Measurements Pupil IOL Custom reports CustomReportTest	OLSEN IOL CALCULATOR Olsen IOL Olsen Toric IOL BARRETT IOL CALCULATOR Barrett Universal II Formula Barrett Toric Calculator Barrett True K Barrett True K Toric Calculator Barrett Rx Formula	Surgeon Generic OUTPUT DEVICES Printing PCL 6 HP Universal Printing PC E C E C E C M C M C M C M C M C M C M C M
ALADDIN REPORT OPTIONS	🗆 Pupil	
● None O Sim-K	O Meridians O Emimeridians Cancel Print	

11.12.2 Custom Reports


If the unit has been provided with Customized Reports they will be available to be selected for printing or exporting in the Printing form. In order to obtain custom reports contact you technical assistance.

REPORT SELECTION		SURGEON
 Aladdin Measurements Pupil IOL Custom report Custom ReportTest 	ts OLSEN IOL CALCULATOR Olsen IOL Olsen Toric IOL BARRETT IOL CALCULATOR Barrett Universal II Formula Barrett Toric Calculator Barrett True K Barrett True K Barrett Rx Formula	Surgeon Generic OUTPUT DEVICES HP Universal Printing PCL 6 USB Drive Export O Both O Network Folder \
ALADDIN REPORT OPTION	<u>S</u>	
☐ 3 zones☐ Ane● None○ Sim		5
Screenshot	Cancel Print	

11.13 Data Exportation

you can select one or more destinations for exportation.

After every measurement, you can export the corresponding reports or xml date made in the current exam.

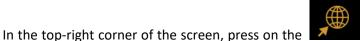

button. The popup of Figure 41 is show where

SEND TO....
DICOM Storage Server
i - base
PhacoOptics®
XML
Cancel Export

Figure 41

Currently available destinations are:

• **DICOM Storage SCP Server**, the selected reports are sent to the designed DICOM storage location according to prior defined settings. Refer to DICOM configuration for further details. Select the desired reports to save on the Storage location.



You will receive the confirmation message in case of successful or unsuccessful transfer.



If Storage Commitment is configured you will receive the confirmation message in case of successful or unsuccessful transfer.



• IMAGEnet 6/IMAGEnet i-base, the selected reports and data are sent to IMAGEnet i-base or IMAGEnet 6 if activated and configured from the settings. If both are enabled the destination is IMAGEnet 6.

Maim	SURGEON		🔎 📥 🗡
OD	Surgeon Generic		-OSI
Data IOI			Olsen
Surgeon Surgeon Generic	IMAGENET 6 Report Export Selection		YL (D) -1.45 ax 8
Target (D) 0 STA Model Spherical Equivalent Cylindrical Power (D) Spherical Power (D) Axis of Placement (1) Expected Refraction	Measurements Pupil IOL Post Refractive IOL Toric IOL Custom reports Custom ReportTest	OLSEN JOL CALCULATOR Olsen JOL Olsen Toric IOL BARRETT JOL CALCULATOR Barrett Universal II Formula Barrett Toric Calculator Barrett True K Barrett True K Barrett Rx Formula	TW (mm) 11.69 JP 0 (mm) 4.11
	Cance	Export	Next

- PhacoOptics[®]: <u>http://www.phacooptics.net/</u>, export biometry data in xml format to a network shared folder, to be imported in PhacoOptics[®] software. Refer to section 13.5 and to PhacoOptics[®] software instruction for further details on configuration.
- **XML:** create XML file with biometry data and IOL calculations (also images of the eye topography) that is exported to the configured network shared folder. Refer to section 13.5.

When the selection of the destination you can press "**Export**" to perform exportation of data to the destinations, or "**Cancel**" to just close.

11.14 IOL calculation

ALADDIN HW3.0 also includes a section for calculation of the intraocular lenses (IOL Calculation).

In order to perform the intraocular lens power calculation, the available power interval, the increments and the calculation constants must be provided for each type of formula and lens. These, however, do not depend solely on the type of lens and the calculation formula, but are also closely linked to the measuring technology and the surgical techniques used. It is strongly recommended that the user optimize the IOL constants in clinical practice and the type of device used for acquisition of the biometric data.

Main Acquisitio		sition	IOL Calc	ulation	n Measurements 💭					
OI		торсо	N DEMO 01,	/01/1950		10/02/.	2015 -	17:55	(DS
Data	IOL Calc	ulation	Toric IOL	Calculation	Post Ref	ractive IO	L	Barrett		Olsen
Biometry					Keratomet	y				
Source	Aladdin	Man	iual		Source	Aladdin		Manual		
AL	23.73	mm		-	кі	40.74	D	8 •	/	
ACD	3.14	mm	1		К2	42.19	D	98 •		
ιτ	4.04	mm		*	CYL	-1.45	D	8 •		
сст	0.544	mm			Index	1.3375	Ur	it D 🔻		
wtw	11.69	mm								
PUPØ	4.11	mm								
Acquisition	Optical									
Lens	Phakic									
Vitreous	Natural									

Figure 42

The screen has various sections, which we will explore in detail below:

- Data
- IOL Calculation
- Toric IOL Calculation
- Post Refractive IOL
- Barrett (Barrett IOL Calculator module)
- Olsen (Olsen IOL Calculator module)

The first time a surgeon enters in IOL Calculation, it will appear the panel shown in Figure 43 that contains the Disclaimer regarding the usage of the IOL Calculation. The Disclaimer will appear every time you enter in IOL Calculation or you want to print a IOL report, unless you check the box below before clicking the OK button. The Disclaimer is also replicated in the IOL settings (see 13.4.1), where you can even reactivate its appearance at every IOL usage.

Disclaimer

This IOL calculator is intended to be an additional tool to help physicians in selecting the appropriate IOL for a particular patient. The tool is intended to be used in concurrency with a comprehensive ophthalmic examination, specific diagnostic tests and measurements for patient who need to perform cataract surgery. The calculator's results are not intended to serve as surgical or medical instruction from TOPCON or be definitive, nor TOPCON can guarantee that the calculator will be accurate in every moment. Medic or Physician who uses this tool must arrive at their own decisions regarding the patients and they are responsible for the refractive outcome. By using this tool, you agree the conditions and hold TOPCON from any claims you may arising out of your use of this calculator.

Don't show this message again

OK

Figure 43

11.14.1 Data

In this section there is a summary of the measurements performed with the instrument.

The screen displayed for the "Data" tab is as shown in Figure 42. As you can see, it has two sections:

- Biometry: detailing the data on the ocular biometry
- Keratometry: detailing the patient's keratometric data

From the **"Source"** field (present in both the Biometry and Keratometry sections, Figure 42), you can choose to use the measurements of the **"Aladdin"** instrument as source, or to enter them manually by pressing the **"Manual"** button. In this case, a panel opens (Figure 44) where you can enter the data manually using the numeric keypad.

BIOMETRY SOURCE MA	ANUAL	7	8	9	CURRENT EYE : OD KERATOMETRY SOURCE MANUAL
AL	23.73 mm	4	5	6	K1 8.28 mm 8 °
ACD LT	3.14 mm 4.04 mm	1	2	3	K2 8.00 mm 98 °
сст	0.544 mm		0	-	Index 1.3375 🔻
WTW PUP	11.69 mm 4.11 mm	~	ESC	ОК	
Acquisition	Optical 🔻				
					Cancel Save

Figure 44

In Figure 44 (Manual Keratometry tab), you can enter values both in Diopters or Millimeters. The values will be automatically recognized according to a specific range. This range goes from 6.75mm to 9.64mm, (from 35Diopters to 50Diopters). Also every mm/D conversion from now on will be performed on the basis of the current index in this section.

In Figure 44 you can manually enter biometry data; knowing that an external instrument does not always work in Optical mode, you can also select acoustic measurement mode.



WARNING: The responsibility for any data entered and checked manually lies exclusively on the user. **WARNING:** Using data from acoustic instruments also means that the constant of every IOL must be

optimized for those kinds of instrument; at present, it is more common to find online only databases of lenses optimized for optical interferometry instruments.

11.14.2 Spherical IOL calculation

The section where you process the data collected for the calculation of the best intraocular lens is divided into three main parts:

The **IOL Calculation** section is fundamental in calculating the total spherical power that will compensate the removal of the crystalline in cataract surgery. This calculation, depending on the case, will be sufficiently precise to guarantee the patient optimal vision.

In certain cases, however, if the corneal astigmatism exceeds a disturbance value, i.e. if faced with a medium or medium-high astigmatism, it is advisable to implant an intraocular lens that also takes this factor into account (refer to the following section).

In the **"Toric IOL Calculation"** section, a surgeon who has established, by means of topography, the need for an astigmatic correction, can use toric lens calculation. The calculation is made on the basis of the ideal spherical IOL and the patient's characteristics. The software then, on the basis of the corneal astigmatism, recalculates the one on the IOL plane, also considering the possible astigmatism induced by surgery.

The Toric IOL is thus chosen, which guarantees post-operation refraction with minimum residual astigmatism.

In the "**Post Refractive IOL**" section you can calculate the intraocular lenses for patients who have already had refractive surgery to correct myopia or hypermetropia.

Maii	n	Acqui	sition	IOL Calc	ulation	Measurements			≜ 🗙			
	OD OT TOPCON DEMO 01/01/1950 10/02/2015 - 17:55 OS											
Data	IOL Ca	alculation	Toric IOL (alculation	Post Refr	active IOL	Barret	tt	Olsen			
Surgeon				Measurements								
Surgeon Ge	eneric		•	AL (mm)	23.73	K1 (D)	40.74	CYL (D)	-1.45 ax 8°			
Target (D)				ACD (mm)	3.14	K2 (D)	42.19	WTW (mm) 11.69			
0				LT (mm)	4.04	CCT (mm)	0.544					
, i												
Oculentis	\bullet	Oculentis	▼	Oculentis 🔹 🔻		Oculentis 🔹 🔻		Oculentis 🔹 🔻				
L-303	•	L-313	•	LS-313 M	LS-313 MF30 V LS-313 M		F30 🔻	LS-313 N	5-313 MF30 🛛 🔻			
Barrett Uni	iversal II 🔻	Barrett Ur	niversal II 🔻	Haigis	igis 🔹 🔻 Barrett Universal II 🔻		SRK/T 🔻					
IOL @ Target 21.97	LF = 1.517 A = 118.300	IOL @ Target 21.77	LF = 1.412 A = 118.100	IOL @ Target 22 .44	A0 = 0.950 A1 = 0.400 A2 = 0.100	IOL @ Target 22.21	LF = 1.640 A = 118.530	IOL @ Target 22.12	A = 118.500			
IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)			
21.00	0.73	21.00	0.58	21.50	0.69	21.00	0.89	21.00	0.80			
21.50	0.35	21.50	0.21	22.00	0.32	21.50	0.52	21.50	0.45			
22.00	-0.02	22.00	-0.17	22.50	-0.04	22.00	0.15	22.00	0.09			
22.50	-0.40	22.50	-0.56	23.00	-0.42	22.50	-0.22	22.50	-0.28			
23.00	-0.79	23.00	-0.95	23.50	-0.79	23.00	-0.60	23.00	-0.64			
									Reset			

Figure 45

Figure 45 above shows the SW environment for entering data for the "IOL Calculation".

In the *"Surgeon"* field, you can choose which surgeon will perform the IOL implant and any customization of the constants or preset of the preferred lenses and formulae will be applied on this basis.

In *"Target"* field the target refractive value for the Post-Op must be inserted.

The *"Measurements"* field summarizes the measurement data.

From the drop-down menu, select the IOL manufacturer and model, as well as the preferred formula with which to calculate the best lens.

Once this data has been entered, the most appropriate lens can be chosen at the discretion of the surgeon.

Pressing "Reset" will reset the initial preset conditions.

11.14.3 Toric IOL calculation

Toric IOL calculation is divided into two main steps. The first one consists on the calculation of the Spherical Equivalent Power; in the second one you can select the toric IOL that produce the best correction.

Main Acquisition		sition	IOL Calc	ulation	Measur	ements			
C	D	👌 торсог	N DEMO 01/0)1/1950		10/02/201	5 - 17:55		OS
Data	IOL Ca	alculation	Toric IOL (alculation	Post Refr	active IOL	Barre	ett	Olsen
Surgeon				Measurements					
Surgeon Gen	eric		T	AL (mm)	23.73	K1 (D)	40.74	CYL (D)	-1.45 ax 8°
Target (D)	SIA (D)	IL (°)		ACD (mm)	3.14	K2 (D)	42.19	WTW (n	nm) 11.69
0	0	and the second division of the second divisio	98	LT (mm)	4.04	CCT (mm)	0.544		
Oculentis		Alcon							
LS-313 Tx	•	Acrysof S	N6AT 🔻		Ţ		•		
Haigis	•	SRK/T			- -		•		
10L @ Target 22.48	A0 = 0.970 A1 = 0.400 A2 = 0.100	IOL @ Target 23.00	A = 119.200	IOL @ Target		IOL @ Target		IOL @ Tar	get
IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D) REF (D)
21.50	0.72	22.00	0.70						
22.00	0.35	22.50	0.35			No.			
22,50	+0.02	23.00	0.00	-					
23.00	-0.39	23.50	-0.35	-		-	-		
23.50	-0.76	24.00	-0.71						

Figure 46

Figure 46 shows the first-step interface that has quite the same structure as the normal IOL calculation. The available toric lenses you can select come from a list of models whose calculation constants have been published by their manufacturer. The user can in case insert new toric manufacturers and/or models inside toric IOL settings section (see 13.4.3)

In addition to choosing the **"Target"**, you need to specify also the **"Surgical Induced Astigmatism (SIA)"** and **"Incision Location (IL)"**. The former identify the astigmatism (in diopters) induced by the incision while the latter identify the surgical incision axis. After having selected the toric IOL model and one of the available formulas, a values table from which to choose the **Spherical Equivalent Power** is obtained. Once you choose a lens, pressing **"Next"** at the bottom right, you enter in the second-step of toric IOL calculation (Figure 47).

OD	В ТОРСС	ON DEMO 01	/01/1950		10/02/20	15 - 17:55		OS	
Data	IOL Calculation	Toric IOL	Calculation	Post Re	fractive IOL	Ba	rrett	Olse	n
Surgical Pre Op	Data	0	Measurements						
SEQ	SIA		AL (mm)	23.73	K1 (D)	40.74	WTW (n	im) 11.69	
23.00		0	ACD (mm)	3.14	K2 (D)	42.19	CYL (D)	-1.45 ax	8 '
Formula	IL		LT (mm)	4.04	CCT (mm)	0.544			
	SRK/T	98	Expected Post O	o Cornea	🗖 A	bulafia-Koch	Correction		
	A = 119.200		K1 (D)	40.74	K2 (D)	42.19	CYL (D)	-1.45 ax	8 '
Toric IOL									
Model	Alcon Acr	Sof SN6AT3	Availab	le Toric Lens Lens	es Res Astigm		OD 105	90 75 60	1
Spherical Eq	uivalent Power (D)	23.00		n.a.	ø		135	1 1	5 30
Cylindrical P	ower (D)	1.50)	of SN6AT2	-0.73 D @ 8		165		15
Spherical Po		22.25	Acrys	of SN6AT3	-0.38 D @ 8		80 -		1-0
			Activa	of SN6AT4	-0.15 D @ 98		_ \	1 /	-
Axis of Place		98	2	IOL Ideal	Designed and	.04	Temporal	1	Nasa
Expected Re	fraction +0.18D	-0.38 D @ 8					ora		

Figure 47

"Measures" and "Surgical Pre Op Data" frames summarize the values used in the first-step calculation.

"Expected Post Op Cornea" frame gives information about the post surgery patient eye Keratometry, taking into account the aforementioned SIA and IL.

If the option "**Abulafia-Koch correction**" is selected the astigmatism is corrected in the "**Expected Post Op Cornea**" by taking into account, together with SIA and IL, the nomogram-based correction. Refer to "<u>Abulafia</u> <u>A, Koch DD, Wang L, Hill WE, Assia EI, Franchina M, Barrett GD: New regression formula for toric intraocular</u> <u>lens calculation. – Journal of Cataract & Refractive Surgery, 2016 – Elsevier</u>". If the Toric IOL calculation is performed using the Abulafia-Koch Correction it is reported in the corresponding Toric IOL report.

As a result, the **"Toric IOL"** frame, immediately below, details the best toric lens computed automatically by the system for the manufacturer and model selected previously in the first-step.

From **"Available Toric Lenses"** table you can choose also a different cylinder value for the lens, based on the Residual Astigmatism you want to achieve (under-correction/overcorrection). In particular, the best toric lens value is shown in the central row and (if available) the ones that under-correct above the central row, the ones that overcorrect below.

At the right side, you can find an image that illustrates the ideal position of the IOL once the implant is in place and in red the incision location angle.

Under the table, the small icon opens the Toric Rotation Misalignment Simulator (Figure 48).


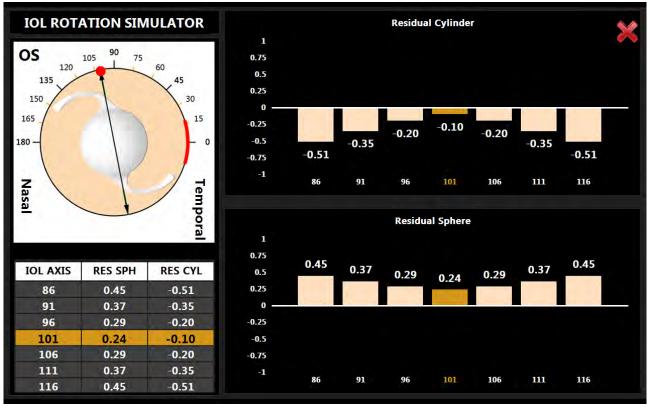


Figure 48

This simulation shows the impact of a wrong axis placement that can occur during surgery and how it affects the residual sphere and cylinder refraction of the patient. The simulator starts at the correct axis placement and displays in the left bottom table a series of misalignment angles close to the ideal one. Selecting a different row you can see the new residual sphere and cylinder at the selected axis.

In Figure 49 is shown the selection of a ten degrees misalignment, with the new axis selected in orange and the correct one in green.



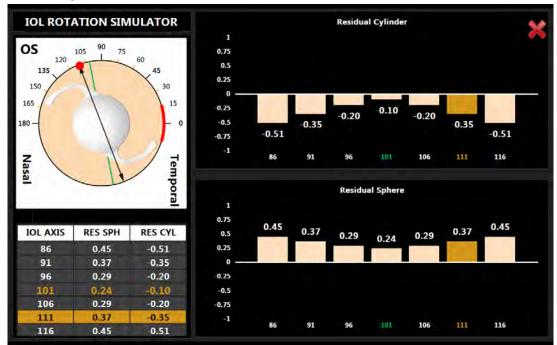

70

11.14.4 Post Refractive IOL calculation

In this section, you can calculate intraocular lenses also for patients who have undergone refractive surgery to correct myopia or hyperopia by using the Camellin-Calossi formula and the Shammas no-history formula (Figure 50).

These formulae are used with patients who have had prior refractive surgery. Each such patient is unique and results may vary widely. You should interpret all IOL power recommendations with caution. In this environment, you need to manually enter certain fundamental data.

Main Acquis		quisitio	ition IOL Calculation		Measurements		A	5	3	
0	р 🚦 ТО	PCON DEM	MO 01	/01/1950		10/02/201	5 - 17:55		OS	
Data	IOL Calculation		Toric IOL Calculation		Post Refractive IOL		Barre	ett	Olsen	
Surgeon Surgeon Gene LT (mm) 4.04	eric			Measurements AL (mm) ACD (mm) LT (mm)	23.73 3.14 4.04	K1 (D) K2 (D) CCT (mm)	40.74 42.19 0.544		-1.45 ax mm) 11.69	
Keratorefractive Surgery Information Form Correction Type				las Available		Calculation Guidelines Guide on the use of the Post Refractive				
Myopic Type RK			□ Camellin Calossi □ Shammas No History			IOL Calculation Formula: Camellin Calossi formula can be used in all types of keratorefractive surgery shown on the screen. Shammas No History can be used only on Myopic or Hyperopic RK, PRK, LASIK and LASEK.				
							_			

Figure 50

The first is the surgeon who performs the operation. As in IOL calculation, the constants may differ from one surgeon to another.

The second data is the "LT", i.e. the crystalline lens thickness that can be modified manually.

Next, select the correction type between the options shown below:

- Myopic
- Hyperopic
- Unknown

If the correction type is unknown:

- it is not possible to select the surgery type
- it is not possible to choose the Shammas No-History formula
- to use the Camellin Calossi formula, you must insert in the Input Data the Pachimetry values and the optic zone diameter.

If the correction type is myopic or hyperopic, you have to select the surgery type performed on the patient from those listed:

- 1. Radial Keratotomy (RK)
- 2. Photo Refractive Keratectomy (PRK)

- 3. Lasik
- 4. Lasek
- 5. LK
- 6. РТК
- 7. Unknown

In case of **Radial Keratotomy, Photo Refractive Keratectomy, Lasik and Lasek** you need to insert in the **"Refractive Change"** frame the correct ametropia type and the correction factor obtained by the operation ("SIRC").

In case of **LK**, **PTK or Unknown** only the Camellin Calossi formula is available and you need to insert the current **Pachymetric data** as well as the diameter of the optical zone to improve the accuracy of the final calculation (Figure 51).

The Unknown option must be selected every time that you don't know the type of surgery or one of the associated information.

For example if you know your patient has undergone Radial Keratotomy or Photo Refractive Keratectomy or Lasik or Lasek but you don't know the SIRC value, select **Unknown** and insert the pachimetry values.

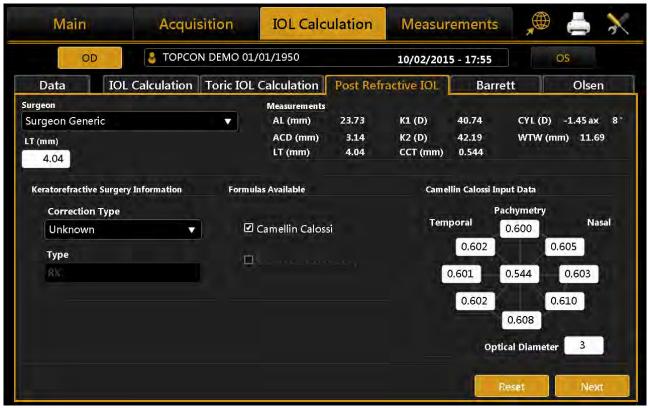


Figure 51

By pressing **"Next"** you move on to the final diagram of the calculation. Here you decide on the **"Target"** and select the lens make and model.

If you highlight the lens selected, the result is memorized and highlighted on the report.

Main	Acquisition	IOL Calculation	Measurements	🔎 📥 💥
OD	E TOPCON DEMO 01/	/01/1950	10/02/2015 - 17:55	OS
Data	IOL Calculation T	oric IOL Calculation	Post Refractive IOL	Barrett Calculator
Surgeon		Measures		
Surgeon Generic		AL (mm) 23.73	K1 (D) 40.74	CYL (D) -1.45 ax 8°
Target (D)		ACD (mm) 3.14	K2 (D) 42.19	WTW (mm) 11.69
0		LT (mm) 4.04	CCT (mm) 0.544	
.ZEISS	▼ Aaren ▼	• (·	v	▼
CT 47S(Acri.Lyc 47S)	Scientific AQUA 4 Y RM		•	•
Camellin Calossi	Shammas	•	▼	
IOL @ Target A = 117.6	00 IOL @ Target A = 118.500) IOL @ Target	IOL @ Target	IOL @ Target
22.26	24.16			
IOL (D) REF (D)	IOL (D) REF (D)	IOL (D) REF (D)	IOL (D) REF (D)	IOL (D) REF (D)
21.50 0.55	23.00 0.85			
22.00 0.19	23.50 0.49			
22.50 -0.17 23.00 -0.54	24.00 0.12 24.50 -0.25			
23.50 -0.90	25.00 -0.63			
			R	eset Back

Figure 52

The final result of the Post-Op calculation is shown in Figure 52 with the suggested lenses highlighted in yellow in each table.

11.14.5 Barrett Calculator

Barrett Calculator is an optional module: ask your distributor for details.

The Barrett Calculator integrates the Barrett IOL Calculator v1.05. The calculation methods are the following:

- Universal Formula II: Barrett Universal II Formula v1.05, for alle yes regardless of axial length
- **Toric Calculator:** Barrett Toric Calculator v1.05, for correction of pre-existing corneal astigmatism with Toric IOLs
- True K: Barrett True K Formula v1.05, for eyes with prior myopic or hyperopic LASIK/PRK/RK
- **True K Toric:** Barrett True-K Toric Calculator v1.05, for eyes with prior myopic or hyperopic LASIK/PRK/RK and corneal astigmatism
- **Rx Formula:** Barrett Rx Formula v1.05, for IOL exchange and piggy back IOLs based on refraction after cataract surgery


All the formulas are based on the Barrett Universal II Formula.

Mair	a)	Acqui	sition		IOL Calc	ulation	Measur	ements	* - *
	OD	5 торсо	N DEMC	01/0	01/1950		10/02/201	5 - 17:55	OS
Data	IOL Cal	ulation	Toric	IOL (Calculation	Post Ref	ractive IOL	Barret	t Olsen
Surgeon					Measurements				a second s
Surgeon Gei	neric		•		AL (mm)	23.73	Kf (D)	40.74	CYL (D) -1.45 ax 8°
Target (D)	0 SIA (D)	0	пс	98	ACD (mm) LT (mm)	3.14 4.04	Ks (D) CCT (mm)	42.19 0.544	WTW (mm) 11.69
		Universal	II Formul	la U	niversal II To	ric T	rue K	True K Tor	ric RX Formula
Model		Oculen	tis LS-31	.3 Tx		Toric Lenses	Res Astigm	Right	Eye so
Spherical I	Equivalent Powe	er (D)	2	2.00				т	35
de la seletaria	l Power (D)			0.00			Carl and	e m 180-	102
Cylindrica		-			Non	and the second sec	the second se	P	
Spherical I	Power (D)		2	2.00	1010-00		-0.85 D @ 12°	P A	
Spherical	Power (D) acement (°)	-	2	102	LS-3	13 T1 -	-0.25 D @ 102°	a a	no ja ca
Spherical	acement (°)	+0.33 D -0		102	LS-3	13 T1 -	-0.25 D @ 102° -0.81 D @ 102°	r a l	no n
Spherical I Axis of Pla	acement (°)	+0.33 D -(102	LS-3	13 T1 - 13 T2 -	-0.25 D @ 102° -0.81 D @ 102°		Hecomorate Ara



11.14.5.1 Universal Formula II (Barrett)

Barrett Universal II Formula v1.05, for all yes regardless of axial length.

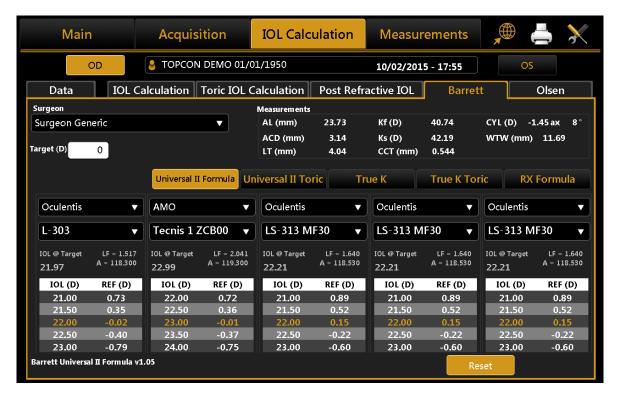

In the *"Surgeon"* field, you can choose which surgeon will perform the IOL implant and any customization of the constants or presetting of the preferred lenses will be applied on this basis.

In *"Target"* field the target refractive value for the Post-Op must be inserted.

The "Measurements" field summarizes the measurement data.

From the drop-down menu, select the IOL manufacturer and model with which to calculate the best lens. Once this data has been entered, the most appropriate lens can be chosen at the discretion of the surgeon. The latter is highlighted in orange. Once selected, the lens will be memorized as the preferred one and will be shown highlighted on the report printout.

Pressing **"Reset"** will reset the initial preset conditions.

11.14.5.2 *Toric Calculator (Barrett)*

Barrett Toric Calculator v1.05, for correction of pre-existing corneal astigmatism with Toric IOLs.

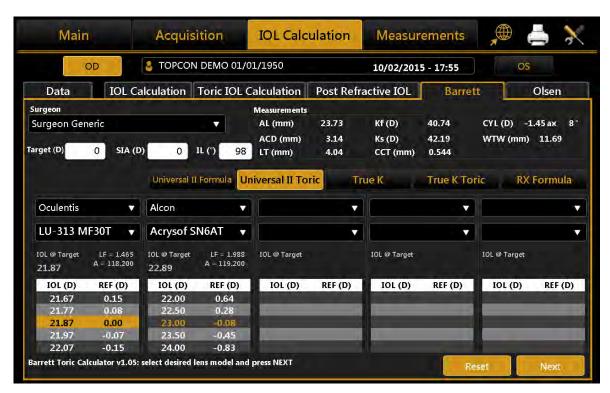


Figure 54

Toric Calculator (Barrett) is divided into two main steps. The first one consists on the calculation of the Spherical Equivalent Power; in the second one you can select the toric IOL that produce the best correction.

Figure 60 shows the first-step interface that has quite the same structure as the spherical IOL calcualtion. The available toric lenses you can select come from a list of models whose calculation constants have been published by their manufacturer. The user can in case insert new toric manufacturers and/or models inside toric IOL settings section (see 13.4.3)

In addition to choosing the **"Target"**, you need to specify also the **"Surgical Induced Astigmatism (SIA)"** and **"Incision Location (IL)"**. The former identify the astigmatism (in diopters) induced by the incision while the latter identify the surgical incision axis.

After having selected the toric IOL model, a values table from which the **Spherical Equivalent Power** is obtained. Once you choose a lens, pressing **"Next"** at the bottom right, you enter in the second-step of toric IOL calculation (Figure 55).

	OD		ТОРСО	N DEMO	D 01/0	01/1950		10/02/201	5 - 17:55		OS
Data	10	IOL Calc	ulation	Toric	IOL (Calculation	Post Ref	ractive IOL	Barre	ett	Olsen
Surgeon					_	Measurements					
Surgeon Ger	neric			•		AL (mm)	23.73	Kf (D)	40.74	CYL (D) -1.45 ax 8°
Target (D)	0	SIA (D)	0	цс	98	ACD (mm)	3.14 4.04	Ks (D) CCT (mm)	42.19 0.544	WTW	(mm) 11.69
					50						
			Universal	II Formu	1	niversal II To	ric I	rue K	True K To	oric	RX Formula
Model					la <mark>U</mark>	niversal II To	ric T	rue K	-		RX Formula
			Alcon Ac	rysof St	la <mark>U</mark>	niversal II To Available		rue K Res Astigm	-	ht Eye	RX Formula
Model	Equiva	ent Powe	Alcon Ac	rysof St	la U	niversal II To Available	Toric Lenses		Rigi		RX Formula
Model Spherical F	Equival Powe	lent Powe r (D)	Alcon Ac	rysof SI	la U N6AT 23.00	niversal II To Available La Non AcrySot	Toric Lenses ens Toric SN6AT2	Res Astigm -0.85 D @ 12° -0.14 D @ 12°	Rigt T m m 190- o r		RX Formula
Model Spherical E Cylindrical	Equival Powe	lent Powe r (D) (D)	Alcon Ac	rysof SI	N6AT 23.00 1.00	niversal II To Available La Non AcrySol AcrySol	Toric Lenses ens Toric SNGAT2 SNGAT3	Res Astigm -0.85 D @ 12°	Rigi		RX Formula

Figure 55

As a result, the **"Toric Calculator"** frame, immediately below, details the best toric lens computed automatically by the system for the manufacturer and model selected previously in the first-step.

From **"Available Toric Lenses"** table you can choose also a different cylinder value for the lens, based on the Residual Astigmatism you want to achieve (under-correction/overcorrection). In particular, the best toric lens value is shown in the central row and (if available) the ones that under-correct above the central row, the ones that overcorrect below.

At the right side, you can find an image that illustrates the ideal position of the IOL once the implant is in place and the incision location angle.

11.14.5.3 *True K (Barrett)*

Barrett True K Formula v1.05, for eyes with prior myopic or hyperopic LASIK/PRK/RK.

Main	Acquis	ition	IOL Calc	ulation	Measur	rements	, the second sec	
OD	🕹 TOPCON	I DEMO 01/0	1/1950		10/02/201	5 - 17:55	0	5
Data IOL Ca	alculation	Toric IOL C	alculation	Post Refra	active IOL	Barret	tt 👘	Olsen
Surgeon			Measurements					
Surgeon Generic			AL (mm) ACD (mm)	23.73 3.14	Kf (D) Ks (D)	40.74 42.19	CYL (D) - WTW (mm)	1.45 ax 8° 11.69
farget (D) 0			LT (mm)	4.04	CCT (mm)	0.544		
HISTORY Correction type:	Oculentis		Oculentis	×.	Oculentis	•	Oculentis	
Correction type:	L-303	•	L-303	•	L-313	•	L-313	
Myopic Lasik 🔻	L 303		L 303	Y	[L 313		L 313	
Pre-Lasik Refraction -4.50	IOL @ Target 23 .48	LF = 1.517 A = 118.300 LF TK = 1.880	10L @ Target 23.48	LF = 1.517 A = 118.300 LF TK = 1.880	IOL @ Target 23.26	LF = 1.412 A = 118.100 LF TK = 1.770	IOL @ Target 23.26	LF = 1.412 A = 118.100 LF TK = 1.770
Post-Lasik Refraction	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
0.25	22.50	0.72	22.50	0.72	22.50	0.57	22.50	0.57
- Automation	23.00	0.35	23.00	0.35	23.00	0.20	23.00	0.20
No history	23,50	-0.02	23.50	-0.02	23.50	-0.18	23.50	-0.18
	24.00	-0.39	24.00	-0.39	24.00	-0.56	24.00	-0.56
TrueK= 40.92 D Corr.= -4.52 D	24.50	-0.77	24.50	-0.77	24.50	-0.94	24.50	-0.94



In *"Target"* field the target refractive value for the cataract Post-OP must be inserted. The *"Measurements"* field summarizes the measurement data.

In the "*HISTORY*" section select the correction type performed in the preceeding Refractive Surgery:

- Myopic Lasik
- Hyperopic Lasik
- Radial Keratotomy

Insert the measured Refraction (**Pre-Lasik Ref., in diopters**) before the Refractive Surgery and the measured Refraction (**Post-Lasik Ref., in diopters**) after the Refractive Surgery, accordingly the selected correction type.

Pre-Lasik Ref. must be negative for Myopic Lasik and Radial Keratotomy corrections, while must be positive for Hyperopic Lasik.

Otherwise select "**No History**" if pre and post refractive surgery measurements are not available in order to obtain an estimate of the correction amount based on the correction type and the eye biometry data.



From the drop-down menu, select the IOL manufacturer and model with which to calculate the best lens. Once this data has been entered, the most appropriate lens can be chosen at the discretion of the surgeon. The latter is highlighted in orange. Once selected, the lens will be memorized as the preferred one and will be shown highlighted on the report printout.

Pressing "Reset" will reset the initial preset conditions.

11.14.5.4 *True K Toric (Barrett)*

Barrett True-K Toric Calculator v1.05, for eyes with prior myopic or hyperopic LASIK/PRK/RK and corneal astigmatism.

Main	Acquis	tion	IOL Calc	ulation	Measur	ements	*	3
OD	STOPCON	DEMO 01/0	1/1950		10/02/201	5 - 17:55	OS	in the second
Data IOL C	alculation	Toric IOL C	alculation	Post Refr	active IOL	Barret	t	Olsen
Surgeon			Measurements					
Surgeon Generic			AL (mm)	23.73	Kf (D)	40.74	CYL (D) -1	L.45 ax 8°
and the state			ACD (mm)	3.14	Ks (D)	42.19	WTW (mm)	11.69
Target (D) 1 SIA (I	0.00	IL (°) 98	LT (mm)	4.04	CCT (mm)	0.544		
	Universal II	Formula Ur	niversal II To	ric Ti	rue K	True K To	ric RX	Formula
HISTORY	Alcon	-		•		•		
Correction type:	Amusto	NCAT						
Myopic Lasik 🔻	Acrysof S	N6AT 🔻			JI	<u> </u>		
Pre-Lasik Refraction -4,50	IOL @ Target 23.03	LF = 1.988 A = 119.200 LF TK = 2.340	IOL @ Target		IOL @ Target		IOL @ Target	
Post-Lasik Refraction	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
0.20	22.00	1.71	_					
□ No history	22.50	1.37	-	_	1	-		_
- Hornordy	23.00 23.50	1.02	100	_				_
TrueK= 40.93 D Corr.= -4.47 D	Contract of the second s	0.32	1.	_	-			_

Figure 57

True K Toric Calculator (Barrett) is divided into two main steps. The first one consists on the calculation of the Spherical Equivalent Power; in the second one you can select the toric IOL that produce the best correction.

Figure 57 shows the first-step interface that has quite the same structure as the spherical IOL calcualtion. The available toric lenses you can select come from a list of models whose calculation constants have been published by their manufacturer. The user can in case insert new toric manufacturers and/or models inside toric IOL settings section (see 13.4.3)

In addition to choosing the **"Target"**, you need to specify also the **"Surgical Induced Astigmatism (SIA)"** and **"Incision Location (IL)"**. The former identify the astigmatism (in diopters) induced by the incision while the latter identify the surgical incision axis.

In the "*HISTORY*" section select the correction type performed in the preceeding Refractive Surgery:

- Myopic Lasik
- Hyperopic Lasik
- Radial Keratotomy

Insert the measured Refraction (**Pre-Lasik Ref., in diopters**) before the Refractive Surgery and the measured Refraction (**Post-Lasik Ref., in diopters**) after the Refractive Surgery, accordingly the selected correction type.

Pre-Lasik Ref. must be negative for Myopic Lasik and Radial Keratotomy corrections, while must be positive for Hyperopic Lasik.

Otherwise select "**No History**" if pre and post refractive surgery measurements are not available in order to obtain an estimate of the correction amount based on the correction type and the eye biometry data.

☑ No History	
TrueK= 41.01 D Corr.= -3.74 D	

After having selected the toric IOL model, a values table from which the **Spherical Equivalent Power** is obtained. Once you choose a lens, pressing **"Next"** at the bottom right, you enter in the second-step of True K toric IOL calculation (Figure 58).

Mair	n	Acqui	sition	IOL C	alculation	Measur	rements	🔎 📥 💸
	OD	👌 торсо	N DEMO 0	1/01/1950		10/02/201	5 - 17:55	OS
Data	IOL Cal	ulation	Toric IO	L Calculatio	on Post Re	fractive IOL	Barrett	Olsen
Surgeon				Measurem	ents			
Surgeon Ger	neric			AL (mm)		Kf (D)		YL (D) -1.45 ax 8°
Target (D)	1 SIA (D)	0.00	n co	ACD (mi 98 LT (mm)		Ks (D) CCT (mm)	42.19 W 0.544	VTW (mm) 11.69
		Universal	II Formula	Universal II	Toric	True K	True K Toric	RX Formula
				True	K 40.93 D Myop	ic Lasik -4.47 D LF	= 2.340	
Model	<u> </u>	Alcon Ac	rysof SN64	Avail	able Toric Lense	s	Right Eye	1
		_		\equiv	Lens	Res Astigm	Right Eye	10
Spherical I	Equivalent Powe	er (D)	23.0				T	
Cylindrical	Power (D)		1.0	00	Non Toric	-0.86 D @ 12°	e m 180-	102
Spherical I	Power (D)		22.5	50	ySof SN6AT2	-0.20 D @ 12°	P let	Add to
Axis of Pla	cement (°)		1(12	ySof SN6AT3	-0.14 D @ 102°	I w	Tas la
Expected r	efraction	+1.12 D -0).20 D @ 1		ySof SN6AT4 IOL ideal t	-0.63 D @ 102° oricity 1.2	5 Sincistar	so Recommended Axis
Barrett True K 1	foric: recommend	led Toric IO	L overview				Back	

Figure 58

As a result, the **"True K Toric"** frame, immediately below, details the best toric lens computed automatically by the system for the manufacturer and model selected previously in the first-step.


From **"Available Toric Lenses"** table you can choose also a different cylinder value for the lens, based on the Residual Astigmatism you want to achieve (under-correction/overcorrection). In particular, the best toric lens

value is shown in the central row and (if available) the ones that under-correct above the central row, the ones that overcorrect below.

At the right side, you can find an image that illustrates the ideal position of the IOL once the implant is in place and the incision location angle.

11.14.5.5 RX Formula (Barrett)

Barrett Rx Formula v1.05, for IOL exchange and piggy back IOLs based on refraction after cataract surgery.

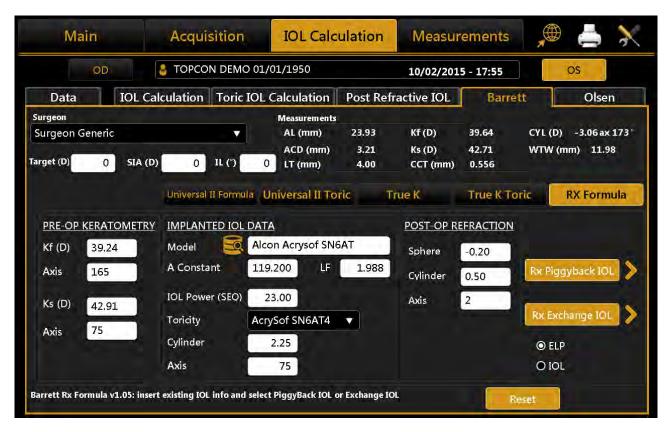

RX Formula (Barrett) is divided into two main steps. The first one consists in inserting all the information regarding the current situation of the patient's eye.

The *"Measurements"* field summarizes the measurement data achieved in the current exam, and they represent the current eye biometry.

In *"Target"* field the target refractive value for the new surgery Post-OP must be inserted. In addition to choosing the **"Target"**, you need to specify also the **"Surgical Induced Astigmatism (SIA)"** and **"Incision Location (IL)"**. The former identify the astigmatism (in diopters) induced by the incision while the latter identify the surgical incision axis.

In the section "*PRE-OP KERATOMETRY*" you can insert the Keratometry data that was measured before the first cataract surgery.

In the section "*IMPLANTED IOL DATA*" you can insert information about the IOL implanted in the first cataract surgery:

- Model: allows to insert the implanted IOL model:
 - selecting from the on-board database of Spheric IOL or Toric IOL lenses, by using the button
 - 3

which opens the selection list

Sph. IOL	Toric IOL
Select manufacturer	
Alcon	- I ,
Select model	
	*


o Inserting manually the description of the implanted IOL



- A Constant/LF: are the calculations constants used for the calculation performed to select the power
 of the IOL implanted in the first cataract surgery, are filled automatically if the lens model is selected
 from the on-board archive or must be inserted manually for manually inserted models.
 Adjust the appropriate Lens Factor/A Constant for IOL by subtracting 0.25 mm from Lens Factor "bag"
 constant if IOL in sulcus.
- IOL Power (SEQ): is the spherical Equivalent power of the Implanted IOL
- **Toricity:** is the cylinder value selection list for the toricity of the implanted IOL, it's automatically populated with a list of cylinder values and submodels if the implanted IOL information has been selected from the on board database



Otherwise if the IOL model information has been inserted manually, this list is locked to "Manual" and the cylinder value must be inserted manually

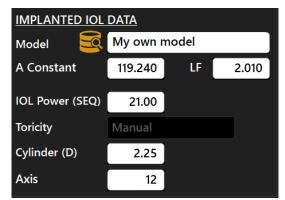

- **Cylinder:** is the cylinder value of the implanted IOL which can be inserted manually or selected from the Toricity list if available
- Axis: is the axis of placement of the implanted toric IOL

In the section "**POST-OP REFRACTION**" you can insert information about the refraction measured after the first cataract surgery during which the Implanted IOL described in the "Implanted IOL data" was applied.

All the mentioned values are required to proceed in the second step of the RX Formula. Then it's possible to proceed to the second step using one of the two available options:

- Rx Piggy Back IOL
- Rx Exchange IOL

Choose between the two options:


- ELP(default) if an error in the predicted ELP is assumed
- IOL for Post Lasik, RK or Low Diopters IOL where predicted ELP is unreliable

The recommended IOL/Piggy Back IOL and Refractive Outcome are calculated according the Barrett Universal II Formula and Barrett Toric Calculator.

ACD (mm) 3.21 Ks (D) 42.71 WTW (mm) 11.98 Target (D) 0 SIA (D) 0 IL (*) 0 LT (mm) 4.00 CCT (mm) 0.556 Universal II Formula Universal II Toric True K True K Toric RX Formula Inversal II Formula Universal II Toric True K True K Toric RX Formula IOL Power Refraction (S.E.) 0.55 D Cylinder Power: IOL Plane 0.50 D - Corneal Plane 1.06 D Predicted refraction (J2 D sph0.14 D @ 92*) Aconst = 119.600, IF = 2.190 Calculated SlA: -1.10 D @ 140* Predicted refraction (S.E.) 0.55 D Predicted refraction (J2 D sph0.014 D @ 92*) Rotate existing IOL (23.00 D g(S.E.) 2.25 D @ 75*) from 75* to 68* Aconst = 119.600, IF = 2.190 Calculated SlA: -1.10 D @ 140* Universal II Formula Universal II Constraitive error: 0.20 D sph. 0.50 D @ 2* Rotate existing IOL (23.00 D g(S.E.) 2.25 D @ 75*) from 75* to 68* Aconst = 10.600 Left Eye Image: Constrained aconstrained	Main	n.		Acqui	sition		IOL Calc	ulation	Measur	ement	s 🔎) 📥 🕅
Surgeon Measurements Surgeon Generic AL (mm) 23.93 Kf (D) 39.64 CYL (D) -3.06 ax 173 Target (D) 0 SIA (D) 0 IL (*) 0 LT (mm) 4.00 CCT (mm) 0.556 Universal II Formula Universal II Toric True K True K Toric RX Formula Nono CSE Optimized constants Optimized constants ACO (SE) 0.500 @ 2* Cylinder Power: 0.20D sph. 0.140 @ 92* Cylinder Power: 0.20D sph. 0.50D @ 2* Optimized constants ACorst = 119.600, IF = 2.190 Calculated SA: -1.10D @ 140* 10L Power Refraction (S.E.) -0.50 (Meniscus) -0.52 D D sh. 0.140 @ 92* Cylinder Power: 0.20D sph. 0.50D @ 2* Corst = 119.600, IF = 2.190 Calculated SA: -1.10D @ 140* 0.50 (Meniscus) -0.29 D D Definition rsi: 7* Octowise for Min. Res. Ast. = -0.38D Corst = 119.600, IF = 2.190 Calculated SA: -1.10D @ 140* 10L Power Res. Cyl. -0.49 D @ 92* -0.15 30 45 60 75 90 105 120 135 150 165 180 0 </th <th>- 40</th> <th>OD</th> <th></th> <th>ТОРСО</th> <th>N DEMO (</th> <th>)1/(</th> <th>)1/1950</th> <th></th> <th>10/02/201</th> <th>5 - 17:55</th> <th></th> <th>OS</th>	- 40	OD		ТОРСО	N DEMO ()1/()1/1950		10/02/201	5 - 17:55		OS
Surgeon Generic Al (mm) 23.93 Kf (D) 39.64 CYL (D) -3.06 ax 173 Farget (D) 0 SIA (D) 0 IL (*) 0 LT (mm) 3.21 Ks (D) 42.71 WTW (mm) 11.98 Intersel II Formula Universal II Formula Universal II Toric True K True K Toric RX Formula IOL Power Refraction (S.E.) 0.55 D Optimized constants Acconst = 119.600, IF = 2.190 IOL Power Refraction (S.E.) 0.55 D Predicted refraction: 0.12 D sph0.14 D @ 92* Cylinder Power: IOL Plane 0.50 D = Correl Optimized constants 0.50 (Meniscus) 0.55 D Predicted refraction: 0.12 D sph0.14 D @ 92* Cylinder Power: IOL Plane 0.50 D = Correl Optimized constants 0.50 (Meniscus) 0.25 D Predicted refraction: 0.12 D sph0.14 D @ 92* Cylinder Power: IOL Plane 0.50 D = Correl Optimized constants 0.50 (Meniscus) -0.29 D Optimized constants Optimized constants Optimized constants 0.50 (Meniscus) -0.29 D Optimized constants Optimized constants 0.50 (Meniscus) -0.29 D Optimized constants Option 100 D Optimized constants <	Data	11	IOL Calc	ulation	Toric IC)L (Calculation	Post Ref	ractive IOL	Bar	rett	Olsen
Target (D) 0 SIA (D) 0 IL (") 0 IT (mm) 3.21 Ks (D) 42.71 WTW (mm) 11.98 Target (D) 0 SIA (D) 0 IL (") 0 IT (mm) 4.00 CCT (mm) 0.556 Universal II Formula Universal II Toric True K True K Toric RX Formula Universal II Formula Universal II Toric True K True K Toric RX Formula $\frac{\text{Rx Piggyback IOL}{0.000 [S.E], 0.50 D @ 2^{2}}$ Cylinder Dower IOL Plane 0.50 – Corneal Plane 1.06 D. Predicted refraction: 0.120 psh0.14 D @ 92^{2} Cylinder Power IOL Plane 0.50 D $= 0.783$ (https://downs.com/doi/10.100 for 10.000 [S.E], 0.250 D @ 2^{2}. Rotate existing refractive error: -0.20D sph. 0.50 D @ 2^{2}. Rotate existing IOL (23.000 LCS; 2.25 D @ 73) from 75' to 68'. Tool 50 (Meniscus) -0.29 D IOL Power Res. Cyl. T-0.0 -0.49 D @ 92^{2} T-0.0 -0.49 D @ 92^{2} T-	Surgeon						Measurements					
Farget (D) 0 SIA (D) 0 IL (") 0 LT (mm) 4.00 CCT (mm) 0.556 Universal II Formula Universal II Toric True K True K Toric RX Formula No 0.00 D (S.E) 0.50 D @ 2' 0.00 D (S.E) 0.50 D @ 2' 0.000 D (S.E) 0.50 D @ 2' 0.50 (Meniscus) 0.55 D Order Average III Toric Total D @ 92' Division refraction: 0.12 D sph0.14 D @ 92' Division refraction: 0.12 D sph0.14 D @ 92' Division refraction: 0.12 D sph0.14 D @ 92' Division refraction: 0.12 P on 75' to 66' Division refraction: 0.12 P on 75' to 166' Division refraction: 0.12 P on 75' to 10' Division refract	Surgeon Ger	neric					AL (mm)	23.93	Kf (D)	39.64	CYL (I	D) -3.06 ax 173°
$\frac{\text{Rx Piggyback 10L}}{0.000 (S.E), 0.50 D @ 2^{2t}}$ $\frac{\text{IOL Power}}{0.55 \text{ D}}$ $\frac{\text{Refraction (S.E.)}}{100 \text{ Meniscus}}$ $\frac{0.55 \text{ D}}{0.50 \text{ (Meniscus)}}$ $\frac{0.55 \text{ D}}{0.50 \text{ (Meniscus)}}$ $\frac{0.29 \text{ D}}{0.50 \text{ (Meniscus)}}$ $\frac{0.29 \text{ D}}{0.50 \text{ (Meniscus)}}$ $\frac{0.29 \text{ D}}{10 \text{ L Power}}$ $\frac{\text{Res. Cyl.}}{1.0 \text{ Cylinder power: 10 \text{ Second 101 C 23.00}}$ $\frac{\text{Res. Cyl.}}{1.0 \text{ Cylinder power: 10 \text{ Second 101 C 23.00}}$ $\frac{\text{Res. Cyl.}}{1.0 \text{ Cylinder power: 10 \text{ Second 101 C 23.00}}$ $\frac{\text{Res. Cyl.}}{1.0 \text{ Cylinder power: 10 \text{ Second 101 C 23.00}}$ $\frac{\text{Res. Cyl.}}{1.0 Cylinder power: 10 \text{ Cylinder p$	arget (D)	0	SIA (D)	0	пс	0	1				WTW	(mm) 11.98
IOL Power Refraction (S.E.) -0.50 (Meniscus) 0.55 D Predicted refraction (S.E.) Predicted refraction (D.2D sph. 6.014 D @ 92° Cylinder Power: IOL Plane 0.50 D - Corneal Plane 1.06 D Predicted refraction (D.2D sph. 6.020 b ph. 6.050 D @ 2° No.50 (Meniscus) 0.55 D Predicted refraction (D.2D sph. 6.014 D @ 92° Rotate existing IOL (23.00 D (S.E.) 2.25 D @ 75°) from 75° to 68° Recom. rot. 7° clockwise for Min. Res. Ast. = -0.38 D Asignatism valDL Axis 0 10L Power Res. Cyl. T-0.0 -0.49 D @ 92° 10.5 0.15 30 45 60 75 90 105 120 135 150 165 180 0 15 30 45 60 75 90 105 120 135 150 165 180				Universal	II Formula	U	niversal II To	ric T	True K	True K	Toric	RX Formula
0.50 (Meniscus) -0.29 D IOL Power Res. Cyl. T-0.0 -0.49 D @ 92° T.0.5 00.411 co.95% T.1.6 0.311 co.95%						E	0.00 D der Power: IOL Pla Predicted refraction xisting refractive e	(S.E.) 0.50 D @ ine 0.50 D ~ Co on: 0.12 D sph. tror: -0.20 D s	2° orneal Plane 1.06 D -0.14 D @ 92° ph. 0.50 D @ 2°	68°	AConst = Calculated	119.600, LF = 2.190 SIA: -1.10 D @ 140°
IOL Power Res. Cyl. T-0.0 -0.49 D @ 92° T.0.5 0.61 D @ 92° T.1.6 0.51 D @ 2° T.1.6 0.51 D @ 2°	0.50 (Meni	scus)	-0.2	29 D	0.		Astig	matism vs IOL Axis		100		T e m
T 1 0 0.31 0 @ 92 -3 -3 -1 -3	IOL Pow	rer	Res.	Cyl.	-0.6- sy -1.2-	_		1	·		190-	2° 0 0
T 1 0 0.31 0 @ 92 -3 -3 -1 -3	T-0.0	t.	-0.49 E) @ 92°	enques	1	/					· ·
1-1-0 -0-221-01-02 Picture Pic	1.0.5		0.14 (-3+	1	5 30 45 60		5 120 135 150	165 180		
			-0.21	0 @ 2	1			102.0010	_	_	Indision Re	commended Axis Current Axis

The recommended TORIC IOL and Axis alignment for the targeted refractive outcome is displayed. The axis that provides the minimum astigmatism for the existing IOL is calculated as well as the rotation in degrees from the current axis of the existing implanted IOL.

The SIA and Optimised Lens Factor/A Constant are provided according to the pre and post op Keratometry and the refractive outcome.

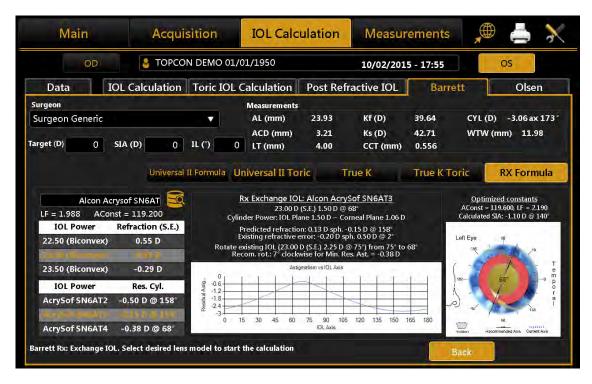

11.14.6 Olsen Calculator

Olsen Calculator is an optional module: ask your distributor for details.

The Dr. Olsen's formula is based on standard paraxial ray tracing using standard keratometry of the anterior surface of the cornea and optical biometry of the eye. The IOL power calculation methodology using ray tracing and improved algorithms to predict the position of the IOL inside the eye is intellectual property protected by international law. It includes a toric calculator which is based on a proprietary method to account for the non-measured posterior surface of the cornea.

11.14.6.1 Olsen Spherical IOL (Dr. Olsen formula)

Olsen Spherical IOL calculator performs calculations of spherical IOLs using the Dr. Olsen's formula.

Mai	n.	Acqui	sition	IOL Calc	ulation	Measur	ements		
	OD .	👌 торсо	N DEMO 01/	01/1950		10/02/201	5 - 17:55	OS	
Data	IOL C	alculation	Toric IOL	Calculation	Post Ref	ractive IOL	Barre	tt	Olsen
Surgeon				Measurements					
Surgeon Ge	neric			AL (mm)	23.93	Kf (D)	39.64	CYL (D) -	3.06 ax 173
				ACD (mm)	3.21	Ks (D)	42.71	WTW (mm)	11.98
Target (D)	0			LT (mm)	4.00	CCT (mm)	0.556	PUP Ø (mm)	4.45
Oculentis	Y	Alcon		Поуа	M	AMO		AIMO	
Oculentis		Alcon	•	Ноуа	•	AMO	<u> </u>	AMO	
L-312	▼	SN60WF	•	iMics Y-60	DH 🔹	Tecnis 1 Z	СВ00 🔻	Tecnis ZA	9003 🔻
IOL @ Target	ACD = 4.20	IOL @ Target	ACD = 4.67	IOL @ Target	ACD = 4.70	IOL @ Target	ACD = 4.92	10L @ Target	ACD = 4.92
21.91	A = 118.5	22.73	A = 119.0	22.15	A = 118.8	23.12	A = 119.4	23.01	A = 119.3
IOL (D)	REF (D)								
21.00	0.67	21.50	0.88	21.00	0.83	22.00	0.80	22.00	0.73
21.50	0.30	22.00	0.53	21.50	0.47	22.50	0.45	22.50	0.37
22.00	-0.07	22.50	0.16	22.00	0.11	23.00	0.09	23.00	0.01
22.50	-0.44	23.00	-0.20	22.50	-0.26	23.50	-0.28	23.50	-0.36
	-0.82	23.50	-0.57	23.00	-0.63	24.00	-0.64	24.00	-0.73
23.00									

In the *"Surgeon"* field, you can choose which surgeon will perform the IOL implant and any customization of the constants or presetting of the preferred lenses will be applied on this basis.

In *"Target"* field the target refractive value for the Post-Op must be inserted.

The "Measurements" field summarizes the measurement data.

From the drop-down menu, select the IOL manufacturer and model. The available lenses you can select come from a list of models whose calculation constants and optical parameters have been validated. The user can in case insert new manufacturers and/or models inside IOL settings section (see 13.4.3).

Once this data has been entered, the most appropriate lens can be chosen at the discretion of the surgeon. The latter is highlighted in orange. Once selected, the lens will be memorized as the preferred one and will be shown highlighted on the report printout.

Pressing "Reset" will reset the initial preset conditions.

11.14.6.2 Olsen Toric IOL (Dr. Olsen formula)

Olsen Toric IOL Calculator (based on Dr. Olsen's formula) is divided into two main steps. The first one consists on the calculation of the Spherical Equivalent Power; in the second one you can select the toric IOL that produce the best correction. The first-step interface has quite the same structure as the spherical IOL calcualtion. The available toric lenses you can select come from a list of models whose calculation constants and optical parameters have been validated. The user can in case insert new toric manufacturers and/or models inside toric IOL settings section (see 13.4.3).

Mair		Acquis	sition		IOL Calcu	ulation	Measur	ements		
	OD I	👌 торсоі	N DEMO	01/	01/1950		10/02/201	5 - 17:55	0	s
Data	IOL Ca	alculation	Toric IQ	OL (Calculation	Post Ref	ractive IOL	Barre	tt	Olsen
Surgeon					Measurements					
Surgeon Ger	heric		-		AL (mm)	23.93	Kf (D)	39.64	CYL (D)	-3.06 ax 173°
Target (D)	0 SIA (D)	0	IL (°)	0	ACD (mm) LT (mm)	3.21 4.00	Ks (D) CCT (mm)	42.71 0.556	WTW (mm PUP Ø (mm	
					Spherical IO		oric IOL		II:	
Oculentis		Alcon								T
LS-313 Tx		Acrysof S	N6AT	•		•				
	ACD = 4.07	IOL @ Target	ACD = 4		IOL @ Target		IOL @ Target		10L @ Target	
10L @ Target 21.64	ACD = 4.07 A = 118.2	22.75	A = 11	9.0						
			A = 11 REF (D		IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
21.64 IOL (D) 20.50	A = 118.2 REF (D) 0.85	22.75 IOL (D) 22.00	REF (D 0.55))	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
21.64 IOL (D) 20.50 21.00	A = 118.2 REF (D) 0.85 0.48	22.75 IOL (D) 22.00 22.50	REF (D 0.55 0.18))	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
21.64 IOL (D) 20.50 21.00 21.50	A = 118.2 REF (D) 0.85 0.48 0.11	22.75 IOL (D) 22.00 22.50 23.00	REF (D 0.55 0.18 -0.18))	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)
21.64 IOL (D) 20.50 21.00	A = 118.2 REF (D) 0.85 0.48	22.75 IOL (D) 22.00 22.50	REF (D 0.55 0.18)	IOL (D)	REF (D)	IOL (D)	REF (D)	IOL (D)	REF (D)

In addition to choosing the **"Target"**, you need to specify also the **"Surgical Induced Astigmatism (SIA)"** and **"Incision Location (IL)"**. The former identify the astigmatism (in diopters) induced by the incision while the latter identify the surgical incision axis.

After having selected the toric IOL model, a values table from which the **Spherical Equivalent Power** is obtained. Once you choose a lens, pressing **"Next"** at the bottom right, you enter in the second-step of toric IOL calculation.

Main		Acqui	sition	IOL Calcu	ulation	Measur	ements	🔎 📥 À
0	D	👌 торсо	N DEMO 01	/01/1950		10/02/201	5 - 17:55	OS
Data	IOL C	alculation	Toric IOL	. Calculation	Post Ref	fractive IOL	Barre	tt Olsen
Surgeon		1.000		Measurements				
Surgeon Gene	eric			AL (mm)	23.93	Kf (D)	39.64	CYL (D) -3.06 ax 173°
and the second second				ACD (mm)	3.21	Ks (D)	42.71	WTW (mm) 11.98
Target (D)	0 SIA (D) 0	IL (°)	LT (mm)	4.00	CCT (mm)	0.556	PUP Ø (mm) 4.45
				Available	and a familie		OS	105 90 75
Model Spherical Equ	uivalent <u>Por</u>		is LS-313 T3 21.50	Ler	15	Res Astigm	OS 135	120 105 90 75 60 45 30
Spherical Equ			21.50	Ler LS-31	15 3 T1 -	1.40 D @ 171°	135 150 165	120 105 75 60 45
Spherical Equ Cylindrical Po	ower (D)		21.50 3.00	Ler LS-31 LS-31	15 3 T1 - 3 T2 -	A PA	135	120 105 75 60 45 30
Spherical Equ Cylindrical Po Spherical Por	ower (D) wer (D)		21.50 3.00 20.00	Ler LS-31 LS-31 LS-31 LS-31	15 3 T1 - 3 T2 - <mark>3 T3 -</mark>	1.40 D @ 171° 0.85 D @ 171°	135 150 165 180 -	
Spherical Equ Cylindrical Po	ower (D) wer (D)		21.50 3.00	Ler LS-31 LS-31 LS-31 LS-31	15 3 T1 - 3 T2 - <mark>3 T3 -</mark> 3 T4	1.40 D @ 171° 0.85 D @ 171° 0.30 D @ 169°	135 150 165	
Spherical Equ Cylindrical Po Spherical Por	ower (D) wer (D) ment (°)	wer (D)	21.50 3.00 20.00	Ler LS-31 LS-31 LS-31 LS-31 LS-31 LS-31	15 3 T1 - 3 T2 - <mark>3 T3 -</mark> 3 T4	1.40 D @ 171° 0.85 D @ 171° 0.30 D @ 169° -0.25 D @ 85°	135 150 165 180 -	120 105 75 60 45 30 15
Spherical Equ Cylindrical Po Spherical Por Axis of Place	ower (D) wer (D) ment (°)	wer (D)	21.50 3.00 20.00 82	Ler LS-31 LS-31 LS-31 LS-31 LS-31 LS-31	15 3 T1 - 3 T2 - <mark>3 T3 -</mark> 3 T4	1.40 D @ 171° 0.85 D @ 171° 0.30 D @ 169° -0.25 D @ 85°	135 150 165 180 -	120 105 175 60 45 30 15 0 15
Spherical Equ Cylindrical Po Spherical Por Axis of Place	ower (D) wer (D) ment (°)	wer (D)	21.50 3.00 20.00 82	Ler LS-31 LS-31 LS-31 LS-31 LS-31 LS-31	15 3 T1 - 3 T2 - <mark>3 T3 -</mark> 3 T4	1.40 D @ 171° 0.85 D @ 171° 0.30 D @ 169° -0.25 D @ 85°	135 165 180- Na <u>59</u>	120 105 175 60 45 30 15 0 15

As a result, the **"Toric Calculator"** frame, immediately below, details the best toric lens computed automatically by the system for the manufacturer and model selected previously in the first-step.

From **"Available Toric Lenses"** table you can choose also a different cylinder value for the lens, based on the Residual Astigmatism you want to achieve (under-correction/overcorrection). In particular, the best toric lens value is shown in the central row and (if available) the ones that under-correct above the central row, the ones that overcorrect below.

At the right side, you can find an image that illustrates the ideal position of the IOL once the implant is in place and the incision location angle.


11.15 Data saving

After performing some acquisitions and eventual IOL calculations, in order to save the data from the examination, click on the "Main" button. As shown in Figure 59, the software will ask the user to confirm the action.



Figure 59

If you want to save the data of the patients without doing an acquisition, you have to press the "main" button when you are in the "acquisition" panel. The system will ask you if you want to save the current patient's details.

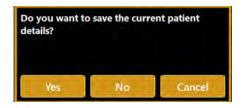

11.16 RX/AL Trend Function

If enabled from the settings (see section 13.6) the RX/AL Trend Function can be accessed from the patient list view.

The RX/AL Trend function allows to review the progression of biometry value of the selected patient over time, based on the data present in the local archive with the possibility to add further data manually.

New	Searc	h	Acquisition		12/07/2019 12:31
Last Name		<u>PATIENT</u>		EXAM LIST	
 ☑ Last Name ◯ ID 		Last Na Name	ΤΟΡϹΟΝ	10/02/2015 17:55:00	New Exam
TOPCON DEMO		Date of	DEMO Birth 01/01/1950		Open
		ID Gender			Delete
		O Male	e O Female	RX/AL Trend	Post Op
		1	2 3 4	5 6 7 8	9 0 ESC
		q	W E	RTYU	ΙΟΡ
		Α	S D F	GHJ	K L →
Edit	Page 1 / 1 Delete	Ζ	x c v	B N M '	←

By pressing the RX/AL Trend button you access the two main actions for the selected patient:

- Data Review, access to the RX/AL Trend data overview section
- New Exam, proceed with a new examination



11.16.1 RX/AL Trend Action Window

The RX/AL Trend Action Window provides access to the two main actions if a patient is selected from the list.

If no patient is selected the two actions are disabled.

The windows provides also access to the RX/AL Trend function specific options. Other options such as Cylinder notation or Display unit depend on the general settings of the Aladdin application.

11.16.2 New Exam Action

You can enter the current refraction measurement data before performing the biometry acquisitions. The data is not compulsory to proceed with acquisition, you can update them afterward in Data Review section.

Vew Exam 👗 Patient M	Ayopia 11/11/2011	05/0	5/2019
	Refraction Da	ata	
 Non-Cycloplegic Cycloplegic 	Objective Subjective	Refraction Over Refraction	
1		and a	
Refraction Sphere:	D	Refraction Sphere:	N
Cylinder:	D	Cylinder:	
Axis:	0	Axis:	
Vertex Distance: 14	.00 mm	Vertex Distance: 14.00 m	nm
		Go To Ac	quisition

If you select Over Refraction option you are prompted also the input of the Contact Lens data.

New Exam 👗 Patient	Myopia 11/11/2011		05/05/2019	X
	Refraction	Data		
○ Non-Cycloplegic ○ Cycloplegic	Objective Subjective	O Refracti O Over Re		
R				1
Refraction		Refrac	tion	
Sphere:	D	Sphere:	D	
Cylinder:	D	Cylinder:	D	
Axis:	•	Axis:	•	
Vertex Distance:	4.00 mm	Vertex Distance:	14.00 mm	
Contact Ler	is	Contact	Lens	
Base Curvature:	D	Base Curvature:	D	
Lens Power:	D	Lens Power:	D	
			Go To Acquisitio	


When pressing on Go to Acquisition a new examination session will be activated and Aladdin will be ready for acquisition.

Using this function it is possible to perform the following acquisition flow:

- 1. Measure Subjective/Objective Refraction externally
- 2. Input Refraction Data in the form
- 3. Perform Topography and Biometry acquisition on Aladdin

11.16.3 Data Review

The data present in the Data Review section is provided by examinations directly performed with Aladdin and stored onboard together with manual data eventually imported or added.

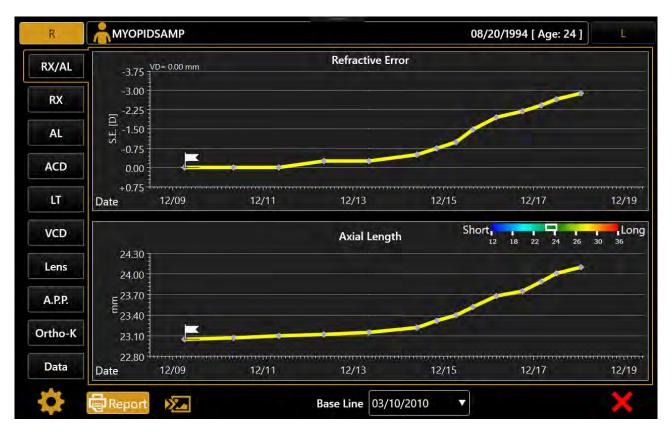

BASELINE

Starting condition of the trend profiles. Can be changed on the fly to any examination date.


MONTH BASE

Windowing period for Variation profile calculation.

NOTES

Notes can be added/edited in relation to each single examination date and reviewed.


ALADDIN - User manual

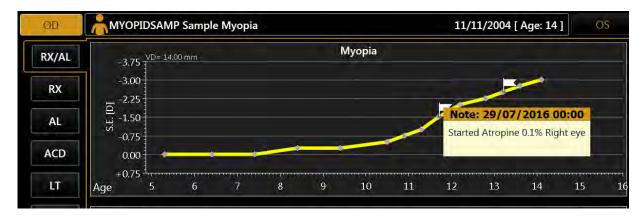





RX/AL: Refraction & Axial Length trend: Combined trend view with Refraction error S.E. and Axial Length.

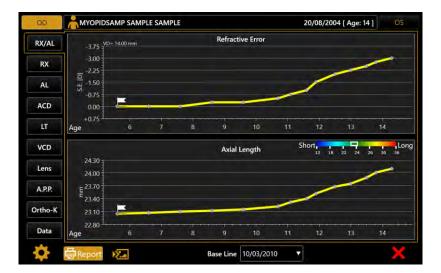

RX: Refraction trend: Refraction S.E. absolute value (at reference Vertex Distance) trend and yearly variation (on 3,6,12 months basis).



AL: Axial Length trend: Absolute Value trend and yearly variation (on 3,6,12 months basis).

OD	MYOPIDS	AMP SAMP	LE SAMPL	E				20/08/200	04 [Age:]	14]	OS
RX/AL					Axial I	ength	SI	nort <mark>.</mark> 12 18	22 24	26 30	Long
RX	24.30 24.00									/	•
AL	E ^{23.70} 23.40							/	-		
ACD	23.10	E	•	-	+	-	-				
LT	22.80 ¹ Age	6	7	8	9	10	11	12	13	14	, , , ,
VCD	+0.40 1			V	ariation AL	/ Year		On a	3	▼ Mont	ths Bas
Lens	+0.30							m		7	
A.P.P.	€+0.20						~	-	\mathbb{N}	-/	•
Ortho-K	+0.10	لتفعر				100	/				
	0.00 1	6	7	-8	9	10	11	12	13	14	

Other Trends:

For other biometry values the absolute value trend and the yearly variation are shown (if the data is available):

- ACD, Anterior Chamber Depth
- LT, Lens Thickness
- VCD, Vitreous Chamber Depth
- APP, Average Pupillar Power (4.5mm diameter)
- Lens Power: estimated from biometry data and refraction data if available

Ortho-K section

View to compare Topography map between two different examinations. Available only for examination performed directly with Aladdin. Compare variations of Kerato-refractive indexes. Useful to verify effects of Ortho-K application.



KERATO REFRACTIVE INDEXES

- 3mm and 5mm Keratometry
- Average Pupillar Power (4.5mm diameter pupil size)
- Asphericity
- Spherical Aberration (Longitudinal Spherical Aberration LSA)
- Curvature Irregularity (SD, Corneal Power std. deviation)
- Asymmetry (SAI index)

DATA

R		SAMP				08	/20/1994 [Age	e: 24 j
RX/AL				R	aw Data		+Add D	ata 🚺 Impor
	DATE	RX SE [D]	AL [mm]	ACD [mm]	LT [mm]	Lens Pow. [D]	A.P.P. [D]	
RX	01/31/2019							1 前
AL	12/20/2018	-3.00	24.10	3.50	4.36	23.74	44.35	1 D
	09/07/2018							1 T
ACD	06/05/2018	-2.75	24.01	3.48	4.36	23.65	44.40	1 d
LT	04/04/2018							🕥 🗇
-	02/03/2018	-2.50	23.89	3.50	4.33	23.98	44.29	S 🗇
VCD	12/02/2017							1 D
- Cord	09/03/2017	-2.25	23.75	3.48	4.34	24.07	44.37	S 🗇
Lens	06/04/2017							8
A.P.P.	02/03/2017	-2.00	23.68	3.50	4.33	23.85	44.50	S 🗊
	11/04/2016							1 🖸
Ortho-K	07/29/2016	-1.50	23.52	3.48	4.32	23.85	44.45	🕥 🛅
Data	05/22/2016							10 1
	03/15/2016	-1 00	23 40	3 47	4 35	23.81	44 35	A m

Raw data table review and management (edit, delete).

Import Button: Import data from template sheet for the current patient only.

Add Button: Add data manually by filling the form.

ROW Buttons

Edit Button: Edit data of the selected row, shows edit form with current values



Delete Button: Delete the selected row

REPORT & SCREENSHOT



REPORT: Create Report with:

- o Refraction and Axial Length trend
- o Data table
- o Notes



SCREENSHOT: create a report page with the current view

Available output destinations:

- Printer (as in the main printing form of the Aladdin application)
- o USB drive export
- o Network Shared folder

11.16.4 RX/AL Trend Option


By pressing on the settings icon you get access to the RX/AL Trend related options.

<u>Time Axis opt</u>	ions			
Date	Age	Fixed Age range:	5	16
<u>Refraction</u>				
Reference Vertex	Distance [mm]	14.00		
<u>Utils</u>				

Time Axis Options

The charts shown in the Data Review section can have the time axis visualized by:

• Date, calendar date of examination data

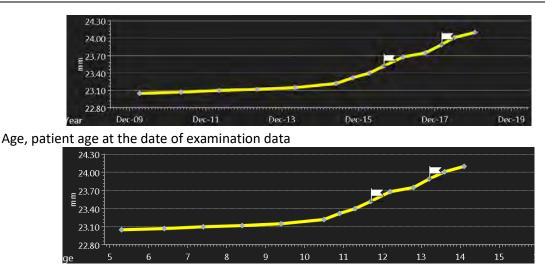

If the Age option is selected for the Time Axis, the Fixed age range can be enabled/disabled and configured.

If the fixed age range is enabled at the least the age range defined by the two related values will be visualized in the time axis. Data eventually exceeding such range will be shown as well.

Refraction Options

The reference vertex distance value will determine the vertex distance (VD) to which all the refraction data will be shown in the charts (transposed if input VD was different).

Utils

• **Import data from** predefined sheet template file. A dialog to select a file from connected USB pen drive will be shown.

In this action can be imported data for:

- The currently active patient (will be matched by ID)
- o Another patient already existing in Aladdin archive (will be matched by ID)
- A patient not present in the Aladdin archive (a new patient will be created if at least ID and date of Birth are specified in the sheet, otherwise the data will be ignored)
- **Export**: create on the connected USB drives an empty sheet template file, to be filled externally and used for importation of previous/already existing data.

DATA IMPORT TEMPLATE SHEET

IMPORTANT: Do not alter the structure of the sheet!



This sheet is to be used for importing automatically a set of data into the RX/AL Trend data collection available for Aladdin.

The sheet can be arranged to contain data for different patients and both Right and Left eye.

Each row of the can contain data for:

- Patient identification
- Eye examination data identification
- Biometry data
- Refraction Data
- Note belonging to the examination date

REQUIRED DATA (for each row)

- Patient ID: patients already present in the Aladdin archive will be matched by this field
- **Patient Date of Birth**: required to import patients which are not already existing in the Aladdin archive
 - (if the DoB is not specified the patient will not be created and the relative data ignored)
- Examination Eye (Dx or Sx only accepted values)
- Examination Date

IMPORTATION RULES

- If the compulsory data is not provided the single sheet row will be ignored
- Patients not already present in the archive
- The Dates (patient DoB and Exam Date) will be parsed according to the *Date Format* selected in the settings of the Aladdin application. Not compatible dates will cause the row to be ignored.
- If in relation to a single exam date/patient ID there are more than one row (typically two for each Exam Date of a patient for Right and Left eye):
 - o Only the first occurrence for each eye will be imported
 - o Only the note (if present) of the first occurrence will be imported

EXAMPLE

Data for one Examination date of Patient with date for both Right and Left eye

	Patient			Exer	noiteni			Refraction	/OverRefraction			Over Sefraction Info		Sice	etry-		Note
ID + LastNam	r Name v	DoB +	Gender +	Eyer	•	Rx VD [mm]	Re Sphere (D)	Re Cyl (D)	Ra Cyl Assi [deg]	Ra type	Eye Status	Lens Power SE (D) Base Curre Inm)	ALImmi	ACD time1	LT (mml		Test note
Panashti		AUXPART		2	10/07/19	12.00	2.62	0.00	0.00	Subjective	Normal	+2.00	25.67	1.00	3.46	A4.37	Note Day Fact Tech Note für anne Prisant Erne das & imported
Pulseit		are overall		H	10/07/13	0.00	4.37	-6.25	172.00	Claircovy	Cucloplegie						ant imported
				0.0	414114	0.10			Net.	delighter							

DATA DESCRIPTION

Patient Identification

	Patient		
Surname 💌	Name 🔻	DoB 🔻	Gender 💌
Surname	Name	11/11/2011	м
Surname	Name	11/11/2011	м
	Surname	Surname Vame Surname Name	Surname Name DoB Surname Name 11/11/2011

• Eye Examination identification Eye identification and Exam date

Exar	Examination										
Eye 🔻	Date 🔻										
Dx	07/05/2019										
Sx	07/05/2019										

Refraction Measurement Data

		Refraction/C	verRefraction			Contact	t Lens
Rx VD	Rx Sphere	Rx Cyl	Rx Cyl Axis	Rx type	Eye Status	Lens Power SE	Base Curve
12	-2.22	0	0				
12	-2.80	0	0				

- o Vertex Distance in mm relative to refraction measurement
- Refraction measurement in preferred notation (positive cylinder notation or negative cylinder notation)
- o Rx type:
 - Subjective or Objective
- Eye status: condition at the moment of refraction measurement
 - Cycloplegia or Normal
- o Contact Lens data (in case of Over Refaction measurement)
 - Lens Power SE in Diopters
 - Base Curve of contact lens in mm

Biometry Data

	Bior	metry	
AL	ACD	LT	APP
25.02	3.6	4.02	45.2
25.15	3.8	4.2	44.7

- Axial Length [mm]
- Anterior Chamber Depth (Epithelium to Anterior Capsule) [mm]
- o Crystalline Lens Thickness [mm]
- APP (Average Pupillar Power) [D]: Corneal Power in the entrance pupil optical zone

12 MEASUREMENTS


All measurements performed during the examination can be reviewed in detail in the "Measurements" section.

There are four types of measurement.

- KER: Keratometry
 - o ZER: Zernike Analysis
- AL: Axial length
- ANT: Anterior Segment sections: CCT, ACD, LT
- PUP: Pupillometry

to which various environments correspond, described in detail in the following sections.

12.1 Topographic map (KER)

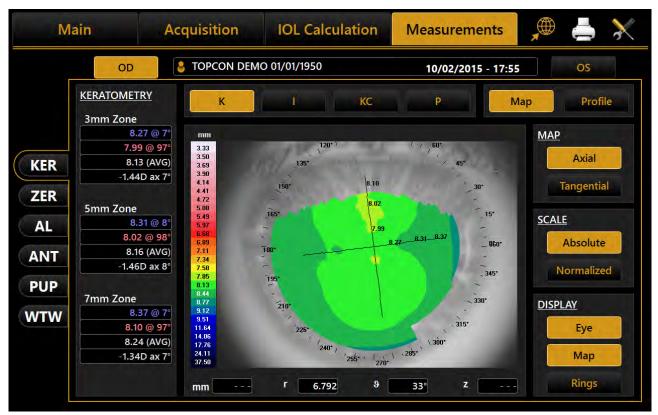



The environment displayed is shown in Figure 60.

Click on the **"OD"** or **"OS"** buttons to display the map of the right or left eye. The R and L buttons are only active if the keratometry of the eye in question has been acquired.

In the right column, you can select the following options:

- Axial or Tangential: axial map or tangential map
- Absolute or Normalized: absolute scale or standardized scale
- Eye, Map, Rings: to display the image of the eye, the map, the rings

Pressing on any point on the map displays the following information:

- Diopters (D)
- Radius (r)
- Meridians (θ)
- Altimetry (z)

The Scale buttons allow to switch between Absolute and Normalized (adjustable) scale color steps. When Normalized is pressed the button is replaced with controls that allows to adjust the color step for the current topography map. Minimum step size is 0.25 D or 0.05 mm depending on the selected measure unit.



Refer to section 13.2 for more settings relative to the topographic map representation. The buttons at the screen top show the topographic map indices (see following paragraphs for details):

- K: Keratometry
- I: Keratorefractive indices
- KC: Keratoconus
- **P**: Pupil

•

12.1.1 Keratometry

Press the *"K"* button to display the keratometric data on the 3 mm, 5 mm and 7 mm zones, as shown in Figure 60 (by settings the zones can be set to 2,4,6 mm).

12.1.2 Keratorefractive indices

Press the "I" button to view the keratorefractive indices:

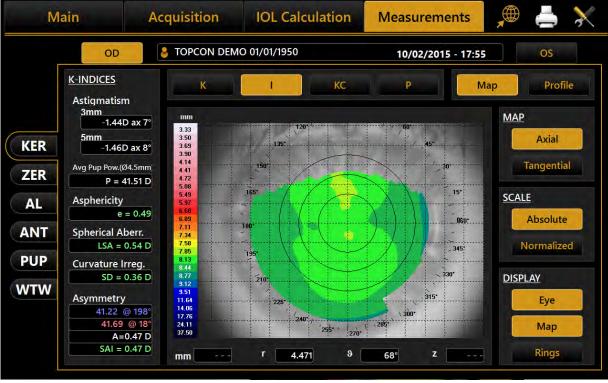


Figure 61

- Astigmatism: Astigmatism at 3 and 5 mm (or 2 and 4 mm)
- Pupil Avg: Average pupil power for a pupil of 4.5 mm
- Asphericity: Asphericity of the cornea at 8 mm diameter
- Spherical Aberration: Longitudinal spherical aberration of a 4.5mm diameter cornea area

- **Curvature Irregularity:** Irregularity of curvature calculated on the standard deviation of the instantaneous readings for a 4.5mm diameter cornea area
- Asymmetry + SAI: Asymmetry between the most curved hemisphere and the flattest one calculated for 4.5mm diameter cornea area and an SAI (Surface Asymmetry Index) that represents the surface asymmetry index of the 4.5mm diameter cornea area.

12.1.3 Keratoconus

Press the *"KC"* button to open Apex keratometry with the following information:

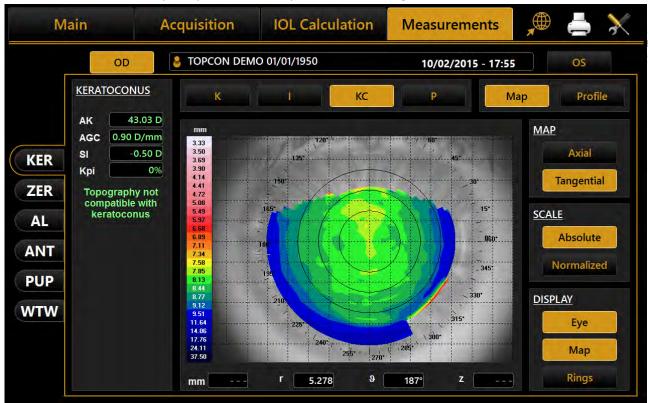


Figure 62

- **AK**: Apical curvature. Represents the power of the cornea in its apex
- AGC: apical gradient of curvature. Represents the corneal power average variations (per unit of length), taking the apical power as reference.
- **SI**: difference between the average power of two circular zones centered in the vertical axis of the ruler and placed in the lower hemisphere and in the upper hemisphere of the cornea respectively.
- *Kpi*: Keratoconus diagnosis probability index.

Based on the combined evaluation of the first three indices with the probability index, there are three different possibilities: topographic picture not compatible with keratoconus (green); suspected keratoconus (yellow); topographic picture compatible with keratoconus (red).

If the topographic picture is compatible with keratoconus or indicates a suspected keratoconus, the numerical values of the geometric parameters of the cone are shown at the bottom of the panel. These are:

- A: area of the keratoconus (mm²)
- *D*: average diameter of the keratoconus (mm)
- *r*, Ø: polar coordinates (mm, °) of the barycentre of the keratoconus in relation to the centre of the *map*
- *RND: circularity factor of the keratoconus*

12.1.4 Pupil

Press the "P" button to open the pupil's indices:

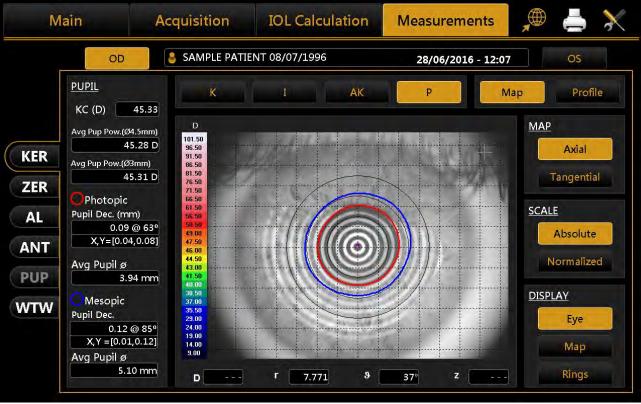


Figure 63

- KC: KC represents the central keratometry in diopters
- Avg Pupil Power: Average pupil power for a pupil of 4.5 mm and 3.0 mm
- Photopic(red outline)
 - **Pupil Dec.:** Pupil decentration in polar and cartesian coordinates from the cornea vertex
 - Avg Pupil Ø: Mean diameter of the pupil
- Mesopic(blue outline)
 - **Pup Dec.**: Pupil decentration in polar and cartesian coordinates from the cornea vertex
 - Avg Pupil Ø: Mean diameter of the pupil

12.1.5 Profile

Press the *"Profile"* button to view the curvature profile along the most curved meridian and the flattest meridian (red and blue).

The difference is displayed in green (Figure 64).

By pressing the arrow buttons, you can vary the flattest and the most curved meridians.

The graph will be modified accordingly.

Pressing the "Map" button, you go back to the topographic map.

	OD SAMPLE PATIENT 01/01/2001 04/10/2013											- 15:05 OS		
199'	1		289°	r.						109	ō	19°	Map. Prof	
-6.0		5.0	-4.0	-3.0	-2.0	-1.0	47.0	1.0	2.0	3.0	4.0	5.0		
							46.0							
			/				45.0	-		-				
		/	-				44.0			/			· ALENYL	
							43.0							
							42.0							
							41.0							
							40.0						218 ⁴	
							39.0 38.0						and the second second	
3							1 30.0						Steepest Meridian	
dK													109	
-6.0		5.0	-4.0	3.0				1.0	20_	3.0	4.0	5.0		
				÷			2.0				· · · · · · ·			
							4.0						Flattest Meridian	
			<u>.</u>	<u>ģ</u>			6.0							
						_	8,0			_			19	

Figure 64

12.2 Zernike

The Zernike module provides a comprehensive view of the wave front aberrations generated by the front surface of the cornea. The results of the Zernike axis are illustrated by means of numerical indices and graphic representations (Figure 65).

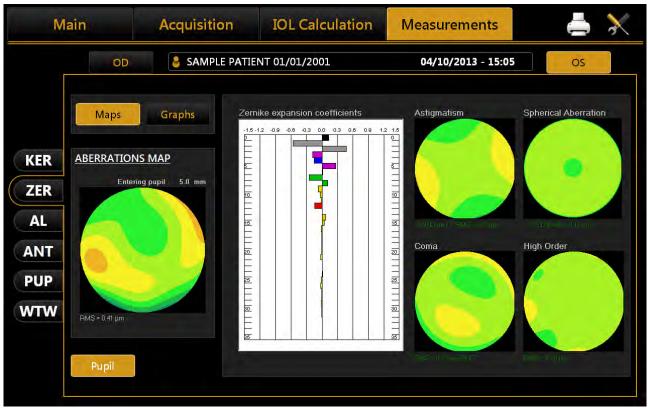


Figure 65

Click on the **"OD**" or **"OS**" buttons to view the results of the Zernike analysis for the right or left eye.

On the left the <u>Aberrations Map</u> is detailed, representing the total aberration that corresponds to the sum of all the aberration components and the RMS value. This allows you to quantify the deviation with respect to an ideal wavefront.

On entering the module, the aberrations map is displayed ("Maps" section):

- <u>Histograms of the Zernike expansion coefficients</u>: each histogram represents the weight of the corresponding polynomial.
- <u>Primary aberrations map:</u>
 - ✓ Astigmatism: the map, the magnitude in diopters, the axis and the RMS value are displayed
 - ✓ Spherical aberration: the map, the quantity of longitudinal spherical aberration in diopters and the RMS value are displayed
 - ✓ Coma: the map, the RMS value and the direction are displayed
 - ✓ High Order: all the components of a higher order than the primaries are grouped; the map and the RMS value are displayed.

Click on "Graphs" at the top left to display the vision quality summary (Figure 66).

This section displays:

- Zernike Coefficient pyramid: represents the numerical value of each coefficient by means of a grey scale; the greater the coefficient, the greater the color contrasts with the pyramid's background.
- **Point Spread Function**: represents the intensity of the wave front in the retina.

- **Spot Diagram**: represents the spatial distribution of the wave front over the retina.
- Visus/Visus Low Contrast: represent the patient's real vision at high and low contrast.

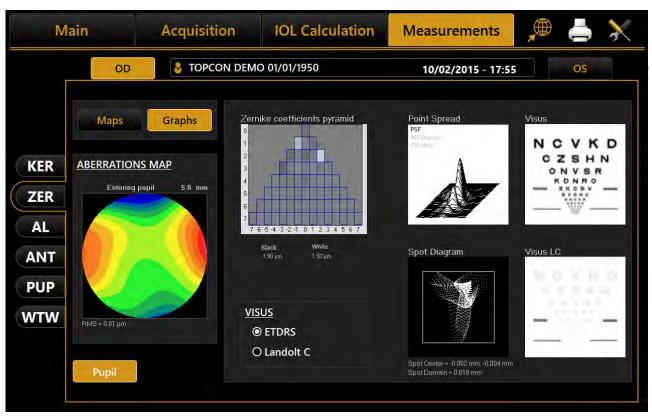



Figure 66

The data displayed refers only to the component induced by the anterior surface of the cornea, not by the eye's entire optical system.

Press the *"Maps"* button to return to the maps display.

The *"Pupil"* button opens a panel (Figure 67) where you can select the diameter of the pupil (in a range between 2 mm and 7.5 mm) to see how the aberrations change with the variation of the pupil diameter.


Double click on each images (except Visus simulation) to review the chart enlarged.

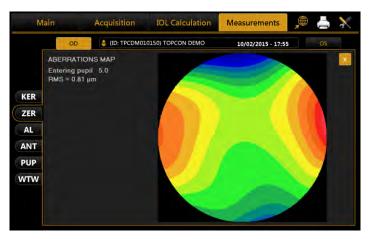

Main		Acquisition	IOL Calculation	Measurements	🔎 📥 💥
	R	S TOPCON E	DEMO 01/01/1950	10/02/2015 - 17:5	5
	Maps	Graphs	Zernike coefficients pyramid	Point Spread PSF PSF Dense	Visus
KER	ABERRATION	<u>S MAP</u>			NCVKD CZSHN onvsr
ZER	7.0 mm	7.5 mm	5 6 7		ZKCSV DVOHO HZOKO
AL	6.0 mm	6.5 mm 5.5 mm	-7 -6 -5 -4 -3 -2 -1 0 1 2 3 4 5 6 7 Black White 1 50 μm 1 50 μm	Spot Diagram	Visus LC
PUP	4.0 mm	4.5 mm			10 C = 10 O 5 = 10 H
WTW	3.0 mm	3.5 mm	<u>VISUS</u> © ETDRS		
	2.0 mm	2.5 mm	O Landolt C	Spot Center = 0.000 mm, -0.002 mm Spot Domain = 0.009 mm	

Figure 67

It is possible to switch between ETDRS and Landolt C Visus simulation view.

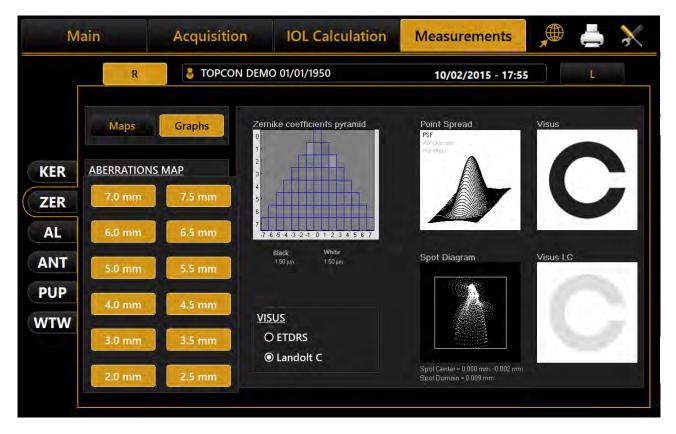

12.3 Axial Lenght (AL)

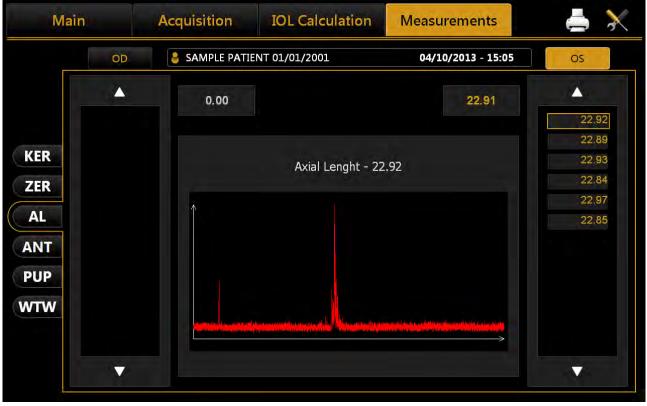


Figure 68

Figure 68 shows an axial length measurement.

In this screen you can select and display the interferometric graph for each measurement and from the left and right columns the measurements performed for the right eye and left eye, respectively. The measurements highlighted in yellow are the ones used to calculate the average axial length and are acceptable with respect to signal/noise. Those highlighted in red are those discarded by the system, for being unacceptable. It is always advisable to repeat a discarded measurement carefully.

Main	Acquisition	IOL Calculation	Measurements	\Rightarrow
O	D 8 TOPCON DEM	10 01/01/1950	10/02/2015 - 17:55	OS
Scan #1	OD CCT 0.544 mm ACD 3.14 mm LT 4.04 mm		OS 0.556 mm 3.21 mm 4.00 mm	CCT ACD LT
ZER AL ANT PUP WTW	OD Scan # 1 ACD 3.14 mm LT 4.04 mm			

12.4 Anterior Segment Sections (ANT, CCT-ACD-LT)



Figure 69 gives an example of ANT (Anterior Segment Section) measurements.

In this screen you can select and display the interferometric graph for each measurement and from the left and right columns the measurements performed for the right eye and left eye, respectively.

As for the axial length, if the instrument does not record good quality signal, or if the data is inconsistent, the acquisition is discarded.


On the top-central section the over-all values for both eyes are reported, while in the column at the left side of the interferometric graph the results for the selected acquisition are reported.



LT results are always with warning signs if the AL measurement for the same eye is not present.

12.5 Pupillometry (PUP)

The pupillometry module allows displaying and analyzing the dynamic and static pupillometry (pupil images acquired in controlled light conditions).

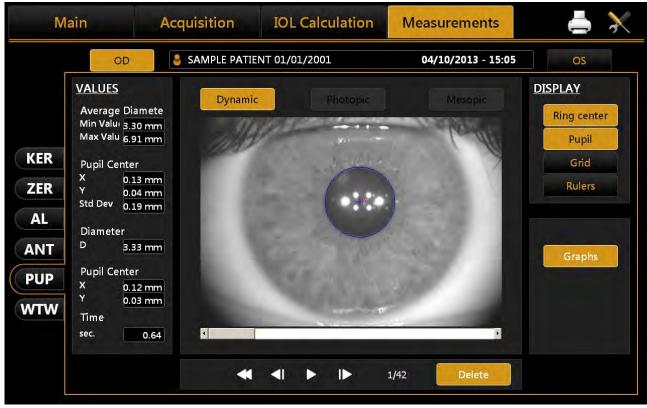

Normally, if the pupillometry is acquired, the software goes into dynamic mode (Figure 70).

Figure 70

Click on **"OD"** or **"OS"** to display the pupillometry of the right or left eye, respectively.

With the patient's eye in view, buttons are located below the home screen. These buttons are used to navigate between the acquired frames. The current frame is shown next to the buttons.

Below a summary of the functions available on this screen:

Display

- Ring Center: Shows the position of the fixation point
- Pupil: Shows the blue ring, which highlights the pupil's edges
- Grid: Shows an overlaid grid
- **Rulers:** Shows calibrated rulers

Sequences

The user can select the sequence of images to be displayed using the buttons at the top:

- **Dynamic** Clicking on the **"Dynamic"** button to display the dynamic pupillometry in the left column, the following information will also be displayed:

- Average: Value of the maximum and minimum pupil diameter measured in all the images acquired during the sequence
- Pupil Center: Cartesian coordinates of the average pupil center and its standard deviation
- Diameter: Pupil diameter for the frame selected
- Pupil Center (frame): Cartesian coordinates of the center of the pupil for the frame selected
- **Photopic, Mesopic** By clicking on the *"Photopic", "Mesopic"* buttons static pupillometry acquisitions will be displayed, with the following information:
 - Value of the average pupil diameter measured in all the images acquired during the sequence.

The other information is the same as that already described for the dynamic pupillometry. The active buttons are those for which at least one acquisition is present.

<u>Delete</u>

Pressing the *"Delete"* button, the system cancels the current pupillometry frame and the data it contains.

<u>Graphs</u>

Pressing the **"Graphs"** button displays the graphs relating to the pupil. This function is explained in the next paragraph.

12.5.1 Graphs


In this section three types of graph are displayed:

- Decentration (Figure 71)
- Latency (Figure 72)
- *Statistics* (Figure 73)

In all these graphs you can select which eye you want to analyze by clicking on "OD" or "OS".

The *"Close"* button closes the graphs.

Decentration

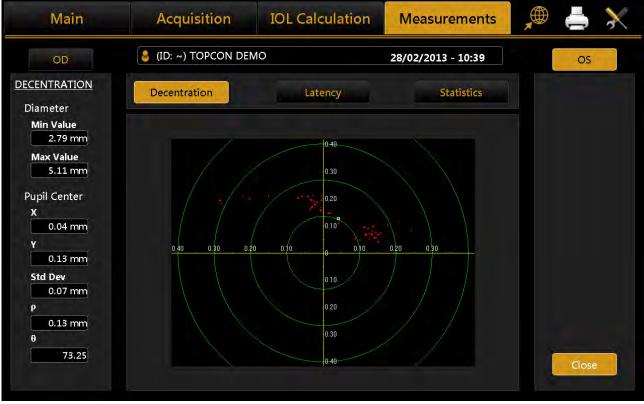



The green concentric circles identify the decentration of the pupil center with respect to the fixation point. The red dots, on the other hand, represent the coordinate variations during acquisition of the dynamic pupillometry.



Figure 72

The graph shows the time in seconds on the abscissa and the pupil diameter in mm on the ordinate, in a scale standardized on the maximum and minimum value recorded. Next the progression of the pupil's diameter over time is represented.

Taking into account that dynamic pupillometry consists of acquiring various images in variable light conditions, from mesopic to photopic and back to mesopic, on the **"Settings"** screen you can set the acquisition times for each mode (explained later). The left column shows the key to the graph.

Red: for acquisition in mesopic light conditions.

Green: to indicate the pupil contraction phase following a change in brightness brought about by the LEDs coming on.

Blue: for the pupil dilation phase following the change from LEDs on to LEDs off.

Remember that these graphs are only available if the acquisition of the dynamic pupillometry has been performed.

Statistics

Main	Acquisition	IOL Calculation	Measurements	\Rightarrow
QD	SAMPLE PATIENT 0	1/01/2001	04/10/2013 - 15:05	OS
	Decentraliz.	Latency	Statistics	
Mean 5 39 mm P10% 3.73 mm P25% 3.87 mm P75% 4.12 mm P90% 4.24 mm Mean 6.24 mm P10% 5.12 mm P25% 5.45 mm P75% 7.57 mm P90% 7.80 mm	900 900 700 800 800 400 400			d# 2575 d# 2575 d # 1030 d # mer
		Рнохоро Меноро		Cloxe

Figure 73

The graph represents the static value of the percentile of the sample for each acquisition in controlled light conditions.

As indicated in the key on the right-hand side and by the values detailed on the left, the red line represents the average value of the sample, the blue frame the value interval between the 25% and 75% percentiles, the green line the value interval between the 10% and 90% percentiles, and the red circle the values outside this interval.

The graph is displayed only if images of the pupil have been acquired in photopic or mesopic conditions.

12.6 White To White (WTW)

The White to White section allows you to view the value of the corneal diameter calculated from limbus.



Figure 74

Clicking on the **Open** button in the **Edit** menu, the user can manually reposition positional indicators in order to refine the diameter measurements.

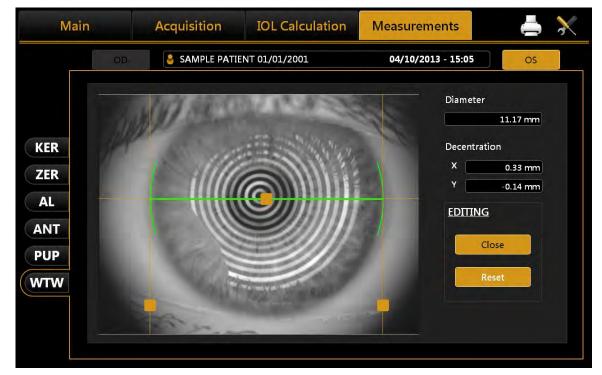


Figure 75

114

Next to the image, obtained by automatic white to white calculation, you can see:

- Corneal diameter;
- **Decentration**: deviation from the center of the iris with respect to the fixation point.

By changing the indicators position also values of corneal diameter and offset of the visual axis x and y are updated.

The Reset button restores all values to the ones obtained by the automatic calculations of system.

13 SETTINGS

cess the "Settings" section, press the ី but	ton.
General Measurements Surgeons	IOL Connectivity Admin
LANGUAGE	KEYBOARD LAYOUT
English ▼ Set Current language English	QWERTY Set Current layout QWERTY
DATE Date Format dd/mm/yyyyy ▼ Set Current Datetime 14/06/2016 20:52 Edit	POINTING DEVICES Mouse Cursor Visibility © Enabled O Disabled
PATIENT REQUIRED FIELDS	OD/OS NOTATION
O ID only	O OD/OS Native Language
	Close

Figure 76

The settings screen is divided into the following categories.

- General
- Measurements
- Surgeons
- IOL
- Connectivity
- Admin

From each settings environment you can close and return to the previous activity by selecting the **"Close"** button.

13.1 General

Refer to Figure 76:

Language: The first time the program is started, the default language set is English and the keyboard layout is "QWERTY".

To change the language settings, select the desired language from those that appear by clicking on the $\underline{\psi}$ button, press **"Set"** to set automatic start with the chosen language. It is suggested to reboot the device to apply all the settings.

Keyboard Layout: To change the keyboard layout, select the desired layout and press "*Set*". You can display the update of the layout in the personal details window ("Main").

Date: Choose the desired date format and press on the **"Set"** button. You can also set the current system date and time by clicking on the "Edit" button.

Pointing Devices: Toggles the mouse cursor (on or off).

OD/OS Notation: Toggles between two different notations, OD/OS will show the Latin notation to indicate which eye is being acquired. The native language will depend on the words for left and right.

Patient Required Fields: Toggles between two different options of required fields for creation of new patient's details. With ID only the ID is the only required field to insert when creating a patient. With this option the patient list is by default shown by ID (can be changed to Surname and Name in the patient list view).

13.2 Measurements

The acquisition settings panel allows you to set parameters for display of the corneal map, the printout and acquisition and display of the pupillometry.

MAP			PUPILLOMETRY	
Type	O Tangential O Millimeters		Display □ Grid ☑ Ruler	 Ring center Pupil
O Absolute	O Normalize	d	NUP OPTION	
O Classic	O ISO	O ISO 2005	MAP OPTION Keratometry	Map Draw
Cylinder Notation			O Sim-K	Meridians
O Positive	O Negative		O Meridians	🗹 3 Zones
Refractive Index			Emimeridians	🗖 Ruler
1.3375 Asphericity ⊚ e O S	F Op	OQ	⊚ 3-5-7 mm) 2-4-6 mm	
			Pupil Decentration Forma	
ACQUISITION Default Sequence	© K-AL-ANT	O K-AL	O Cartesian 🛛 🛛 P	'olar

Figure 77

13.2.1 Map

Туре

Select a map type:

- Axial
- Tangential

Scales

Select a scale measure:

- Diopters
- Millimeters

Select a scale type:

- Absolute
- Normalized

Select a scale color map (see paragraph 13.2.5 for details):

- Classic
- *ISO*¹
- ISO 2005²

Cylinder Notation

Select the type of cylinder notation:

- Positive
- Negative

Refractive Index

Select the refractive index to work with. You can choose from 5 indices:

- 1.3315
- 1.3320
- 1.3360
- 1.3375
- 1.3380

Be careful because changing the index will result in a differently calculated Keratometry and Topographic Map. Change this value at your discretion.

Asphericity

Select an asphericity unit of measure:

- e
- SF
- p
- Q

13.2.2 Acquisition

Select the default Biometry sequence acquisition:

- **K-AL-ANT:** the default acquisition sequence is Keratometry->Axial Length->Anterior Segment (CCT-ACD-LT)
- **K-AL:** the default acquisition sequence is Keratometry->Axial Length. Anterior Segment (CCT-ACD-LT) can be performed by pressing on the relative button

13.2.3 Pupillometry

Select one or more items with which to personalise the display of the pupil images:

- Grid
- Ruler
- *Ring center*: the center of the pupil (in blue) and the fixation point (in red) will be displayed
- *Pupil*: the outline of the pupil will be displayed in blue

¹ ISO 19980:2012(en) Ophthalmic instruments — Corneal topographers

² ISO 19980:2005(en) Ophthalmic instruments — Corneal topographers

13.2.4 Map Option

Map Design

Select one or more items with which to personalise the map display:

- Meridians
- 3 Zones
- Ruler
- Grid

Keratometry

Select one of the keratometric indices:

- Sim-K
- Meridians
- Emimeridians
- **3-5-7 mm or 2-4-6 mm**, selects the 3 Zones diameters to which the Meridians or Emimeridians values are displayed

Pupil decentration format


Display the pupil center decentration in Cartesian or polar coordinates in the application and Aladdin report

- Cartesian
- Polar

13.2.5 Topography Map Color scale description

Previous software versions have always used the same color scale for Topography maps (absolute and normalized), this is called "Classic" color scale. This software version introduces the possibility to switch the scale to the ISO³ and ISO 2005⁴ color scale, by going to Measures section of the settings. The color scale option selected affects any topographic map drawing in the ALADDIN HW3.0 application and in the printed reports (also the custom ones).

The Classic and ISO absolute color scale are shown in the following figure for the same topography map.

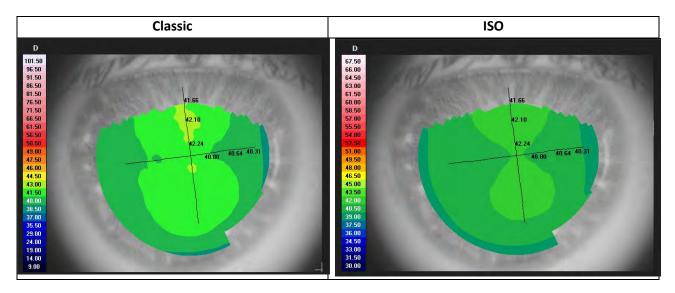

³ ISO 19980:2012(en) Ophthalmic instruments — Corneal topographers

⁴ ISO 19980:2005(en) Ophthalmic instruments — Corneal topographers

	Classic	D	101.50	96.50	91.50	86.50		76.50	71.50	66.50	61.50	56.50	50.50	49.00	47.50	46.00	44.50	43.00	41.50	40.00	38.50	37.00	35.50	29.00	24.00	19.00	14.00	9.00
Absolute scale	ISO	D	67.50	66.00	64.50	63.00	61.50	60.00	58.50	57.00	55.50	54.00	52.50	51.00	49.50	48.00	46.50	45.00	43.50	42.00	40.50	39.00	37.50	36.00	34.50	33.00	31.50	30.00
	ISO 2005	D	65.50	64.00	62.50	61.00	59.50	58.00	56.50	55.00	53.50	52.00	50.50	49.00	47.50	46.00	44.50	43.00	41.50		38.50	37.00	35.50				29.50	28.00
	Classic	D	44.75	44.50	44.25		43.75	43.50	43.25	43.00	42.75	42.50	42.25	42.00	41.75	41.50	41.25	41.00	40.75	40.50	40.25	40.00	39.75	39.50	39.25		38.75	38.50
Normalized (adjustable) scale	ISO	D	44.75	44.50	44.25	44.00	43.75	43.50	43.25	43.00	42.75	42.50	42.25	42.00	41.75	41.50	41.25	41.00	40.75	40.50	40.25		39.75	39.50	39.25	39.00	38.75	38.50
	ISO 2005	D	44.75		44.25		43.75	43.50	43.25	43.00	42.75	42.50	42.25	42.00	41.75	41.50	41.25	41.00	40.75	40.50	40.25	40.00	39.75	39.50	39.25	39.00		38.50

13.3 Surgeons

The **"Surgeons"** panel allows you to create different user profiles.

General	Measurements	Surgeons	IOL	Connectivity	Admin
SURGEON DATA					
New Surgeon Surgeon Demo Surgeon Gene Surgeon Samp	eric		Last Name Name Address		
			Date Of Birth	Sex O Male	O Female
Add	Edit	Delete			
					Close
		Fi	aure 78		

Figure 78

Select the surgeon in the left column to display the data. Press the **"Edit"** button to modify the data entered.

Press the "Add" button to add a new surgeon.

Once you have entered/modified the data, press the **"Back"** or **"Save"** buttons, respectively, to cancel or save the data.

13.4 IOL

General	easurements	Surgeons	IOL	Connectivity	Admin
General		Preset	IOL List		
Surgeon Surgeon Generic		•	Disclaimer This IOL calculator is inte		
Calculation Formulas Haigis Hoffer Q Holladay I SRK II SRK/T Enable Abulafia-Kod 	Barrett Univ Correction	versal II	physicians in selecting the The tool is intended to be comprehensive ophthalm and measurements for pa- surgery. The calculator's or medical instruction fro can guarantee that the ca Medic or Physician who decisions regarding the pa- refractive outcome. By us hold TOPCON from any ca this calculator.	e used in concurrency wi nic examination, specific atient who need to perfor results are not intended om TOPCON or be definit alculator will be accurate uses this tool must arrive patients and they are res sing this tool, you agree	ith a diagnostic tests rm cataract to serve as surgical tive, nor TOPCON e in every moment. e at their own ponsible for the the conditions and
Post Op Calculation Fo © Camellin Calos:					
 Currentin Current Shammas 			Visualize the disclaimer	at each IOL usage	Save
					Close

Figure 79

Here you can configure various options and presets for the IOL environment (Figure 79) associated with the surgeon selected, divided into three different environments described in detail below:

- General
- Preset
- IOL list

13.4.1 General

General environment (Figure 79) displays the terms and conditions of use of the IOL calculation section and lets you choose which formulas will be activated. You can reactivate the appearance of the Disclaimer for the selected surgeon at every IOL usage by checking the box below it and then saving.

13.4.2 Preset

General	Preset		OL List	
Surgeon		IOL	Toric IOL	Post Refractive IOL
Surgeon Generic	•			
		Olsen Sph IOL	Olsen Toric IOL	
Target Refraction				
0				
Manufacturer	Manufacturer	Manufacturer	Manufacturer	Manufacturer
Oculentis 🔹	Oculentis 🗸	Oculentis 🔹	Oculentis 🔹	Oculentis 🔹
Model	Model	Model	Model	Model
L-303 🔻	L-313 🔻	LS-313 MF30 🗸	LS-313 MF30 🛛 🔻	LS-313 MF30 🗸
Formula	Formula	Formula	Formula	Formula
Barrett Universal II 🔻	Barrett Universal II 🔻	Haigis 🔻	Barrett Universal II 🔻	SRK/T 🔻
			Res	et Save

Figure 80

Preset environment (Figure 80) is composed of 5 main sections, where each surgeon can set:

- Under **"IOL"** the preset for the **IOL Calculation** (included the Barrett Universall II spherical section of Barrett module)
- Under **"Toric"** the preset for the **Toric IOL Calculation** (included the Barrett Universall Toric and True K toric section of Barrett module)
- Under **"Post Refractive IOL"** the preset for the **Post Refractive IOL** (included the Barrett TrueK section of Barrett module)
- Under "Olsen Sph. IOL" the preset for Olsen Spherical IOL calculation section of Olsen module
- Under "Olsen Toric IOL" the preset for Olsen Toric IOL calculation section of Olsen module

Clicking on **"Save"** button the selected settings are saved. This function will be particularly useful during the IOL calculation, when each setting is reloaded each time as preset.

Pressing the "Reset" button, previous selections are reset, deleting all presets associated with the surgeon.

"IOL" and "Post Refractive IOL" sections have the same screen, the "Toric" one is presented in Figure 81.

rgeon				
		IOL	Toric IOL	Post Refractive IOL
urgeon Generic	•			
		Olsen Sph IOL	Olsen Toric IOL	
Target Refraction:	SIA	Incision Location		
0	0	Steep Axis	• Fixed Value	
Manufacturer	Manufacturer	Manufacturer	Manufacturer	Manufacturer
Oculentis	Oculentis	Oculentis	Oculentis	Oculentis
Model	Model	Model	Model	Model
LS-313 Tx 🔻		▼	▼	▼
Formula	Formula	Formula	Formula	Formula
Haigis 🔻	▼	▼	▼	
oric Calculator Selectior	1			
	ntor 💿 Ocule	entis Calculator	Res	set Save

Figure 81

Toric preset uses a different set of lenses (with only toric ones) and requires some additional settings with respect to "IOL" and "Post Refractive IOL" sections.

In particular, you can specify the "Surgical Induced Astigmatism (SIA)" induced by the surgeon and the "Incision Location" used during the surgery. You can decide to set the "Incision Location" at the Steep Axis of the Keratometry or at a Fixed Value to be specified.

The "Toric Calculator Selection" is to choosen between:

- "Generic Toric Calculator", that allows to define the lenses preset from a fully customizable collection of models and manufacturers.
- "Oculentis Calculator", that allows to limit the collection of available lenses to Oculentis manufacturer only. In this case the manufacturer selection is locked both in the preset set-up and in the Toric IOL calculation environments.

With the Oculentis Toric Calculator option selected you have access to the additional following functionalities:

- Review the Toric IOL alignment image with the specific Oculentis aspect also in the Toric calculation interface
- Print the Oculentis specific order forms also using the Barrett Universal II Toric Calculator (if activated)
- Get the Lentis Comfort Toric (LS-313MF15Tx) Plausibility Checkup Report printed out together with the Order Form

This setting is per Surgeon.

In the "Olsen Sph. IOL" and "Olsen Toric IOL" preset the formula is locked to Olsen formula.

General		innit		IOL List								
Surgeon Generic			205	- 1	Toric ICt.	i in	ia liiveaved	au i				
Target Refraction	5A 0		Main		Acquisition	,	IOL Calco	ulation	Measure	ments		X
Manufacturer	Manifecturer			0D	S TOPCON D. O	1/01/1	950		10/02/2015	-17:55	05	
Oculentia	Oculentia		Da	64 ()	IOL	Calcula	dian	Table 3	Of Orlightness		Post Refca	tive 101
Model (U-313 T	Model		Surgeon Gen	eric		2	Measures AL (mm)	21.71	N3 (0)	40.74	CR. (0) 1	45ax 98
Pormula Halgis +	Pormula		Terget (D)	SLA (D) 0	11 (T) 22		ACD (mm) LT (mm)	3.14 4.04	K2 (D) CCT (mm)	42.19 0.544	WTW press	11.69
			Oculentis		Oculentis		Oculentis		Oculentis		Oculentis	
			10-313 7					•				
		-	Haigis			•						
			22.19	40-1405 41-3408 42-3108	Kik 9 Tarpet		SPA IN Terget		All of Second		it) in farger	
		i i	XOL (D) 21522	REF (D) 0415	IOL (D) HE	F (D)	IOL (D)	HEF (D)	30L (D)	413 (D)	10L (0)	REF (D)
			22.05 12.05	0.07	-	-	-		-		-	
			22.29	-0.54				_	-	-		

Figure 82

13.4.3 IOL list

In this section (Figure 83) you can manage IOL spherical and torical lenses list. You can change the available manufacturers and models by adding, deleting or editing them.

For each lens you can display and edit the constants used in each formula.

The two main environments, the "Spherical" and the "Torical" ones, can be selected clicking on the corresponding button and have a similar layout.

13.4.3.1 IOL Spherical List

The **"Spherical"** layout is shown in Figure 83.

On the left side you can find the list of manufacturers, in the center their related spherical models while in the right side the calculation constants of the selected lens.

ALADDIN - User manual

General	Preset IOL List	
Suraeon		
Surgeon Generic	Spherical Torical	Export Post Op
MANUFACTURER AND MODEL	Acri.Lens 11C	. 118.000
.ZEISS	Acri.Lens 11C	nu A
1stQ	Acri.Lens 12C SRI	KII A 118.900
Aaren	Acri.Lyc 45LC SR	K/T A 118.900
Alcon	Acri.Lyc 45S	
AMO		fferQ pACD 5.460
ARGONOPTICS AUROLAB	AT LISA 809M(AT LISA 366D) AT LISA tri839MP	lladay SF 1.720
Bausch&Lomb		igis a0 1.210
Biotech	CT 47S(Acri.Lyc 47S)	igis a0 1.210
Ciba		igis a1 0.400
Corneal	OT ACOUNTA 400M (Apri Ornert 46LO)	igis a2 0.100
Croma	CT ASPHINA 509M (Acri.Smart 36A)	
Curamed	CT ASPHINA 509M (India) Car	mellin Calossi A 118.000
DISTRA	CT ASPHINA 603P (XI_Stabi ZO)	ammas A 118.000
		rrett LF 1.850



List of functions for IOL manufacturer column:

- Add: add a new manufacturer not present in the current list
 - o Insert the manufacturer name
 - o Insert the model name
 - o Insert the kind of formula and constant
 - Insert the value of the constant, "Hoffer Q pACD" in the case below (other constants will automatically be converted).

ADD				
Manufacturer				
TEST				
Model				_
TEST				
Constant				
Hoffer Q		pACD	4	
		Betk		
	E	igure 84		

- Edit: edit the name of the current manufacturer in the list
 - **Delete:** delete a manufacturer. Please note that this function will also delete every IOL associated to the current manufacturer.

List of functions for IOL model column:

- **Add:** add a new IOL model to the current manufacturer:
 - o Insert the model name
 - o Insert the kind of formula and constant
 - Insert the value of the constant, "Haigis AO" in the case below (other constants will automatically be converted).

Manufacture	5 .			
ZEISS				
Model				
TEST				
Constant				
Halgis		Au	1.3	



- **Edit:** edit the name of the current IOL
- Delete: delete the selected IOL
- **Calculation Constants History:** gives information on any change of calculation constants or optical advanced parameters values as shown in Figure 86.

The possible sources of edit are "Manual" (constants values changed manually by the user), "ULIB" (constants values changed after an ULIB update) and "Restore" (constants values restored by the user to a previous version). To restore a previous version you need to select the version to be restored and then click on the yellow right-arrow and then click on save button.

culentis"L 312		Selecte	a l		Curr	ent
93/10/2017 10:54:07 93/10/2017 10:53:51 93/10/2017 10:53:41 18/07/2017 11:13:41	Manual Restore Manual ULIB	SRKII A SRK/T A HofferQ pACD Holladay SF Haigis a0 Haigis a1 Haigis a2 Camellin Calossi A Shammas A Barrett LF	118.700 118.500 5.261 1.500 -2.476 0.046 0.300 118.000 118.000 1.622	•	SRKII A SRK/T A HofferQ pACD Holladay SF Haigis a0 Haigis a1 Haigis a2 Camellin Calossi A Shammas A Barrett LF	118.700 118.500 5.261 1.500 -2.476 0.046 0.300

Figure 86

With the "**Ulib**" button you can import .zip files of ULIB (User Group for Laser Interference Biometry) format. Download the file and copy it to the root (main card) of an empty <u>FAT32-formatted USB pen</u>

- 1. Insert the USB pen in the ALADDIN HW3.0 device.
- 2. Click on the *"Ulib"* button on the "IOL list" panel.
- 3. Select which data to import among the list of source tables available in the downloaded package.



13.4.3.2 IOL Spherical model Advanced paremeters (Olsen)

To use an IOL model with the Olsen formula, advanced optical parameters for the IOL model must be defined.

Advanced IOL model parameters are not contained in the Ulib database and they will not be changed or added to any IOL model when performing the Ulib update.

Aladdin by default provides these parameters for a subset of the IOL models. The user, accordingly to his knowledge and needs, can edit these parameter or add them for a model that didn't have them by default.

To edit or insert the advanced IOL parameters for an existing (or newly created) IOL model select the desired

IOL model and press the edit button

	EDIT	Constants	Olsen
ure	Manufacturer	Optic Material: Other	
AN	Oculentis	Refractive Index*	1.461
AIN EXe		Avg. Thickness'	1.04
ya Imi		Sperical Aberration	0.00
/lee	Model	A Const	118.5
ns.	and a second s	JOL ACD Const (Olsen)	4.20
BI	1-312	Optic Configuration? Use	r Defined
)-" edi		Front Radius (mm) Back	Radius (mm)
edi R∕		11.32	1.32
orc		* Required Field	
de			
i l			
Č.			
-		Bark	Save

Optic Material and Refractive index: the user can select the optic material that correspond to the selected IOL model to get the Refractive index field value filled automatically. Otherwise, the user can select "Other" from the list to specify a different refractive index. The refractive index value is used in the Olsen formula to perform the IOL calculations. It doesn't affects other formulae.

Aladdin provides standard refractive index values for the following materials:

- Acrylic
- Hydrophobic Acrylic
- PMMA
- Hydrophilic Acrylic
- Silicone
- Collamer
- HEMA

It is recommended to obtain the exact value from the manufacturer; this information is usually provided in the IOL model datasheet.

Avg. Thickness, Front Radius and Back Radius: these values must be specified with respect to the same IOL power. It is recommended to specify these values with respect to an average IOL power value of 22D.

Spherical Aberration: [optional] the spherical aberration value must be provided by the manufacturer according to the spherical aberration correction applied to the IOL model. Typical values are -0.27 μ m for a full correction. Zero value is for a model that does not add or subtract the spherical aberration to the eye. If the value is left empty, a 0.1 μ m value of spherical aberration add to the eye is assumed.

A Const: optimized A Const value used for SRK/T formula on the selected IOL model. Can be inserted to get the IOL ACD Const (Olsen) automatically calculated from this value if a directly optimized value for IOL ACD Const (Olsen) is not available.

IOL ACD Const (Olsen constant): if a directly optimized value for IOL ACD Const (Olsen) is available can be inserted and used without the need of inserting a corresponding A Const value.

Optic Configuration: this value is to indicate the relationship between the front IOL curvature and the back IOL curvature. Available values in the list are:

- Biconvex 1:1, front radius is positive back radius is negative
- Biconvex 2:1, front radius is positive back radius is negative
- Biconvex 3:1, front radius is positive back radius is negative
- Biconvex 1:2, front radius is positive back radius is negative
- Biconvex 1:3, front radius is positive back radius is negative
- PlanoConvex, front radius is 0 back radius is negative
- ConvexoPlano, front radius is positive back radius is 0
- ConvexoConcavo, front and back radius are positive values

13.4.3.3 IOL Torical List

The "Torical" layout is shown in Figure 87.



Figure 87

List of functions for IOL manufacturer column:

Add: add a new manufacturer not present in the current list

- Insert the manufacturer name
- o Insert the model name
- o Insert the kind of formula and constant
- Insert the value of the constant, "Hoffer Q pACD" in the case below (other constants will automatically be converted)
- Define the "Sphere Power Range", inserting the minimum, the maximum and the step of the spherical power of the lens
- Choose the "Cylinder Definition" of the lens, with "Sub models" or "Cylinder Range Based".
 - If you select "Sub models" (Figure 88) you can add a list of sub models each of them with a different cylinder value of Toricity, using the IOL button under the table (add, edit and delete).
 - If you select "Cylinder Range Based" (Figure 89), you have to insert the minimum, the maximum and the step of the cylinder value in order to define the toricity range of the lens.

anufacturer NEW MANUFACTURER		Calculation Constan	ts
lodel/Series		SRK/T	*
NEW MODEL		Value A	118.2
ylinder Definition			- and the second
Submodels		Constants	Ölsen
phere Power Range		SRK/T A	118.2
MIN MAX 10 40	0.5	HofferQ pACD	5.01
ric Lens Submodel List		Holladay SF	1.23
Lens SUBMODEL TO	Toricity 0.75	Haigis a0	0.78
SUBMODEL TI	1,5	Haigis al	0.4
		Haigis a2	0.1
		Barrett LF	1,465

Figure 88

and a statements		
NEW MANUFACTURER	Calculation Constant	5
Model/Series	SRK/T	
NEW MODEL	Value A	118.2
Cylinder Definition		LIGIT
Cylinder Range Based 🔹	Gonstants	Olsen
Sphere Power Range		118.2
MIN MAX STEP 10 40 0.5	SRK/T A HofferQ pACD	5.01
Cylinder Power Range MIN MAX STEP	Holladay SF	1.23
0.75 7 0.5	Haigīs a0	0.78
	Haigis al	0.4
	Haigis a2	0.1
	Barrett LF	1.465

Figure 89

- Edit: edit the name of the current manufacturer in the list
- **Delete:** delete a manufacturer. Please note that this function will also delete every IOL associated to the current manufacturer.

List of functions for IOL model column:

- **Add:** add a new IOL model to the current manufacturer, with the same procedure described in the "Add Toric Manufacturer" section above.
- Edit: edit the name and the properties of the current IOL
- Delete: delete the selected IOL
- **View Properties:** visualize properties, calculation constants and list of sub models of the lens (or cylinder range)
- **Calculation Constants History:** gives information on every change of calculation constants values. The layout is similar to the one in Figure 86.

The possible sources of edit are "Manual" (constants values changed manually by the user) and "Restore" (constants values restored by the user to a previous version). To restore a previous version you need to select the version to be restored and then click on the yellow right-arrow and then click on save button.

13.4.3.4 IOL Toric model Advanced paremeters (Olsen)

Refer to section 13.4.3.2 for details.

13.5 Connectivity

This panel (Figure 90) allows you to configure all the settings relative to network connectivity with external softwares or storage destinations.

General Measurements Surgeons	IOL Connectivity Admin
NETWORK FOLDER Current Network Folder	EXPORT TO EXTERNAL SW Connection to external software settings
Username	 Perform Exportation when Saving Configure Don't Ask for Confirmation when Exporting
Password Verify	IMAGEnet 6 © Enabled O Disabled Current Web Server Address
XML EXPORT Export to shared folder	http://10.1.1.66/IMAGEnet/Exam
Enable O Disable	DICOM Config. DICOM
i-base INTEGRATION © Enabled O Disabled	Enable DICOM Pending request
	Exportation Mode
Current i-base IP Address	Encapsulated Pdf Document
	Close

Figure 90

13.5.1 Network folder configuration

The "NETWORK FOLDER" panel allows the user to configure and use a remote network folder to store ALADDIN HW3.0 reports.

That resource will then become selectable as a destination in the report's print form.

In order for ALADDIN HW3.0 to be able to connect to the remote network folder, you must configure ALADDIN HW3.0 setting the correct access credentials for the remote resource.

Configuration parameters:

• **Network folder path:** the path to access the network folder location (without trailing backslashes) eg.

\\10.0.0.81\path_to\AladdinSandbox

\\TopconNetwork\path_to\AladdinSandbox

• Username: specify the domain name if needed

eg.

TopconDomain\username

• **Password:** for the specified username

When you click on the "Configure" button the system starts searching for the network resource. This procedure may take some time depending on the network. Failure or success to connect to the network resource is reported as shown in Figure 91. Connection failure may be due to unreachable resource path or to wrong credentials.

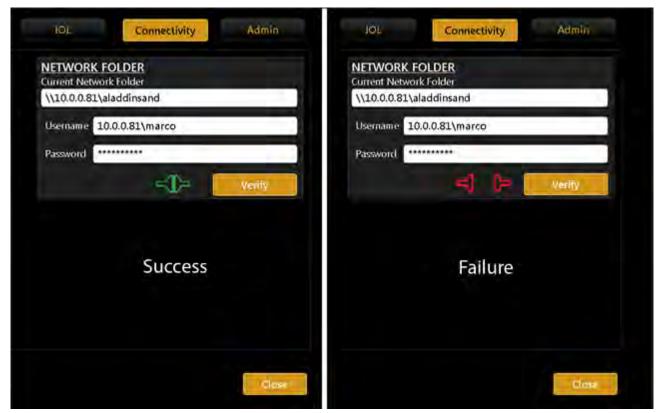


Figure 91

13.5.2 XML Export

Enables/Disables XML option for exporting XML data of the exam to the network folder by the export window.

13.5.3 IMAGEnet i-base software

ALADDIN HW3.0 can receive and transfer data to Topcon IMAGEnet i-base through a wireless or LAN network. IMAGEnet i-base is activated by clicking on the Enabled Option and by providing the IP address of the machine we want to connect to. By clicking on the Configure button (Figure 94) the user is presented with a list of IP Addresses of the machines running IMAGEnet i-base software that are reachable from ALADDIN HW3.0. Once selected the proper IP(Figure 92), ALADDIN HW3.0 is ready to exchange data with IMAGEnet i-base machine.

NETWORK FOLDER	SERVER LIST	-		
Current Network Folder	192.168.56.1			
\\server\path_to_serve			Search	1 million (1)
			Search	Configure
Username domain\user				
Password *******				
XML EXPORT				m
Export to shared folder	4			
Enable				Configure DICOM
BASE INTEGRATION				
Enabled				Pending request
© Enableu (Patient List Max Size	
Current IBase IP Addres:			100	/ Capture Image ▼
192.16				
			Close	
				Close
				and a start

13.5.4 IMAGEnet 6 Server software

ALADDIN HW3.0 can receive and transfer data to Topcon IMAGEnet 6 Server through a wireless or LAN network. IMAGEnet 6 Server is activated by clicking on the Enabled Option and by providing the IP address of the external server we want to connect to. Once we have selected the proper IP, ALADDIN HW3.0 is ready to exchange data with the IMAGEnet 6 Server machine (Figure 93).

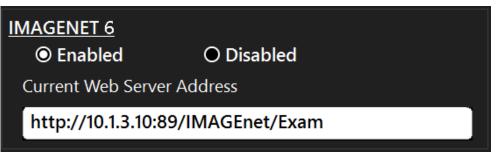


Figure 93

13.5.5 Export to External Software settings

•	
EXPORT TO EXTERNAL SW	
Connection to external software settings	
Perform Exportation when Saving	Configure
	coningate
Don't Ask for Confirmation when Exporting	

In this panel it is possible to control two main aspects of the device behaviour in relation to the exportation function.

- *"Perform Exportation when Saving"*, allows to be prompted automatically the exportation form once saving a new or existing exam.
- *"Don't ask for Confirmation when Exporting"*, allows to skip the selection of targets when performing the exportation while saving or when pressing on the export button

It is also possible to access configuration ("Configure" button) of the other external software destinations which are not present in the other panels of this section.

By activating the checkbox of an export destination, this is included in the targets of the exportation function.

General	Measurements	Surgeons	IOL	Connectivity	Admin
NETWORK FOLDE	older		EXPORT TO EXTERN		
Username visia\a	tenti SEND TO arrigu 🛛 DICOM St	orage Server		ع ۱ Exporting	Configure
Password ******	** 🛛 i-base	ics®	Configure	isabled	
XML EXPORT Export to shared for © Enable					
I-BASE INTEGRAT	ION		Clo		Config. DICOM Pending request
Current i-base If		nfigure	Encapsulated Pd	If Document	
					Close

PhacoOptics®: <u>http://www.phacooptics.net/</u>

In order to get exam data to exported to PhacoOptics[®] software, running on an external PC, it is necessary to configure a network path and access credentials for the target of exportation. Refer to PhacoOptics[®] manuals on how to get the data exported from ALADDIN HW3.0 inside the application.

NETWORK FOLDER Current Network Fold VISIAOPTHLAP, Username visia \arr Password ******* Username domain\username Password ******* VML EXPORT Export to shared folde O Enabled IBASE INTEGRATION © Enabled Current IBase IP Aq	General	Measurements	Surgeons	IOL	Connectivity	Admin
XML EXPORT Export to shared folde Configure O Enable Infigure DICOM IBASE INTEGRATION Inding request O Labled Inding request	Current Network Fol \\VISIAOPTHLAP Username visia\a	Current Network	Folder tination			Configure
IBASE INTEGRATION Inding request Image: State of the state of th	Export to shared fol		***	-(1)	Configure	Infigure DICOM
	● Enabled				CLOSE	

13.5.6 DICOM

The DICOM panel of Connectivity section allows to set the needed parameters for the connections to the available DICOM services.

The available services are:

- **Modality Worklist**, The DICOM Modality Worklist service provides a list of imaging procedures that have been scheduled for performance by the acquisition device.
- **Patient Root Query,** This enables the device to find patient's details from a DICOM server.
- **Storage,** The DICOM Store service is used to send images or other persistent objects (structured reports, etc.) to picture archiving and communication system (PACS) or workstation.
- **Storage Commitment** The DICOM Storage Commitment service is used to confirm that an image has been permanently stored by a device.

For each services the needed parameter are:

- Remote Application Entity (AE) title
- Remote IP address
- Remote connection port

The "Local Application Entity title" is the identifier name through which the device presents itself to the servers.

The "*N*-EVENT Report node port" is the port at which the device is able to receive N-EVENT REPORTS for storage commitments (default is 115).

Local IP Address
10.1.1.12
Local N-EVENT Node Port
115
Query/Retrieve SCP
Application Entity Title
DVTK_QR_SCP
IP Address
10.1.3.10
Port
106
Connection Test Storage Commitment SCP
Application Entity Title
DVTK_STRC_SCP
IP Address
localhost
Port
105
Connection Test

The connectivity to the defined server can be tested using the "C-ECHO" function activated by the relative "*Connection Test*" button. The result of connection test is shown by the green or red icon.

Application AE	
WL_SCP_AE_TITLE	
IP Address	
10.1.1.30	1.2
Port	
107	117
- <u>×</u> -	Connection Test

In order to configure properly the full DICOM workflow it could be necessary to perform some operations or configurations on the server's side. In order to do this contact the System Administrator.

The DICOM module of ALADDIN HW3.0 is describe in detail in its **DICOM Conformance Statement**. Visit <u>http://www.topconmedical.com/conformance.cfm</u> to download it.

13.6 Admin

General Measurements Surgeons	IOL Connectivity Admin
INFORMATION S/N: ALA20-HW4VMMA S/V: 1.9.0 ALPHA2 REPORT Image: Constraint of the second sec	CALIBRATION Check
Saving Exam After Print O Enable O Disable Header	APPLICATION Privacy
Topcon Europe Medical BV Output Format	Backup_Restore Write Filter is ON
Output Filename Convention	Upgrade Disable Filter
Standard O Alternative	Close App Admin Mode
Enable RX/AL Trend Function	Remote Assistance
	Close

Figure 94

This is the instrument's administration panel (Figure 94).

It provides certain information on the system: serial number (S/N) and software version (S/V). The *"Check"* button starts the calibration check procedure.

• Check the calibration

See the paragraph describing the procedure <u>Checking the calibration</u>.

It is absolutely essential to check the calibration when the device has been transported from one place to another and when it has suffered an impact or thermal shocks.

A It is recommended to check the measurements every day when turning on the device.

The "Enable RX/AL Trend Function" gives access to the function for reviewing biometry data over time in charts from the exam list view. See section 11.16 for further details.

The "Application" frame manages the behaviour of the integrated software:

- **Upgrade** → Updates the integrated software
- **Backup** → Starts the backup procedure on a USB driver
- Close App → Closes the application
- **Privacy** → Refer to section 13.7

13.6.1 Report

In the "Report" panel of Admin section it is possible to set a Custom header for all the reports, as well as setting the output format of reports exported to network folder. The available formats are: Pdf, Jpeg, Bmp, Tiff, Png.

It is also possible to set the automatic saving of the exam after printing.



Output filename convention

Option to choose between two filename conventions for the reports exported to network folder:

• Standard:

patientID[_patientSurname][_patientName][_patientDoB]_reportName[_eyelabel]_ExamD ateTime[_progressiveNumbert].extension

• Alternative:

patientID[_patientDoB]_Aladdin_eyelabel[_progressiveNumber].extension

eyelabel is:

- R, right eye
- L, left eye
- B, both eyes
- N, no eye related

13.6.2 Remote Assistance

In case you need remote Assistance the ALADDIN HW3.0 application integrates Teamviewer QS (Quick Support) pre-installed.

- There is NO NEED to disable the Write Filter protection
- Ensure to have available internet access for your ALADDIN HW3.0 unit
- Go to Settings -> Admin and press "Remote Assistance".
- Wait for the Teamviewer Window to open
- Communicate to the Remote Operator the ID shown under "Your ID" and wait for the incoming connection
- The password is masked, the operator knows it already

INFORMATION		CALIBRATION			
S/N: 00000000 S/V:	1.8.0	ener 201			
REPORT		OPCON			
Saving Exam After Print O Enable	AL	ADDIN			
Header	Alinw Remnte Cont			Privacy	
Topcon Europe Medical BV	Hease send your 10 to the the incoming connection .		tore	Write Filter is ON	
Output Format	vera El 13	835 401 341		Disable Filter	
Output Filename Convention	Ready to connect (secur			product inter	
	trans Sedera minera con	Cancel	P	Admin Mode	
				mole Assistance	
		1		and the second second	

To turn off manually the Remote Assistance you can close the Teamviewer window or press "Remote Ass. OFF"

If you get one of the following windows please check your internet connection to ALADDIN HW3.0 or contact your IT staff.



13.6.3 Updating the integrated software

In this section is described the software upgrade procedure from one version to the following version. To update the software, perform the following operations:

- 1. Unpack the update packet in the root (main card) of an empty FAT32-formatted USB external drive.
- 2. Switch on the ALADDIN HW3.0.
- **3.** Cancel the calibration check (Figure 95).

New	Search	Search Server	Acquisition	19/06/2016 18:43 🗙
PATIENT DETAILS				
Last Name		ID		
Name	WARNIN	G: Check calibration before sta	arting the measurements	Ok Clear
Date of Birth (dd/mn	10 March 10		nce	
				0
🗲 ТОРСС	1	2 3 4 5		0 ESC
ALADD		QWER	TV	
	A	S D F	GH	L.
	z	X C V	BNM	·

Figure 95

4. Click on the settings icon (Figure 96).

New	Search	ĩ	S	earc	h Se	rver		Acqu	uisiti	on	19,	/06/2016 1	11:01
PATIENT DETAILS													
Last Name					ID		_			-			0k
Name					Gende							F	Clear
Date of Birth	(dd/mm/yyyy)					Male Female							P
										-	-		
		1	2	3		4	5	6	7	8	3	9 0	ESC
	ΡΟΠ		Q	W		R	T		Y	U	I	0	Р
- ALADDIN					-	~	54	-					
	JUIN	A	1 5	S	D	F	G	G		J	K	L	Ę.

Figure 96

5. Click on the "Admin" tab (Figure 97).

ALADDIN - User manual

INFORMATION	CALIBRATION
S/N: S/V:	Check
00000000 1.5.0 Beta	
PNLREPORT	APPLICE
Saving Exam After Print	
O Enable	INF IS ON
Header	
Topcon Europe Medical BV	
and an and a second	Upgrade Disable Filter
Output Format Pdf	
	Close App Admin Mode
	Remote Assistance

Figure 97

6. Insert the USB stick with the "Aladdin upgrade" files in one of the USB ports (Figure 98).



Figure 98

7. Click on the "Upgrade" button (Figure 99).

ALADDIN - User manual

INFORMATION		CALIBRATION	
S/N:	S/V:	Check	
00000000	1.5.0 Beta		
PNLREPORT		APPLICATION	
Saving Exam After Pi	rint		
O Enable	O Disable	Write Filter is ON	
Header		Backup	
Topcon Europe	Medical BV	Upgrade Disable Filter	
Output Format	Pdf 🔹		
		Close App Admin Mode	
		Remote Assistance	
	1000		
	1.00	Land Contraction of the second	

Figure 99

8. Click on "**Ok**" to reboot the system and start the upgrade (Figure 100).

Check
PLICATION
Privacy
ure. The system will now Write Filter is ON
Disable Fitter
Close App Admin Mode
Remote Assistance

Figure 100

9. After the upgrade you will see the message of Figure 101 on the screen.



Figure 101

- **10.** The system will reboot and starts the "Aladdin application".
- **11.** After restart, the software updates the system; this operations could take some minutes please don't restart machine during this procedure (Figure 102).

Processing]
	1




12. If you see the message of Figure 103, please switch off your ALADDIN HW3.0 and turn it on again. Your ALADDIN HW3.0 should now work fine.



Figure 103

13. Your ALADDIN HW3.0 is upgraded. Please check in the settings, "Admin" tab if the S/V is now the new one (Figure 104).

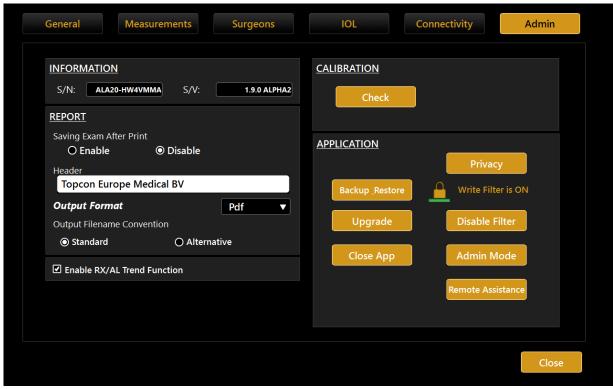



13.6.4 Backup & Restore

It's recommended to perform a backup to have a safety copy of every patient stored data. Depending on the expected size of the entire archive, we suggest to use an external USB drive or have available enough space on a designed network shared folder.

Pressing on the "Backup _Restore" button the Backup and Restore utility is opened.

It is possible to use an external USB drive or a network shared folder to backup and restore data.

Backup Contents

With this utility it is possible to backup:


- Local Exams archive: the complete list of patients and exams that are currently stored in the local database of the machine
- **Surgeon's presets and IOL collections:** Surgeons list with all the data associated with each one, such as default IOL lens presets, IOL collections (customized constants or manually added IOL models)
- Application user settings: interface settings such as visualization options, display units, scales, network settings, report header
- Machine calibration: internal calibration parameters of the machine to be stored for safety.

The machine calibration can be restored only by the technician to the same device from which the backup was made.

Backup Procedure

To perform the backup connect the desired USB external storage device or switch to Network Folder option.

- 1. Select the desired destination:
 - USB, select the desired partition from the list of the available

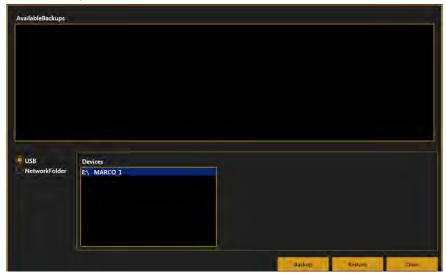


• Network folder, define (if not already defined) the desired network shared folder destination, check the connectivity.

	dalone_74150001_2016-06-20-16-56-33		
USB	NetworkFolder		
NetworkFolder	\\Server\Path_To_Backup_Folder		
	Username		
	Domain\User		
	Password		
	Password	=	Configure

2. Press "Backup" button, wait for the procedure to complete. Press "OK" to confirm the operation or "No" to do not perform the backup.

AvailableBackups AladdinBackup_Stan	dalone_0000000	00_2016-06-21-10-51-32			
AladdinBackup_Stan	dalone_0000000	00_2016-06-21-10-52-48			
				i	
		Are you sure you want to backup	2	I	
				_	
			Ok .	No	
USB	Devices		-		
NetworkFolder	E\ MARCO				
					 _

Do not turn off the device or unplug the power supply and ethernet cable while performing this operation. Do not unplug USB devices if USB is the backup destination. This may take several minutes depending on the exam archive size.

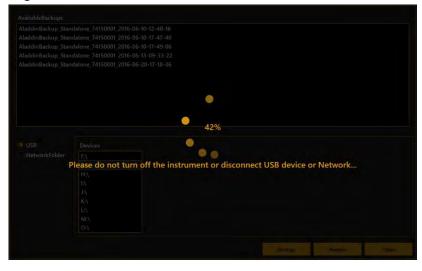

3. After the procedure is completed, a new entry appears in the list of available backup data with the following naming convention:

AladdinBackup_Standalone_<device serial number>_<backup date>

Restore Procedure

The machine calibration backup CANNOT be restored to a different machine from the original one.

If attempting to restore a different machine backup to the current machine the calibration restore is skipped. Ask to technical assistance in order to restore the calibration.

- 1. Select the source of the backup image to restore (USB or Network Folder). Navigate the list of detected backup images available at the selected source.
- 2. Select from the list the backup image you want to restore and press "Restore" button. Press "OK" to confirm the operation.

AvailableBackups							
AladdinBackup_Stan	dalone 741500	01 2016-06-10-12-44	8-16				
AladdinBackup_Stan							
AladdinBackup_Stan							
AladdinBackup_Stan							
AladdinBackup_Stan	dalone_741500	01_2016-06-20-17-18	8-36				
AladdinBackup_Stan	dalone_741500	01_2016-06-20-17-19	9-12				
● USB ● NetworkFolder	Devices F:\ G(A) H:\ I:A J:A J:A L:A M:A O:\ O:\	Are you sure yo	ou want to restore th	Dk	up?		
					Backup	Restore	Close

3. The Restore procedure starts.

Do not turn off the device or unplug the power supply and ethernet cable while performing this operation. Do not unplug USB devices if USB is the backup source.

	•	42%			
Please do not turn	off the instru	ment or disco	nnect USB device o	r Network	

- In sequence, it is prompted to confirm if restoring or not each kind of backup content.
 Press "OK" if you want to restore the content or press "No" to skip the restore of the mentioned content.
 - (1) Restore machine calibration files.

The calibration can be restored only to the device from which the backup was made originally.

This operation is allowed only to authorized technicians. Contact your distributor for assistance. The machine calibration restore is skipped if the procedure is performed by the user.

(2) Restore Application User Settings.

The current content on the machine will be overwritten.

Do you want to restore Application	user settings?	
	Ok	No


(3) Restore Surgeons settings.

The current content on the machine will be overwritten.

Do you want to restore Surgeons da	ita?	
	Ok	No

(4) Restore Local Exams archive.

The backup content will be added to the current content if not already present. This may take several minutes depending on the number of exams in the archive.

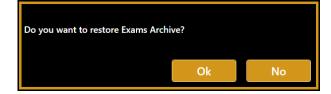

13.6.5 Shut down

Press the "*Close App*" button to close the application and return to the Windows desktop. You will be asked to confirm this operation.

Press the stand-by button to shut down the device.

13.7 Privacy & Security Settings

If the password login is kept enabled and a password has been already defined it is necessary to input the defined login password to operate the device.

PRIVACY AND SECURITY	
In this section you can configure your data protection options within this device. The main option is the definition of the user password to access the device. For further details on the data protection options please refer to sections 4.8 and 13.7 of the user manual.	
Use Password Login to operate the device Configure Password A PASSWORD NOT CONFIGURED!	
Auto LOG OFF / Screen Saver Timeout [min]:	
☐ Hide Patient Details in reports and exported exam packages	
□ Anonymize filenames of exported reports and exam packages	
Password protect exported PDF reports (use login password)	
Hide Patient Names and Disable Actions when in assistance mode	
Close	

13.7.1 Password protected Login

The options allows to protect with password the access to the Aladdin on-board application operations and data. By default the option is enabled but requires the password to be configured.

When the option is enabled and the password configured the login is prompted when (refer to section 11.1):

- the Aladdin on-board application starts
- the auto log OFF timeout occurs (if enabled)
- the lock button in the main top bar is pressed

Press on "Configure Password" button to define the preferred password.

It is recommended to use a strong password. The password is case-sentitive. No password strength checks are applied.



It is recommended to apply a password aging policy.

Type the desired password in the password field and type it again in the confirm password field.



13.7.2 Screen Saver / Auto Log OFF

The option allows to set automatic locking of the Aladdin on-board after a configurable timeout.

By default the option is enabled with timeout set to 15 minutes.

When no interaction with the device is performed during the timeout the screen saver will be displayed. When the screen saver is closed the login screen (if enabled and configured) will be displayed.



13.7.3 Hide Patient Details in reports and exported exam packages

The option, if enabled, allows to mask the patient details in the header of the printed or exported reports. Only the used patient ID will be reported normally. Date Of Birth

🚝 ΤΟΡΟΟΓ		Тор	con Europe Medical bv
Patient	*****N ***O	Surgeon	Surgeon Name
Patient ID	QREWVBQ	Exam Date	11/15/2018 - 16:57

If for assistance reasons the exam data is exported the patient details will be replaced with indications of the source device, only the patient ID will be maintained.

13.7.4 Anonymize filenames of exported reports and exam packages

If the option is selected the report output filename convention will be forced to be the alternative one, refer to section 13.6.1 for further details.

REWVBQ_20181115_165700_Aladdin_B.pdf
QREWVBQ_20181115_165700_Aladdin_B_01.pdf
QREWVBQ_20181115_165700_Aladdin_B_02.pdf
QREWVBQ_20181115_165700_Aladdin_B_03.pdf
QREWVBQ_20181115_165700_Aladdin_B_04.pdf

13.7.5 Password protect exported PDF reports


This option is available if the login password is enabled and configured. If this option is enabled the PDF reports exported to network shared folder will be accessible only by typing the user password as defined to access the Aladdin on-board application.

13.7.6 Hide Patient Names and Disable Actions when in assistance mode

The options, used in combination with password protected login, allows to make the patient details and related actions non-accessible when interacting with the Aladdin on-board application in technician/service mode. To switch the application to technician mode press on the lock button to lock the application. Let the technician staff to login to application with the technician password, the Application will switch to Technician Mode. To switch the application back to user mode press again the lock button and enter the user password.

User Mode

In normal user mode the patient details and related actions are regularly accessible.

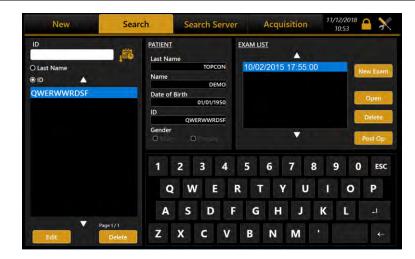

Technician Mode

In technician Mode the patients details are not accessible and delete/edit actions are not available.

When exporting exam packages in this mode the patient details are automatically anonymized.



14 CHANGING THE FUSES

Step 1 Open the fuse box cover using a screwdriver

Step 2 Take out the fuse box (use a screwdriver to release it)

Step 3

Remove the blown fuse from its seat and replace it with an identical one, as indicated in the table below and on the instrument label.

Step 4 Push the fuse box carefully back into position



Figure 105

🛆 It is mandatory to use fuses only with the indicated characteristics

The use of undersized fuses can cause the interruption of power to the device during normal working conditions. In this case there is no risk to the user, nor for the patient, but the device turns off at inopportune moments, and this can cause data loss

The use of oversized fuses can lead to damage to the internal electronics of the device due to current overload for non-interruption by fuses. In this case you do not identify risks to the patient, but on the user or damage to the device and then stop working and possible data loss.

Fuse type	Fuse value
20 x 5 mm	T 2.5 A L 250 V anti-surge

15 TECHNICAL SPECIFICATIONS



No modification of this equipment is allowed

NOTE: the manufacturer shall provide, upon request, circuit diagrams. the list of components, descriptions, calibration instructions or other information that will assist the technical assistance personnel in the repair of parts of the device specified by the manufacturer as repairable by the technical support staff.

NOTE: For the isolation of the device from the supply mains power, the device is provided with a removable power cable.

GENERAL INFORMATION					
FEATURE	MEA	ASUREMENT METHOD	LIGHT SOURCE		
	Placido disk				
	Keratometry	24 rings equally distributed in a			
	conus	43D sphere			
Corneal	Analysed points	Over 100,000 points			
topography -	Measured points	Over 6,000 points	Red LED Type1 and Type 2		
Keratometry		Up to 9.8 mm on a sphere of 8mm			
	Corneal coverage	radius (42.2 diopters with			
		n=1.3375)			
	Focusing system	Guided focus			
Pupillometry	Image analysis		Infrared and white light LED		
AL					
ACD	Low cohoronco inte	rforomotry on optical fibor	Slad @ 820mm		
ССТ	Low-conerence inte	erferometry on optical fiber	Sled @ 830nm		
LT					

OPTICAL RADIATIONS						
FEATURE	LIGHT SOURCE		WAVELENGHT	POWER ON EYE		
Central fixation LED	Yellow greer	n LED	572 nm	< 0.01 mW		
Illumination of Placido disk	Red LED Typ	e1	633 nm	<0.02 mW		
for topographic analysis	Red LED Typ	e2	615 -630 nm	<0.02 mW		
		Blue	473 nm	0.03 cd		
Pupillometric analysis	White LED	Green	532 nm	0.005 cd		
		Red	630 nm	0.008 cd		
	IR LED		780 nm	<0.4 mW		
AL measurement	Sled		830 nm	< 0.7 mW		
AL measurement	IR LED (crown)		770 nm	< 0.1 mW		
	Sled		830 nm	< 0.7 mW		
CCT-ACD-LT measurement	Red LED Typ	e1	633 nm	<0.02 mW		
	Red LED Typ	e2	615 -630 nm	<0.02 mW		

CAUTION - The light emitted from this device is potentially hazardous. The longer the duration of exposure, the higher the risk of ocular damage. Exposure to the device light when using it at maximum intensity will infringe the safety indication after a 60-minutes use. ALADDIN HW3.0 has a series of LEDs of various types and powers installed. All the characteristics are detailed in the Technical Specifications section in this manual. The LED groups comply with the emission limits for the Group 2 instruments according to the standard ISO 15004-2.

INFORMATION ON MEASUREMENTS				
	MEASURE	MEASURING RANGE	DISPLAY RESOLUTION	IN VIVO REPEATABILITY
	Curve radius	5.00 – 12.00 mm	0.01 mm	±0.02 mm
Keratometry	ry Curve Radius in Diopter (D) (n=1.3375) 28.00 - 67.		0.01 D	±0.12 D
Axial Length		15.00 – 38.00 mm	0.01 mm	±0.016 mm
Anterior Chamber Depth		1.50 – 6.50 mm 0.01 mm		±0.04 mm
Central Corneal	۲hickness	0.300 – 0.800 mm	0.001 mm	±0.02 mm
Lens Thickness	Phakic eye	1.50 – 6.50 mm	0.01 mm	±0.06 mm
	Pseudo-phakic eye	0.50 – 3.50 mm	0.01 mm	±0.06 mm
Pupil dimension		0.50 – 10.00 mm	0.01 mm	±0.05 mm
Limbus (White-T	o-White)	8.00 – 14.00 mm	0.01 mm	±0.05 mm

ENVIRONMENTAL CONDITIONS				
	IN USE	STORAGE	TRANSPORT	
Temperature	10 - 40° C	-20 - 70° C	-20 - 70° C	
Relative humidity	8-75% (non condensing)	8-75% (non condensing)	8-75% (non condensing)	
Atmospheric pressure	800-1060 hPa	700-1060 hPa	700-1060 hPa	

ELECTRICAL DATA		
Power suppl	у	AC 100-240V 50/60 Hz
Power consumption		<100 VA
Fuse Type Value	Туре	20 x 5 mm
	Value	T 2.5 A L 250 V anti-surge

	MECHANICAL SPECIFICATIONS		
	Device	Device Packaged	
Width	320mm	600mm	
Height	490mm	800mm	
Length	470mm	710mm	
Weight	18 kg	29kg	

ON-BOARD PC SPECIFICATIONS		
Operating system	WINDOWS 10 64bit	
Processor	Intel [®] Celeron [®]	
RAM	4GB	
Hard disk	At least 500GB SATA + mSATA SSD 32GB	
External connections	LAN integrated, 2x USB	

-

16 PERFORMANCE TESTS

BRIEF SUMMARY OF PERFORMANCE TESTS AND RESULTS

A prospective, single site clinical study comparing the performance of the ALADDIN HW3.0 with LENSTAR LS 900 was conducted in 66 eyes (1 eye for each enrolled subject). This study evaluated the agreement and precision in the subsequent endpoints:

Data Type	All Devices
	Axial Length (AL)
	Cylinder Axis (AX)
	Anterior Chamber Depth (ACD)
Quantitative	Lens Thickness (LT)
Measurements	Central Corneal Thickness (CCT)
	White-to-white distance (WTW)
	 Keratometry at the flattest meridian (K1_D)
	 Keratometry at the steepest meridian (K2_D)

The related unit of measurements are:

- Axial length: millimeters;
- Cylinder Axis: degrees;
- Anterior Chamber Depth: millimeters;
- Lens Thickness: millimeters;
- Central Corneal Thickness: millimeters;
- White-to-white distance: millimeters;
- Keratometry: diopters.

Analysis of Agreement

The measurements for the Agreement analysis arise from the first three acceptable measurement types acquired from each device.

For each endpoint (AL, AX, ACD, LT, CCT, WTW, K1_D, K2_D), the analysis of agreement was performed for every ALADDIN HW3.0 unit per LENSTAR LS 900 unit configuration for each operator separately.

The analysis was carried out estimating the mean difference, its standard deviation, 95% Limits of Agreement, Bland-Altman scatter and difference plot using Bland-Altman method for multiple observations on an individual (Bland and Altman, 1999).

The overall analysis of agreement is reported in the table below.

The endpoint means with their standard deviations were provided too.

Details about the Bland-Altman analysis: model terms were estimated via REML assuming exchangeable replications and a device * subjects interaction term.

The agreement between instruments is summarized in Table 1:

	ALADDIN HW3.0		LENSTAR	LENSTAR LS 900		Agreement parameters		
	ALADDIN	ALADDI	LENSTAR	LENSTAR	Difference	Difference	Lower	
Endpoint	Mean	N SD	Mean	SD	Mean	SD	LoA	Upper LoA
AL[mm]	24.04	1.43	24.05	1.43	-0.010	0.024	-0.058	0.038
AX [°]	76.73	78.26	76.56	77.56	0.174	6.176	-12.177	12.525
ACD[mm]	3.67	0.38	3.67	0.39	0.004	0.040	-0.076	0.085
LT[mm]	3.67	0.27	3.65	0.27	0.021	0.054	-0.088	0.129
CCT[mm]	0.555	0.035	0.555	0.035	0.000	0.004	-0.008	0.008
WTW[mm]	12.27	0.34	12.30	0.36	-0.032	0.068	-0.167	0.103
K1[D]	43.16	1.43	43.17	1.44	-0.010	0.105	-0.221	0.201
K2[D]	44.26	1.52	44.26	1.53	-0.006	0.178	-0.362	0.349

Table 1. Agreement between ALADDIN HW3.0 and LENSTAR LS 900.

Note to the Table 1:

- data in table concern the summarized comparison produced by 3 operators using 3 different Aladdin and 3 different Lenstar LS 900 units.

The results of the study demonstrate that the Axial Length (AL), Cylinder Axis (AX), Anterior Chamber Depth (ACD), Lens Thickness (LT), Central Corneal Thickness (CCT), White-to-white distance (WTW) and Keratometry (K1_D and K2_D) measurements of the ALADDIN HW3.0 are substantially equivalent to those of the predicate device.

Analysis of Precision

The measurements for the Precision analysis arise from the first three acceptable measurement types acquired from each device.

Analysis of precision, for each endpoint, will be performed on ALADDIN and LENSTAR separately using a REML method for repeated measures within subject including the subsequent terms: device identifier (A1, A2, A3 for ALADDIN and L1, L2, L3 for LENSTAR), operator identifier (1, 2 and 3), subject identifier (1 to 66) and every two-way interactions and the three-way interaction; a standard variance component matrix will be used as covariance structure of R-side matrix. The related results will be used to estimate: - the repeatability SD as the square root of model MSE; the repeatability limit as repeatability SD multiplied by 2.8; the repeatability coefficient of variation as (repeatability SD / mean) multiplied 100; - the reproducibility SD as the square root of (device variance component estimate + operator variance component estimate + device*subject variance component estimate + device*operator variance component estimate + device*operator variance component estimate + device*operator variance component estimate + model MSE); the reproducibility Ilmit as reproducibility limit as reproducibility coefficient of variation as (repeatability SD / mean) multiplied 100; - the reproducibility SD as the square root of (device variance component estimate + operator*subject variance component estimate + device*operator variance component estimate + device*operator variance component estimate + device*operator variance component estimate + model MSE); the reproducibility limit as reproducibility SD multiplied by 2.8; the reproducibility coefficient of variation as (reproducibility SD multiplied by 2.8; the reproducibility coefficient of variation as (reproducibility SD multiplied by 2.8; the reproducibility coefficient of variation as (reproducibility SD / mean) multiplied by 2.8; the reproducibility coefficient of variation as (reproducibility SD / mean) multiplied 100.

Variance component estimates were provided for every model term.

To assess if the variability of repeated measures within a subject is fairly constant across the range of results:

1) plot of the standard deviation of repeated results within eye (vertical axis) versus the mean of repeated results (horizontal axis) was provided for each endpoint;

2) Spearman correlation coefficients and related *p*-values for hypothesis H_0 : $\mathbb{B}_{\text{spearman}} = 0$ between mean and standard deviation were provided for each endpoint.

Analysis described in previous point 1 and 2 was provided for ALADDIN HW3.0 and LENSTAR LS 900 separately.

Details about REML analysis: model terms will be estimated by using sum of squares.

The results of the analysis of repeatability and reproducibility are summarized in Table 2:

Repeatability and reproducibility.								
DEVICE	Endnoint	Overall Mean	R	EPEATABILI	ТҮ	REF	RODUCIBI	.ITY
DEVICE	DEVICE Endpoint	Overall Mean	SD	Limit	% COV	SD	Limit	% COV
ALADDIN HW3.0	AL[mm]	24.04	0.020	0.056	0.084	0.024	0.068	0.100
ALADDIN HW3.0	AX [°]	76.73	2.859	8.004	3.726	2.918	8.170	3.803
ALADDIN HW3.0	ACD[mm]	3.67	0.026	0.073	0.708	0.026	0.074	0.721
ALADDIN HW3.0	LT[mm]	3.67	0.031	0.086	0.833	0.032	0.090	0.878
ALADDIN HW3.0	CCT[mm]	0.555	0.005	0.013	0.837	0.005	0.013	0.858
ALADDIN HW3.0	WTW[mm]	12.27	0.066	0.184	0.536	0.066	0.186	0.541
ALADDIN HW3.0	K1_D [D]	43.16	0.077	0.217	0.179	0.082	0.230	0.191
ALADDIN HW3.0	K2_D [D]	44.26	0.121	0.339	0.274	0.127	0.355	0.286
LENSTAR LS 900	AL[mm]	24.05	0.014	0.040	0.059	0.015	0.041	0.061
LENSTAR LS 900	AX [°]	76.56	5.304	14.852	6.928	5.312	14.873	6.938
LENSTAR LS 900	ACD[mm]	3.67	0.036	0.100	0.974	0.036	0.102	0.991
LENSTAR LS 900	LT[mm]	3.65	0.052	0.145	1.418	0.053	0.147	1.440
LENSTAR LS 900	CCT[mm]	0.555	0.004	0.011	0.722	0.004	0.011	0.726
LENSTAR LS 900	WTW[mm]	12.30	0.065	0.183	0.532	0.066	0.184	0.533
LENSTAR LS 900	K1_D [D]	43.17	0.114	0.320	0.265	0.118	0.329	0.272
LENSTAR LS 900	K2_D [D]	44.26	0.186	0.522	0.421	0.191	0.535	0.432

Table 2
Repeatability and reproducibility.

Note to Table 2:

- Repeatability includes variation due to measurement error;

- Reproducibility includes variations due to the device, the operator, the interaction between device and subject, the interaction between operator and subject, the interaction between device and operator, the interaction between device, operator and subject, and measurement errors;

Repeatability % COV = (Repeatability SD / abs(overall mean))*100;

- Reproducibility % COV = (Reproducibility SD / abs(overall mean))*100.

The result of the repeatability and reproducibility study demonstrate that ALADDIN HW3.0 is substantially equivalent for both repeatability and reproducibility to the predicate device.

The analyses of reproducibility and repeatability shows that ALADDIN HW3.0 and LENSTAR LS 900 performances are substantially equivalent for the assessed endpoints.

17 DECLARATION OF CONFORMITY

DECLARATION OF CONFORMITY/ Dichiarazione di conformità:

Manufacturer/*Fabbricante:*

Name/ <i>Nome</i> :	VISIA imaging S.r.I.
Address/Indirizzo:	Via Martiri della Libertà, 95/e
	52027 San Giovanni Valdarno (AR) – ITALY
me of device/Nome del	dispositivo:

Name of device/Nome del dispositivo:

BIOMETER and CORNEAL ANALYSER integrated

Marca/ <i>Brand</i> :	TOPCON
Model/ <i>Modello</i> :	ALADDIN

I, the undersigned, hereby declare that the aforementioned devices comply with Directive 93/42/EEC (implemented in Italy by Legislative Decree no. 46/97) and subsequent amendments and additions (Directive 2007/47/EC – implemented in Italy by Legislative Decree no. 037/10) for Class IIa equipment.

Io, sottoscritto, con la presente dichiaro che i dispositivi specificati sopra sono conformi alla Direttiva 93/42/CEE (recepita in Italia con D.Lgs 46/97) e successive modifiche e integrazioni (Direttiva 2007/47/CE – recepita in Italia con D.Lgs 037/10) per i dispositivi di Classe IIa.

Alessandro Foggi Managing director *Legal representative*

18 APPENDIX A: INSTALLING AN EXTERNAL PRINTER


18.1 Getting drivers and transferring them to ALADDIN

It is recommended to download the latest drivers for the printer and this should be done through an external PC. Download the latest drivers in the Technical Support section of the printer manufacturer site. Download the drivers, which usually come in a zipped folder. Unzip it and keep the *.inf* file somewhere easily accessible (for example "C:\Drivers"). Copy them to the <u>root directory of an USB FAT32 formatted pen driver</u> that will later be connected to the ALADDIN USB port.

Now go back to ALADDIN and disable the "Write Filter" carrying out the following steps:

18.2 Disabling the Write Filter

- Disable the Write Filter in the Admin tab of Settings panel and Confirm. Settings→Admin→ Disable Filter Press OK
- 2. The machine will restart automatically.

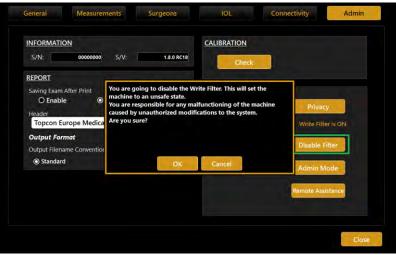



3. Upon restarting, the Aladdin application will warn about the unsafe state of the system.

Press Cancel to avoid the restarting of the machine.





Rev. 19 EN of 10/07/2019

- Close the Aladdin application:
 Settings→Admin→ Close App Press OK
- 5. You will have access to the Desktop of the machine.
- 6. Connect the USB pen with the downloaded drivers to Aladdin.

INFORMATION		CALIBRATION	
S/N: 92150003	S/V:	Check	
IBASE INTEGRATION	á.	ADDITCATION	_
O Enabled Current IBase IP Add	caused by unauthorized	my malfunctioning of the machine modifications to the system.	Write Filter is OFF
IMAGENET O Enabled	O Disabled	Close App	Enable Filter Admin Mode
Current Web Server A			Remote Assistance
http://10.1.3.10:89,	/IMAGEnet/Exam	1	



18.3 Installing a local printer (USB)

The following instructions are for a Windows 10 system.

Close the software using the procedure explained before in step 2. Right click on the Windows start button and select Control Panel.

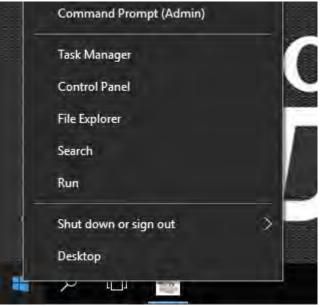


Figure 109

In the control Panel select "Devices and Printers".

All Control Panel Items		-				
↑ 💽 > Control Pane	I > All Control Panel Items >	マ ひ Search Control Panel	م			
Adjust your computer's setti	ngs	View by: Small icons	•			
禮 Administrative Tools	AutoPlay	Backup and Restore (Windows 7)				
Real BitLocker Drive Encryption	Dir Management	Credential Manager				
Pate and Time	Default Programs	🔒 Device Manager				
Devices and Printers	🔲 Display	Base of Access Center				
File Explorer Options	File History	Flash Player (32-bit)				
A Fonts	- HomeGroup	ipsi IDS Camera Manager				

Select the **Add Printer** function.

P Devices and Print	15			- 0	×	
< - ^ ₹	« All Control Panel Items > Devices an	d Printers 🤉 😺	C	Search Devices and Printers	p	
Add a device	dd a printer			E +	6	
V Devices (2)						٨
ALH32-GNZFDO	Generic Non-PnP Monitor					
✓ Printers (4)						

Figure 110

The **Add Printer Wizard** will appear. The printer could be automatically detected. If yes follow the steps, if not press on "The printer that I want isn't listed".

Choose a dev Select a device	ice or printer	to add to	o this PC				
				PI3BB6A0 (H rinter	P LaserJet	MFP M7	25)

Figure 111

 \sim

Select "Add a local printer or network printer with manual settings". Click Next.

	\leftarrow	🖶 Add Printer		
		Find a printer by other options		
		 My printer is a little older. Help me find it. Select a shared printer by name 		Browse
		Example: \\computername\printername or http://computername/printers/printername/.printer		biowaciii
		○ Add a printer using a TCP/IP address or hostname		
		○ Add a Bluetooth, wireless or network discoverable printer		
		Add a local printer or network printer with manual settings		
			Next	Cancel
		Figure 112		
Press Next to the followi	ng	screen.		

÷	🖶 Add Printer		×
	Choose a printer port		
	A printer port is a type of con	nection that allows your computer to exchange inform	nation with a printer.
	Use an existing port:	LPT1: (Printer Port)	~
	○ Create a new port:		
	Type of port:	HP Universal Print Monitor	\sim
			Next Cancel

You now need to identify the make and model of the printer. Click on the Have Disk... button.

÷	🖶 Add	Printer					
	Install	the printer driver					
				e list. Click Windows Update to see more mo installation CD, click Have Disk.	dels.		
	Manuf Brothe Canon		^	Printers Brother Color Leg Type1 Class Driver Brother Color Type3 Class Driver			^
	Dell Eastma	n Kodak Company	~	Brother Color Type4 Class Driver Brother Generic Jpeg Type1 Class Driver			~
		s driver is digitally signer me why driver signing i		Windows Update		ve Disk	
					Next	Cano	cel

Figure 113

Click the **Browse** button and select the folder where you extracted the drivers (.inf file) for this printer. Once you've done that, click **OK**.

Install Fr	om Disk	×
	Insert the manufacturer's installation disk, and then make sure that the correct drive is selected below.	ОК
		Cancel
	Copy manufacturer's files from:	
	A:\ ~	Browse



Select your printer model from the **Printers** list and then click the **Next** button.

Printers	
📮 HP Color LaserJet 1600 Class Driver	
HP Color LaserJet 2500 PCL6 Class Driver	
🔄 HP Color LaserJet 2500 PS Class Driver	
Figure 115	



It is now very important to re-enable the Write Filter. Follow the instructions under <u>"Re-</u> <u>Enabling the Write Filter</u>", otherwise there is the risk of damaging the device.

18.4 Installing a network printer (LAN)

Do not change the Ethernet settings of "Local Area Connection 2" or "**Reserved**" adapter.

If specific configuration for LAN network is needed it the settings can be modified accordingly on "Local Area Connection" or "External" Ethernet adapter.

Connect ALADDIN to an external network through the LAN port using an Ethernet cable. Go to control panel and select devices and printers, refer to previous section. The printer could be automatically detected. If yes follow the steps, if not press on "The printer that I want isn't listed".

Add a device		_		×
Choose a device or printer to ac	ld to this PC			
	NPI3BB6A0 Printer	(HP LaserJet MFP	M725)	
The printer that I want isn't listed		Next	Ca	ancel

Figure 116

Select "Add a printer using a TCP/IP address or hostname", and then select "Standard TCP/IP Port" from the Type of port drop-down menu. Click Next.

÷	🖶 Add Printer
	Find a printer by other options
	○ My printer is a little older. Help me find it.
	○ Select a shared printer by name
	Browse
	Example: \\computername\printername or http://computername/printers/printername/.printer
	Add a printer using a TCP/IP address or hostname
	O Add a Bluetooth, wireless or network discoverable printer
	O Add a local printer or network printer with manual settings
	Next Cancel



Enter the **Printer Name or IP address** in the "Printer Name" textbox. **Port Name** will automatically be entered. You can leave the default. Click **Next**. There will be a small delay while your computer configures the port.

÷	🖶 Add Printer					
	Type a printer hostname o	or IP address				
	Device type:	TCP/IP Device V				
	Hostname or IP address:	10.1.1.36				
	Port name:	10.1.1.36				
	Query the printer and automati	ically select the driver to use				

Figure 118

Next Cancel

From here on, follow the steps for automatic driver selection or select manually the drivers as described for <u>"Installing a local printer (USB)"</u>.

18.5 Re-Enabling the Write Filter

 Open the Aladdin application. The system will warn about the unsafe state. Wait for the automatic restart or press OK. The system will restart enabling the write filter. Once restarted, the system will be in a safe state. 	The preservation and analytic structure of the structure
OR	
	Figure 119
 From inside the application if already open go to: 	General Measurements Surgeons JOL Report Admin INFORMATION CALIBRATION S/N: S/N: Check IBASE INTEGRATION APPLICATION O Enabled Ø Disabled
Settings→Admin→ Enable Filter	Current IBase IP Address 10.1.1.12 Configure Upgrade Upgrade
 After restarting the system will be in a safe state. 	IMAGENET O Enabled Admin Mode O Enabled O Disabled Remote Adsistance Current Web Server Address http://10.1.3.10:89/IMAGEnet/Exam Close App Close Close App Close App Close Close App Close App Close Close App Close App
	Figure 120

19 APPENDIX B: BIBLIOGRAPHY

Please refer to the following literature references on the formulae (in case of specific questions please contact Visia Imaging):

- Haigis: http://www.augenklinik.uni-wuerzburg.de/uslab/ioltxt/haid.htm
- HofferQ: HOFFER KJ: The Hoffer Q formula: A comparison of theoretic and regression formulas. J Cataract Refract Surg, 19:700-712, 1993; ERRATA 20:677, 1994
- Reply: Errata in printed Hoffer Q formula. Journal of Cataract & Refractive Surgery, Volume 33, Issue 1, Pages 2-3, January 2007, Authors: Kenneth J. Hoffer, MD
- Holladay: HOLLADAY JT, PRAGER TC, CHANDLER TY, MUSGROVE KH, LEWIS JW, RUIZ RS: A threepart system for refining intraocular lens power calculations. J Cataract Refract Surg, 14:17-24, 1988
- SRKII: RETZLAFF J: A new intraocular lens calculation formula, Am Intra-Ocular Implant Soc J 6:148-152, 1980
- SRK/T: RETZLAFF J, SANDERS DR, KRAFF MC: Development of the SRK/T intraocular lens implant power calculation formula. J Cataract Refract Surg 16 (3):333-340, 1990
- Barrett Universal II: Barrett Graham D.: An improved universal theoretical formula for intraocular lens power prediction. Journal of Cataract & Refractive Surgery, Volume 19, Issue 6, 713-720
- Olsen: Olsen Thomas et al.: C constant: New concept for ray tracing–assisted intraocular lens power calculation. Journal of Cataract & Refractive Surgery, Volume 40, Issue 5, 764 – 773, 2014
- Olsen Thomas: United States Patent US8657445B2, System and method for determining and predicting IOL power in situ, February 2014
- Olsen: Olsen T., Corydon L., Gimbel H.: Intraocular lens power calculation with an improved anterior chamber depth prediction algorithm. Journal of Cataract & Refractive Surgery, Volume 21, Issue 3, 313-9, 1995

Correction of corneal radii/corneal refraction after corneal refractive surgery:

- HOLLADAY JT: IOL calculations following RK. Refract Corneal Surg 5(3):203, 1989
- HOFFER KJ: Intraocular lens power calculation for eyes after refractive keratotomy. J Refract Surg 11:490:493, 1995

Calculation of phakic implants:

- vd HEIJDE GL, FECHNER PU, WORST JGF: Optische Konsequenzen der Implantation einer negativen Intraokularlinse bei myopen Patienten. Klin MB1 Augenheilk 192:99-102, 1988
- HOLLADAY JT: Refractive power calculations for intraocular lenses in the phakic eye. Am J Ophthalmol 116:63-66, 1993
- HAIGIS W: Biometry in complicated situations, 9th Conv. of DGII 1995, Rochels et al (Hrsg.), Springer, 17-26, 1996

Relations between ultrasound and optical biometer calculation constants:

• RETZLAFF J, SANDERS DR, KRAFF MC (1990): Lens Implant Power Calculation - A manual for ophthalmologists & biometrists, 3rd edition, Slack Inc, Thorofare NJ, USA

- HAIGIS W, LEGE B, MILLER N, SCHNEIDER B: Comparison of immersion ultrasound biometry and partial coherence interferometry for IOL calculation according to Haigis, Graefes Arch Clin Exp Ophthalmology (2000) 238:765-773
- HOLLADAY, JT: International intraocular lens implant registry 2003. J Cataract Refract Surg (2003) 29:176-197
- HAIGIS W: Relations between optimized IOL constants. Symposium on Cataract, IOL and Refractive Surgery of the American Society of Cataract and Refractive Surgery (ASCRS), Philadelphia, PA, USA, June 1-5, 2002, Abstracts, p.112, 2002
- •

Intraocular lens power calculation AFTER corneal refractive surgery:

- Camellin-Calossi: M. Camellin, MD; A. Calossi, Optom "A new formula for intraocular lens power calculation after refractive Corneal Surgery", Journal of Refractive Surgery, vol. 22 Feb. 2006.
- Shammas No-history: SHAMMAS H.J., SHAMMAS M.C: "No-history method of intraocular lens power calculation for cataract surgery after myopic laser in situ keratomileusis", J Cataract Refract Surg 2007; 33:31–36 Q 2007 ASCRS and ESCRS.
- Shammas No-history: SHAMMAS H.J., SHAMMAS M.C., GARABET A., KIM J.H., SHAMMAS A., LABREE L.: Correcting the Corneal Power Measurements for Intraocular Lens Power Calculations After Myopic Laser In Situ Keratomileusis" - American Journal of Ophthalmology (Impact Factor: 4.02). 10/2003; 136(3):426-32.
- Shammas No-history: SHAMMAS H.J., SHAMMAS M.C., HILL W.E.: Intraocular lens power calculation in eyes with previous hyperopic laser in situ keratomileusis" J Cataract Refract Surg 2013; 39:739–744 Q 2013 ASCRS and ESCRS.

Toric IOL calculation:

- HB Fam, KL Lim: Meridional analysis for calculating the expected spherocylindrical refraction in eyes with toric intraocular lenses. Journal of Cataract & Refractive Surgery, 2007 Elsevier
- N Alpins: Astigmatism analysis by the Alpins method. Journal of Cataract & Refractive Surgery, 2001 Elsevier
- G Savini, KJ Hoffer, M Carbonelli, P Ducoli: Influence of axial length and corneal power on the astigmatic power of toric intraocular lenses Journal of Cataract & ..., 2013 Elsevier
- JT Holladay, TV Cravy, DD Koch: Calculating the surgically induced refractive change following ocular surgery. Journal of Cataract & Refractive Surgery, 1992 Elsevier
- Abulafia A, Koch DD, Wang L, Hill WE, Assia El, Franchina M, Barrett GD: New regression formula for toric intraocular lens calculation. Journal of Cataract & Refractive Surgery, 2016 Elsevier
- Alpins N, Barrett GD, Hansen MS, Berdahl JP, Hardten DR, Holladay JT. Innovative toric IOL calculators and how to use them: Barrett Toric Calculator. Cataract Refract Surg Today Europe. May 2015 supplement

ALADDIN HW3.0



Please specify the following when contacting your local supplier regarding questions about this product:

- Product name: Aladdin HW3.0
- Software version and Serial Number: as written in Settings > Admin section
- Period of use: Please inform us of the date of installation
- Defective condition: Please provide us with as much detail as possible

Aladdin HW3.0

User Manual - rev. 19 10/07/2019

Published by:

VISIA imaging S.r.l.



MANUFACTURER VISIA imaging S.r.I. Via Martiri della Libertà, 95/e 52027 San Giovanni Valdarno (AR) Italy